

1390 (Modified)
2000

U.S. DEPARTMENT OF COMMERCE PATENT AND TRADEMARK OFFICE

ATTORNEY'S DOCKET NUMBER

TRANSMITTAL LETTER TO THE UNITED STATES

PA-9912

DESIGNATED/ELECTED OFFICE (DO/EO/US)

U.S. APPLICATION NO. (IF KNOWN, SEE 37 CFR

CONCERNING A FILING UNDER 35 U.S.C. 371

To be assigned 09/914604

INTERNATIONAL APPLICATION NO.

INTERNATIONAL FILING DATE

PRIORITY DATE CLAIMED

PCT/GB00/00916

March 10, 2000

March 12, 1999

TITLE OF INVENTION

Genetic Analysis

APPLICANT(S) FOR DO/EO/US

Michael A. Reeve, Nicholas I. Workman, and Luis Martin-Parras

Applicant herewith submits to the United States Designated/Elected Office (DO/EO/US) the following items and other information:

1. ☒ This is a **FIRST** submission of items concerning a filing under 35 U.S.C. 371.
2. ☐ This is a **SECOND** or **SUBSEQUENT** submission of items concerning a filing under 35 U.S.C. 371.
3. ☒ This is an express request to begin national examination procedures (35 U.S.C. 371(f)). The submission must include items (5), (6), (9) and (24) indicated below.
4. ☒ The US has been elected by the expiration of 19 months from the priority date (Article 31).
5. ☒ A copy of the International Application as filed (35 U.S.C. 371 (c) (2))
 - a. ☒ is attached hereto (required only if not communicated by the International Bureau).
 - b. ☒ has been communicated by the International Bureau.
 - c. ☐ is not required, as the application was filed in the United States Receiving Office (RO/US).
6. ☐ An English language translation of the International Application as filed (35 U.S.C. 371(c)(2)).
 - a. ☐ is attached hereto.
 - b. ☐ has been previously submitted under 35 U.S.C. 154(d)(4).
7. ☐ Amendments to the claims of the International Application under PCT Article 19 (35 U.S.C. 371 (c)(3))
 - a. ☐ are attached hereto (required only if not communicated by the International Bureau).
 - b. ☐ have been communicated by the International Bureau.
 - c. ☐ have not been made; however, the time limit for making such amendments has NOT expired.
 - d. ☐ have not been made and will not be made.
8. ☐ An English language translation of the amendments to the claims under PCT Article 19 (35 U.S.C. 371(c)(3)).
9. ☒ An oath or declaration of the inventor(s) (35 U.S.C. 371 (c)(4)).
10. ☐ An English language translation of the annexes of the International Preliminary Examination Report under PCT Article 36 (35 U.S.C. 371 (c)(5)).
11. ☐ A copy of the International Preliminary Examination Report (PCT/IPEA/409).
12. ☒ A copy of the International Search Report (PCT/ISA/210).

Items 13 to 20 below concern document(s) or information included:

13. ☒ An Information Disclosure Statement under 37 CFR 1.97 and 1.98.
14. ☐ An assignment document for recording. A separate cover sheet in compliance with 37 CFR 3.28 and 3.31 is included.
15. ☒ A **FIRST** preliminary amendment.
16. ☐ A **SECOND** or **SUBSEQUENT** preliminary amendment.
17. ☐ A substitute specification.
18. ☐ A change of power of attorney and/or address letter.
19. ☐ A computer-readable form of the sequence listing in accordance with PCT Rule 13ter.2 and 35 U.S.C. 1.821 - 1.825.
20. ☐ A second copy of the published international application under 35 U.S.C. 154(d)(4).
21. ☐ A second copy of the English language translation of the international application under 35 U.S.C. 154(d)(4).
22. ☒ Certificate of Mailing by Express Mail
23. ☒ Other items or information:

copy of this transmittal letter for charging purposes
return postcard

APPLICATION NO. (IF KNOWN, SEE 37 CFR

INTERNATIONAL APPLICATION NO.

ATTORNEY'S DOCKET NUMBER

To be assigned 4604

PCT/GB00/00916

PA-9912

24. The following fees are submitted:

CALCULATIONS PTO USE ONLY

BASIC NATIONAL FEE (37 CFR 1.492 (a) (1) - (5)) :

- ☐ Neither international preliminary examination fee (37 CFR 1.482) nor international search fee (37 CFR 1.445(a)(2)) paid to USPTO and International Search Report not prepared by the EPO or JPO **\$1000.00**
- ☒ International preliminary examination fee (37 CFR 1.482) not paid to USPTO but International Search Report prepared by the EPO or JPO **\$860.00**
- ☐ International preliminary examination fee (37 CFR 1.482) not paid to USPTO but international search fee (37 CFR 1.445(a)(2)) paid to USPTO **\$710.00**
- ☐ International preliminary examination fee (37 CFR 1.482) paid to USPTO but all claims did not satisfy provisions of PCT Article 33(1)-(4) **\$690.00**
- ☐ International preliminary examination fee (37 CFR 1.482) paid to USPTO and all claims satisfied provisions of PCT Article 33(1)-(4) **\$100.00**

ENTER APPROPRIATE BASIC FEE AMOUNT =**\$860.00**

Surcharge of **\$130.00** for furnishing the oath or declaration later than ☐ 20 ☐ 30 months from the earliest claimed priority date (37 CFR 1.492 (e)).

\$0.00

CLAIMS	NUMBER FILED	NUMBER EXTRA	RATE
Total claims	24 - 20 =	4	x \$18.00
Independent claims	3 - 3 =	0	x \$84.00
Multiple Dependent Claims (check if applicable).			<input type="checkbox"/>

\$72.00**\$0.00****\$0.00****TOTAL OF ABOVE CALCULATIONS =****\$932.00**

- ☐ Applicant claims small entity status. (See 37 CFR 1.27). The fees indicated above are reduced by 1/2.

\$0.00**SUBTOTAL =****\$932.00**

Processing fee of **\$130.00** for furnishing the English translation later than ☐ 20 ☐ 30 months from the earliest claimed priority date (37 CFR 1.492 (f)).

\$0.00**TOTAL NATIONAL FEE =****\$932.00**

Fee for recording the enclosed assignment (37 CFR 1.21(h)). The assignment must be accompanied by an appropriate cover sheet (37 CFR 3.28, 3.31) (check if applicable).

\$0.00**TOTAL FEES ENCLOSED =****\$932.00**

Amount to be:

refunded

\$

charged

\$

- a. ☐ A check in the amount of _____ to cover the above fees is enclosed.
- b. ☒ Please charge my Deposit Account No. **500-588** in the amount of **\$932.00** to cover the above fees. A duplicate copy of this sheet is enclosed.
- c. ☒ The Commissioner is hereby authorized to charge any additional fees which may be required, or credit any overpayment to Deposit Account No. **500-588**. A duplicate copy of this sheet is enclosed.
- d. ☐ Fees are to be charged to a credit card. **WARNING:** Information on this form may become public. **Credit card information should not be included on this form.** Provide credit card information and authorization on PTO-2038.

NOTE: Where an appropriate time limit under 37 CFR 1.494 or 1.495 has not been met, a petition to revive (37 CFR 1.137(a) or (b)) must be filed and granted to restore the application to pending status.

SEND ALL CORRESPONDENCE TO:

Royal N. Ronning, Jr.
Amersham Pharmacia Biotech, Inc.
800 Centennial Avenue
Piscataway, New Jersey 08855

(732) 457-8423

SIGNATURE

Royal N. Ronning, Jr.

NAME

32,529

REGISTRATION NUMBER

August 28, 2001

DATE

Information upon delivery.
Use telephone:

914604



EL717698470US

EL717698470US



UNITED STATES POSTAL SERVICE®


POST OFFICE TO ADDRESSEE

IN (POSTAL USE ONLY)

Code: _____
Day of delivery: _____
Time: _____
Address: _____
City: _____
State: _____
Zip: _____
Weight: _____
Dimensions: _____
Special services: _____
Insurance: _____
Signature: _____

DELIVERY (POSTAL USE ONLY)


Delivery attempt: _____
Employee signature: _____
Day of delivery: _____
Time: _____
Address: _____
City: _____
State: _____
Zip: _____
Weight: _____
Dimensions: _____
Special services: _____
Insurance: _____
Signature: _____
☐ WAIVER OF SIGNATURE (Domestic Only): Additional merchandise insurance is required. If waiver of signature is requested, I wish delivery to be made without obtaining signature of addressee or addressee's agent (if delivery employee judges that article can be left in secure location) and I authorize that delivery employee's signature constitutes valid proof of delivery.
NO DELIVERY ☐ Weekend ☐ Holiday

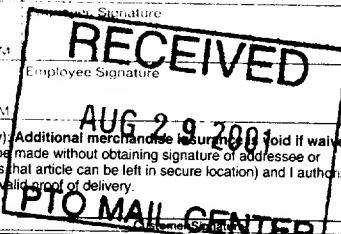

IF YOU WANT TO FILE A CLAIM FOR DAMAGE OR LOSS OF CONTENTS, YOU MUST PRESENT THE ARTICLE, CONTAINER, AND PACKAGING TO THE POST OFFICE FOR INSPECTION.

FROM: (PLEASE PRINT) _____ PHONE: () _____

TO: (PLEASE PRINT) _____ PHONE: () _____

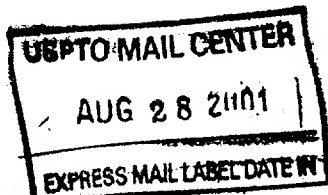

FOR PICKUP OR TRACKING CALL 1-800-222-1811 www.usps.com **EMS**

Addressee Copy
Label 11-F August 2000

EMCA or government agency, fill in account number on label.

move label backing and adhere over these instructions.
sure to remove the Express Mail number from the backing and
in for your records.

1-800-222-1811 for our convenient pick-up service
one low fee no matter how many pieces) or drop off your
Express Mail package at a post office, or an Express Mail box.

Rec'd PCT/PTO 31 JAN 2002

PA-9912

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Application of:	M. Reeve, et al.	Group Art Unit:	To be assigned
Serial Number:	09/914,604	Examiner:	To be assigned
Filing Date:	August 28, 2001		
Title:	Genetic Analysis		

SECOND PRELIMINARY AMENDMENT

Assistant Commissioner for Patents
Box PCT
Washington, D.C. 20234

Sir:

Please consider the following amendments and remarks in connection with the prosecution of the captioned application.

In the Specification


At the end of the written description, before the claims, please insert the Sequence Listing attached hereto.

Remarks

Applicants have amended the specification by insertion of a sequence listing. The sequence listing has also been submitted under separate cover in an initial computer readable form as well as on a paper copy in response to the "Notification of Missing Requirements Under 35 U.S.C. 371 in the United States Designated/Elected Office (DO/EO/US)" having a mailing date of October 25, 2001.

Early and favorable consideration is requested.

Respectfully submitted,

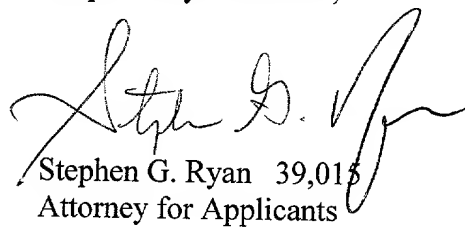
Stephen G. Ryan 39,015
Attorney for Applicants

Amersham Biosciences
800 Centennial Avenue
P. O. Box 1327
Piscataway, New Jersey 08855-1327

Tel: (732) 457-8071
Fax: (732) 457-8463

SEQUENCE LISTING

<110> Reeve, Michael Alan
Workman, Nicholas Ian
Martin-Parras, Luis

<120> Genetic Analysis

<130> PA9912

<140> 09/914,604

<141> 2001-08-28

<150> PCT/GB00/00916

<151> 2000-03-10

<150> EPO 99301933.0

<151> 1999-03-12

<160> 32

<170> PatentIn Ver. 2.1

<210> 1

<211> 3954

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Plasmid

<400> 1

```
tgcgcgcttt cggatgatgac ggtgaaaacc tctgacacat gcagctcccg gagacgggtca 60
cagcttgtct gtaagcggat gccgggagca gacaagcccg tcagggcgcg tcagcgggtg 120
ttggcggttg tcggggctgg cttaactatg cggcatcaga gcagattgta ctgagagtgc 180
accatatgcg gtgtgaaata ccgcacagat gcgtaaggag aaaataccgc atcaggcgcc 240
attcgccatt caggctgcgc aactgttggg aaggcgatc ggtgcgggccc tcttcgctat 300
tacgccagct ggcgaaaggg ggatgtgctg caaggcgatt aagttgggta acgccagggt 360
tttcccagtc acgacgttgt aaaacgacgg ccagtgaatt cgagctcggg accggggccc 420
ccctcgaggc cgacgggtatc gataagcttg atcgagctg gtaatccgga cggccgcgtc 480
gaagatgttg ggggtgttgta acaatatcga ttccaattca gggggggcca cctgatatcc 540
tttgtattta attaaagact tcaagcggtc aactatgaag aagtgttcgt cttcgtccca 600
gtaaggatcc gcactttgaa ttttgtaatc ctgaagggat cgtaaaaaa gctcttcttc 660
aaatctatac attaaagcga ctcgaaatcc acatatcaaa tatccgagt tagtaaaca 720
tccaaaaccg tgatggaatg gaacaacact taaatgtaca ccctggtaat ccgttttaga 780
atccatgata ataattttct ggattatttg taattttttt tgcacgttca aaattttttg 840
caaccctttt ttggaacaaa acactacggg aggctgcgaa atgttcatac tgttgagcaa 900
```

ttcacgttca ttataagctt ttcactgcat acgacgattc tgtgatttgt attcagccca 960
 tatcgtttca tagcttctgc caaccgaacg gacatttctga agtattccgc gtacgtgatg 1020
 ttcacctcga tatgtgcata tgtaaaagca attgttccag gaaccagggc gtatctcttc 1080
 atagccatgg aatacgctt tttcagtgtt gcgatgctaa tgccgttaca aataatccga 1140
 gcaccaagaa tggctgcgcg ctgctctggt acttgacgtc gtatttgacg gggtccttga 1200
 gaaagtattt aaactggaac acaatctgag gaatgatcaa agcaaccaac gccaacgcac 1260
 aataactagt gcaataccaa gacctccaa taatagcacc cagacttggt taataacctc 1320
 tggctctgat attgctccag atggaattgg acgatatggc tcattaattg cgtcgatata 1380
 tctatcatac cagtcgttga ttgtctgtgt atagccagta agacaaggac cagacatcat 1440
 catgcaaaga atcgcttaag ccttcttggt cctttatgag gatctctctg attttcttg 1500
 cgtcgagttt tccggttaaga cctttcggtc cttcgtccac aaacacaact cctccgcgca 1560
 acttttctgc ggttggtact tgactggcca cgtaatccac gatctctttt tccgtcatcg 1620
 tctttccgtg ctccaaaaca acaacggcgg cgggtccgga ttaccagctg cgatcaagct 1680
 tatcgatacc gtcgacctcg acctgcaggc atgcaagctt ggcgtaatca tggcatagc 1740
 gtttctctgt gtgaaattgt tatccgctca caattccaca caacatacga gccggaagca 1800
 taaagtgtaa agcctggggt gcctaattgag tgagctaaact cacattaatt gcgttgcgct 1860
 cactgcccgc tttccagtcg ggaaacctgt cgtgccagct gcattaatga atcggccaac 1920
 gcgcggggag aggcggtttg cgtattgggc gctcttccgc ttcctcgctc actgactcgc 1980
 tgcgctcggt cgttcggctg cggcgagcgg tatcagctca ctcaaaggcg gtaatacgtt 2040
 tatccacaga atcaggggat aacgcaggaa agaactgtg agcaaaaggc cagcaaaagg 2100
 ccaggaaccg taaaaaggcc gcgttgctgg cgtttttcca taggctccgc cccctgacg 2160
 agcatcacia aaatcgacgc tcaagtcaga ggtggcgaaa cccgacagga ctataaagat 2220
 accaggcggt tccccctgga agctccctcg tgcgctctcc tgttccgacc ctgccgctta 2280
 ccggatacct gtccgccttt ctcccttcgg gaagcgtggc gcttctcat agctcacgct 2340
 gtaggtatct cagttcgggt taggtcgttc gctccaagct gggctgtgtg cacgaacccc 2400
 ccgttcagcc cgaccgctgc gccttatccg gtaactatcg tcttgagtcc aaccggtaa 2460
 gacacgactt atcgccactg gcagcagcca ctggtaacag gattagcaga gcgaggtagt 2520
 taggcgggtg tacagagttc ttgaagtggg ggcctaacta cggctacact agaaggacag 2580
 tatttggtat ctgcgctctg ctgaagccag ttaccttcgg aaaaagagtt ggtagctctt 2640
 gatccggcaa acaaacacc gctggtagcg gtgggttttt tgttgcaag cagcagatta 2700
 cgcgcagaaa aaaaggatct caagaagatc ctttgatctt ttctacgggg tctgacgctc 2760
 agtggaaacga aaactcacgt taagggattt tggatcatgag attatcaaaa aggatcttca 2820
 cctagatcct tttaaattaa aaatgaagtt ttaaataaat ctaaagtata tatgagtaaa 2880
 cttgggtctga cagttaccaa tgcttaatca gtgaggcacc tatctcagcg atctgtctat 2940
 ttogttcatc catagttgcc tgactccccg tctgttagat aactacgata cgggaggggt 3000
 taccatctgg cccagtgct gcaatgatac cgcgagaccc acgctcacgg gctccagatt 3060
 tatcagcaat aaaccagcca gccggaaggg ccgagcgcag aagtggctct gcaactttat 3120
 ccgcctccat ccagtctatt aattgttgcc ggaagctag agtaagtagt tccgcagtta 3180
 atagtttgcg caacgttggt gccattgcta caggcatcgt ggtgtcacgc tctcgtttt 3240
 gtatggcttc attcagctcc ggttcccaac gatcaaggcg agttacatga tccccatgt 3300
 tgtgcaaaaa agcgggttagc tccttcggct ctccgatcgt tgtcagaagt aagttggccg 3360
 cagtgttatc actcatggtt atggcagcac tgcataattc tcttactgtc atgccatccg 3420
 taagatgctt ttctgtgact ggtgagtact caaccaagtc attctgagaa tagtgtatgc 3480
 ggcgaccgag ttgctcttgc ccggcgtcaa tacgggataa taccgcgcca catagcagaa 3540
 ctttaaaagt gctcatcatt ggaaaacgtt cttcggggcg aaaactctca aggatcttac 3600
 cgctgttgag atccagttcg atgtaaccca ctctgacacc caactgatct tcagcatctt 3660
 ttactttcac cagcgtttct ggtgagcaa aaacaggaag gcaaaatgcc gcaaaaaagg 3720
 gaataagggc gacacggaaa tgttgaatac tcatactctt cctttttcaa tattattgaa 3780

gcatttatca gggttattgt ctcatgagcg gatacatatt tgaatgtatt tagaaaaata 3840
aacaaatagg gggtccgcgc acatttcccc gaaaagtgcc acctgacgtc taagaaacca 3900
ttattatcat gacattaacc tataaaaaata ggcgtatcac gaggcccttt cgtc 3954

<210> 2

<211> 37

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 2

cgcagctggg aatccggacg cccgcgtcga agatggt 37

<210> 3

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 3

cgcagctggg aatccggacc cgcgcgctt gttggt 36

<210> 4

<211> 2946

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Plasmid

<400> 4

atgaccatga ttacgccaaag ctctaatacg actcactata gggaaagctt cccgacgtct 60
caggctaatag ttggcccacc gacgttccac gatggggcgc tottaagggc ttagaccctc 120
gtcgggagta tttctgtgat ctggcgacac tcacgcgaga agtcattacc ggcgatatga 180
attcactggc cgtcgtttta caacgtcgtg actgggaaaa ccctggcgtt acccaactta 240
atcgccttgc agcacatccc cctttcgcca gctggcgtaa tagcgaagag gcccgccaccg 300
atcgcccttc ccaacagttg cgcagcctga atggcgaatg ggaaattgta aacgttaata 360
ttttgttaaa attcgcgtta aatttttgtt aaatcagctc attttttaac caataggccg 420

```

aaatcggcaa aatcccttat aaatcaaaag aatagaccga gataggggtg agtgttgttc 480
cagtttgtaa caagagtcca ctattaaaga acgtggactc caacgtcaaa gggcgaaaaa 540
ccgtctatca gggcgatggc ccactacgtg aaccatcacc ctaatcaagt tttttggggg 600
cgaggtgccg taaagcacta aatcgggaacc ctaaagggag ccccgattt agagcttgac 660
ggggaaagcc ggcgaacgtg gcgagaaagg aagggaagaa agcgaaagga gcgggcgcta 720
gggcgctggc aagtgtagcg gtcacgctgc gcgtaaccac cacaccgccc gcgcttaatg 780
cgccgctaca gggcgctca ggtggcactt ttcggggaaa tgtgcgcgga acccctattt 840
gtttattttt ctaaatacat tcaaataatgt atccgctcat gagacaataa ccctgataaa 900
tgcttcaata atattgaaaa aggaagagta tgagtattca acatttccgt gtcgccctta 960
ttcccttttt tgcggcattt tgccttccctg tttttgtca cccagaaacg ctggtgaaag 1020
taaaagatgc tgaagatcag ttgggtgcac gagtgggtta catcgaactg gatctcaaca 1080
gcgtaagat ccttgagagt tttcgccccg aagaacgttt tccaatgatg agcactttta 1140
aagtcttgct atgtggcgcg gtattatccc gtattgacgc cgggcaagag caactcggtc 1200
gccgcataca ctattctcag aatgacttgg ttgagtactc accagtcaca gaaaagcatc 1260
ttacggatgg catgacagta agagaattat gcagtgtcgc cataaccatg agtgataaca 1320
ctgcggccaa cttacttctg acaacgatcg gaggaccgaa ggagctaacc gcttttttgc 1380
acaacatggg ggatcatgta actcgccttg atcggtggga accggagctg aatgaagcca 1440
taccaaacga cgagcgtgac accacgatgc ctgtagcaat ggcaacaacg ttgcgcaaac 1500
tattaactgg cgaactactt actctagctt cccggcaaca attaatagac tggatggagg 1560
cggataaagt tgcaggacca cttctgcgct cggcccttcc ggctggctgg tttattgctg 1620
ataaatctgg agccggtgag cgtgggtctc gcggtatcat tgcagcactg gggccagatg 1680
gtaagccctc ccgtatcgta gttatctaca gcacggggag tcaggcaact atggatgaac 1740
gaaatagaca gatcgctgag ataggtgcct cactgattaa gcatttgtaa ctgtcagacc 1800
aagtttactc atatatactt tagattgatt taaaacttca tttttaattt aaaaggatct 1860
aggtgaagat cttttttgat aatctcatga caaaaatccc ttaacgtgag ttttcgttcc 1920
actgagcgtc agaccccgta gaaaagatca aaggatcttc ttgagatcct ttttttctgc 1980
gcgtaatctg ctgcttgcaa acaaaaaaac caccgctacc agcgtggtt tgtttgcggg 2040
atcaagagct accaactctt tttccgaagg taactggctt cagcagagcg cagataccaa 2100
atactgtcct tctagtgtag ccgtagttag gccaccactt caagaactct gtagcaccgc 2160
ctacatacct cgctctgcta atcctgttac cagtggctgc tgccagtggc gataagtcgt 2220
gtcttaccgg gttggactca agacgatagt taccggataa ggcgcagcgg tcgggctgaa 2280
cgggggggtt gtgcacacag cccagcttgg agcgaacgac ctacaccgaa ctgagatacc 2340
tacagcgtga gctatgagaa agcgccacgc ttcccgaagg gagaaaggcg gacaggatatc 2400
cggtaagcgg cagggtcgga acaggagagc gcacgaggga gcttccaggg ggaaacgcct 2460
ggtatcttta tagtcctgtc ggggttctgcc acctctgact tgagcgtcga tttttgtgat 2520
gctcgtcagg ggggcggagc ctatggaaaa acgccagcaa cgcggccttt ttacggttcc 2580
tggccttttg ctggcctttt gctcacatgt tctttcctgc gttatccctt gattctgtgg 2640
ataacggtat taccgccttt gagtgagctg ataccgctcg ccgcagccga acgaccgagc 2700
gcagcagctc agtgagcgag gaagcggaag agcgcccaat acgcaaaccg cctctccccg 2760
cgcgttggcc gattcattaa tgcagctggc acgacaggtt tcccgactgg aaagcgggca 2820
gtgagcgcaa cgcaattaat gtgagttagc tcactcatta ggcaccccag gctttacact 2880
ttatgcttcc ggctcgtatg ttgtgtggaa ttgtgagcgg ataacaattt cacacaggaa 2940
acagct
2946

```

<210> 5

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 5

gttttcccag tcacgacgt

19

<210> 6

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 6

tccggacgtc tcaggctaatt gtt

23

<210> 7

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 7

tgtaacgaca cattgctgga tacc

24

<210> 8

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 8

atataactct cgctccttga taac

24

<210> 9
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 9
aggcgtctga ggctgcggct atgg 24

<210> 10
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 10
aaccgcgtcgc gacgagagtc taag 24

<210> 11
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 11
gatatacgtg atatattttg attg 24

<210> 12
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 12
gatcgggtatc cagcaatgtg tcggttaca 28

<210> 13
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 13
agctgttatc aaggagcgag agttatat 28

<210> 14
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 14
catgccatag ccgcagcctc agacgcct 28

<210> 15
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 15
ctagcttaga ctctcgtcgc gacgggtt 28

<210> 16
<211> 28
<212> DNA
<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 16

ttaacaatca aaatatatca cgtatatc

28

<210> 17

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 17

tgtaacgaca cattgctgga taccgatcc

29

<210> 18

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 18

atataactct cgtccttga taacagctt

29

<210> 19

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 19

aggcgtctga ggctgcggct atggcatgg

29

<210> 20
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 20
aaccgcgcgc gacgagagtc taagctagt

29

<210> 21
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 21
gatatacgtg atatattttg attgttaag

29

<210> 22
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 22
gcgctagggg tccttactgg gacgaagacg aa

32

<210> 23
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 23
gcgctagga tccttactgg gacgaagacg aa

32

<210> 24
<211> 93
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 24
cccggggat cctcgtttta ttgggccgag ttttggtccg tagtgcttgg ttagatatgc 60
ttatgttcac aaaatcatcc ttgtacagaa ttc 93

<210> 25
<211> 93
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 25
cccggggat cctcgtttta ttgggccgag ttttggtccg tagtgcatgg ttagatatgc 60
ttatgttcac aaaatcatcc ttgtacagaa ttc 93

<210> 26
<211> 93
<212> DNA
<213> Artificial Sequence

<220>
<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 26
cccggtgta cacaaaagtt tacctgaaga acgtgggggg tcgtgcctgg tcttgcgta 60
cctggtctca ggagaggggc cccatgggaa ttc 93

<210> 27
<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 27

ctactgatcg gatccccg

18

<210> 28

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 28

aaacgacggc cagtgaat

18

<210> 29

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 29

ctactgatcg gatccccg

18

<210> 30

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 30

aaacgacggc cagtgaat

18

<210> 31
<211> 20
<212> DNA
<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 31
ggccgagttt tgggccgtag

20

<210> 32
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: Synthetic
Oligonucleotide

<400> 32
gtcttgcgtc acctggcttc ag

22

PA-9912

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Application of: M. Reeve, et al. Group Art Unit: To be assigned
Serial Number: To be assigned Examiner: To be assigned
Filing Date: August 28, 2001
Title: Genetic Analysis

FIRST PRELIMINARY AMENDMENT

Honorable Assistant Commissioner of Patents
Box Patent Application
Washington, D.C. 20231

Sir:

Please consider the following amendments and remarks in connection with the prosecution of the captioned application, which is a filing under 35 U.S.C. § 371 and claims priority to international application number PCT/GB00/00916 filed March 10, 2000. This application also claims priority to application number 99301933.0 filed in the European Patent Office on March 12, 1999.

In the Claims

Please amend page 155, line 1 as follows:

[CLAIMS]

What is claimed is:

Please cancel claims 3, 12, 13, 19, 20, and 26, without prejudice.

Please amend claim 1 as follows:

1. (once amended) [A]In a method of providing a mixture of DNA fragments enriched in fragments that are characteristic of a phenotype of interest, [by]which method includes providing affected DNA in fragmented form and providing unaffected DNA in fragmented form, [which method comprises]the improvement comprising:
 - a) mixing the fragments of the affected DNA and the fragments of the unaffected DNA under hybridising conditions to form hybrids;
 - b) recovering a mixture of hybrids that contain mismatches;
 - c) recovering fragments of the affected DNA from the mixture of hybrids that contain mismatches[;and optionally repeating steps a), b) and c) one or more times].

Please amend claim 2 as follows:

2. (once amended) The method of claim 1 wherein the affected DNA is pooled DNA of one or more individuals who show the phenotype of interest, and the unaffected DNA is pooled DNA of one or more individuals who do not show the phenotype of interest.

Please amend claim 6 as follows:

6. (once amended) The method of [any one of claims 1 to 5]claim 1, wherein step b) is performed by use of a mismatch-binding protein.

Please amend claim 7 as follows:

7. (once amended) The method of [any one of claims 1 to 6]claim 1, wherein either the fragments of the affected DNA or the fragments of the unaffected DNA are tagged by one member of a specific binding pair, and step c) is performed by using the other member of the specific binding pair.

Please amend claim 9 as follows:

9. (once amended) The method of [any one of claims 1 to 8]claim 1,
[wherein]further comprising subjecting the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest[, is subjected] to self-hybridisation [followed by recovery of]to form duplexes and subsequently recovering the perfectly matched duplexes.

Please amend claim 10 as follows:

10. (once amended) The method of [any one of claims 1 to 9]claim 1,
[wherein]further comprising mixing the mixture of DNA fragments enriched in
fragments that are characteristic of the phenotype of interest[, is mixed] with an
excess of the fragments of the affected DNA under hybridisation conditions[,]to
form duplexes and subsequently recovering the followed by recovery of]
perfectly matched duplexes.

Please amend claim 11 as follows:

11. (once amended) The method of [any one of claims 1 to 10]claim 1, wherein each
of the affected DNA and the unaffected DNA is provided in fragmented form by
digestion with from 4 to 7 six-cutter restriction endonuclease enzymes together
with from 0 to 50 four-cutter restriction endonuclease enzymes.

Please amend claim 15 as follows:

15. (once amended) The method of [claim 13 or]claim 14, wherein steps a to f) are
repeated using each different subset of r restriction endonuclease enzymes to give
 $(n!)/[(n-r)!r!]$ different arrays.

Please amend claim 16 as follows:

16. (once amended) The method of [any one of claims 13 to 15]claim 14, wherein the n restriction endonuclease enzymes are selected from 4-cutters and 5-cutters and 6-cutters.

Please amend claim 17 as follows:

17. (once amended) The method of [any one of claims 13 to 16]claim 14, wherein the n is 3 to 10 and r is 2 to 4.

Please amend claim 21 as follows:

21. (once amended) [The]A set of arrays [of claim 19 or claim 20]produced by the method of claim 14, derived from a set of n = 6 six-cutter restriction endonuclease enzymes which are *BamHI*; *Bsr GI*; *Hind III*; *NcoI*; *SpeI*; and *AflIII*.

Please amend claim 22 as follows:

22. (once amended) [The]A set of arrays [of claim 19 or claim 20]produced by the method of claim 14, derived from the set of n = 6 six-cutter restriction

endonuclease enzymes which are *EcoRI*; *BspHI*; *BglIII*; *XbaI*; *Acc65I*; and *ApaLI*.

Please amend claim 23 as follows:

23. (once amended) A nucleic acid characterisation method which comprises presenting to [the]a set of arrays [of any one of claims 19 to 22]produced by the method of claim 14 a nucleic acid fragment of interest under hybridisation conditions, and observing a pattern of hybridisation.

Please amend claim 28 as follows:

28. (once amended) The double-stranded DNA molecule of claim 27, wherein[the following criteria are satisfied]:
- a) inter-fragment length differences are greater for larger fragments;
 - b) all possible fragments are unambiguously resolvable by electrophoresis from one another;
 - c) size gaps between bands comprising different numbers of inter-restriction-site units are larger than size gaps between bands comprising the same number of inter-restriction-site units;
 - d) the size gaps and size spread from the largest to the smallest fragment are electrophoretically compatible.

Please add new claim 29 as follows:

29. (new) The method of claim 1 which further comprises steps a), b) and c) one or more times.

Please add new claim 30 as follows:

30. (new) The method of claim 14 wherein one of the restriction enzyme specific primers is tagged.


Remarks

Claims 1-28 are pending in the instant application. Applicants have amended claims 1, 2, 6, 7, 9, 10, 11, 15, 16, 17, 21, 22, 23, and 28 to more fully conform with U.S. practice and to delete multiple dependencies. Applicants have cancelled, without prejudice, claims 3, 12, 13, 19, 20, and 26. Applicants have also added new claims 29 and 30. A version of the claims marked up to show the amendments, as well as a clean version of the claims encompassing the amendments, is attached hereto.

Applicants respectfully assert that all amendments are fairly based on the specification, and respectfully request their entry.

Applicants believe that the claims, as amended, are in allowable form, and earnestly solicit the allowance of claims 1, 2, 4-11, 14-18, 21-25, and 27-30.

Respectfully submitted,

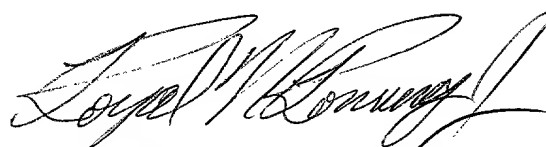

Royal N. Ronning, Jr. 32,529
Attorney for Applicants

Amersham Pharmacia Biotech, Inc.
800 Centennial Avenue
P. O. Box 1327
Piscataway, New Jersey 08855-1327

Tel: (732) 457-8423
Fax: (732) 457-8463

Claims (marked-up version showing amendment(s))

[CLAIMS]

What is claimed is:

1. (once amended) [A]In a method of providing a mixture of DNA fragments enriched in fragments that are characteristic of a phenotype of interest, [by]which method includes providing affected DNA in fragmented form and providing unaffected DNA in fragmented form, [which method comprises]the improvement comprising:
 - a) mixing the fragments of the affected DNA and the fragments of the unaffected DNA under hybridising conditions to form hybrids;
 - b) recovering a mixture of hybrids that contain mismatches;
 - c) recovering fragments of the affected DNA from the mixture of hybrids that contain mismatches[;and optionally repeating steps a), b) and c) one or more times].
2. (once amended) The method of claim 1 wherein the affected DNA is pooled DNA of one or more individuals who show the phenotype of interest, and the unaffected DNA is pooled DNA of one or more individuals who do not show the phenotype of interest.

6. (once amended) The method of [any one of claims 1 to 5]claim 1, wherein step b) is performed by use of a mismatch-binding protein.
7. (once amended) The method of [any one of claims 1 to 6]claim 1, wherein either the fragments of the affected DNA or the fragments of the unaffected DNA are tagged by one member of a specific binding pair, and step c) is performed by using the other member of the specific binding pair.
9. (once amended) The method of [any one of claims 1 to 8]claim 1,
[wherein]further comprising subjecting the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest[, is subjected] to self-hybridisation [followed by recovery of]to form duplexes and subsequently recovering the perfectly matched duplexes.
10. (once amended) The method of [any one of claims 1 to 9]claim 1,
[wherein]further comprising mixing the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest[, is mixed] with an excess of the fragments of the affected DNA under hybridisation conditions[,]to form duplexes and subsequently recovering the[followed by recovery of] perfectly matched duplexes.

11. (once amended) The method of [any one of claims 1 to 10]claim 1, wherein each of the affected DNA and the unaffected DNA is provided in fragmented form by digestion with from 4 to 7 six-cutter restriction endonuclease enzymes together with from 0 to 50 four-cutter restriction endonuclease enzymes.
15. (once amended) The method of [claim 13 or]claim 14, wherein steps a to f) are repeated using each different subset of r restriction endonuclease enzymes to give $(n!)/[(n-r)!r!]$ different arrays.
16. (once amended) The method of [any one of claims 13 to 15]claim 14, wherein the n restriction endonuclease enzymes are selected from 4-cutters and 5-cutters and 6-cutters.
17. (once amended) The method of [any one of claims 13 to 16]claim 14, wherein the n is 3 to 10 and r is 2 to 4.
21. (once amended) [The]A set of arrays [of claim 19 or claim 20]produced by the method of claim 14, derived from a set of n = 6 six-cutter restriction endonuclease enzymes which are *BamHI*; *Bsr GI*; *Hind III*; *NcoI*; *SpeI*; and *AflIII*.
22. (once amended) [The]A set of arrays [of claim 19 or claim 20]produced by the

method of claim 14, derived from the set of $n = 6$ six-cutter restriction endonuclease enzymes which are *EcoRI*; *BspHI*; *BglII*; *XbaI*; *Acc65I*; and *ApaLI*.

23. (once amended) A nucleic acid characterisation method which comprises presenting to [the] a set of arrays [of any one of claims 19 to 22] produced by the method of claim 14 a nucleic acid fragment of interest under hybridisation conditions, and observing a pattern of hybridisation.
28. (once amended) The double-stranded DNA molecule of claim 27, wherein[the following criteria are satisfied]:
- a) inter-fragment length differences are greater for larger fragments;
 - b) all possible fragments are unambiguously resolvable by electrophoresis from one another;
 - c) size gaps between bands comprising different numbers of inter-restriction-site units are larger than size gaps between bands comprising the same number of inter-restriction-site units;
 - d) the size gaps and size spread from the largest to the smallest fragment are electrophoretically compatible.
29. (new) The method of claim 1 which further comprises steps a), b) and c) one or more times.

30. (new) The method of claim 14 wherein one of the restriction enzyme specific primers is tagged.

2014-04-04 14:04:04

Claims (clean version encompassing amendments)

What is claimed is:

1. (once amended) In a method of providing a mixture of DNA fragments enriched in fragments that are characteristic of a phenotype of interest, which method includes providing affected DNA in fragmented form and providing unaffected DNA in fragmented form, the improvement comprising:
 - a) mixing the fragments of the affected DNA and the fragments of the unaffected DNA under hybridising conditions to form hybrids;
 - b) recovering a mixture of hybrids that contain mismatches;
 - c) recovering fragments of the affected DNA from the mixture of hybrids that contain mismatches.
2. (once amended) The method of claim 1 wherein the affected DNA is pooled DNA of one or more individuals who show the phenotype of interest, and the unaffected DNA is pooled DNA of one or more individuals who do not show the phenotype of interest.
4. The method of claim 1, wherein the affected DNA is DNA of one individual who shows the phenotype of interest, and the unaffected DNA is pooled DNA of a

complete set of ancestors who do not show the phenotype of interest.

5. The method of claim 1, wherein the affected DNA is DNA from cells of an individual that show the phenotype of interest, and the unaffected DNA is DNA from cells of the individual that do not show the phenotype of interest.
6. (once amended) The method of claim 1, wherein step b) is performed by use of a mismatch-binding protein.
7. (once amended) The method of claim 1, wherein either the fragments of the affected DNA or the fragments of the unaffected DNA are tagged by one member of a specific binding pair, and step c) is performed by using the other member of the specific binding pair.
8. The method of claim 7, wherein the fragments of the unaffected DNA are tagged with biotin, and step c) is performed by use of immobilised streptavidin.
9. (once amended) The method of claim 1, further comprising subjecting the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest to self-hybridisation to form duplexes and subsequently recovering the perfectly matched duplexes.

10. (once amended) The method of claim 1, further comprising mixing the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest with an excess of the fragments of the affected DNA under hybridisation conditions to form duplexes and subsequently recovering the perfectly matched duplexes.
11. (once amended) The method of claim 1, wherein each of the affected DNA and the unaffected DNA is provided in fragmented form by digestion with from 4 to 7 six-cutter restriction endonuclease enzymes together with from 0 to 50 four-cutter restriction endonuclease enzymes.
14. A method of making a set of arrays of fragments of DNA of interest, which method comprises;
- a) selecting, from a set of n restriction endonuclease enzymes, a subset of r restriction endonuclease enzymes;
 - b) digesting genomic DNA with the subset of r enzymes;
 - c) ligating to the resulting fragments restriction-enzyme-cutting-site-specific adapters with unique polymerase chain reaction amplifiable sequences;
 - d) splitting the resulting fragments into r^2 aliquots;
 - e) amplifying each aliquot with two restriction-enzyme-specific primers;
 - f) forming an array of the r^2 aliquots of non-tagged amplicon strands; and

g) repeating steps a) to f) using one or more different subsets of r restriction endonuclease enzymes.

15. (once amended) The method of claim 14, wherein steps a to f) are repeated using each different subset of r restriction endonuclease enzymes to give $(n!)/[(n-r)!r!]$ different arrays.
16. (once amended) The method of claim 14, wherein the n restriction endonuclease enzymes are selected from 4-cutters and 5-cutters and 6-cutters.
17. (once amended) The method of claim 14, wherein the n is 3 to 10 and r is 2 to 4.
18. The method of claim 17, wherein $n = 6$ and $r = 3$.
21. (once amended) A set of arrays produced by the method of claim 14, derived from a set of $n = 6$ six-cutter restriction endonuclease enzymes which are *BamHI*; *BsrGI*; *Hind III*; *NcoI*; *SpeI*; and *AfIII*.
22. (once amended) A set of arrays produced by the method of claim 14, derived from the set of $n = 6$ six-cutter restriction endonuclease enzymes which are *EcoRI*; *BspHI*; *BglII*; *XbaI*; *Acc65I*; and *ApaLI*.

0914504-013103

23. (once amended) A nucleic acid characterisation method which comprises presenting to a set of arrays produced by the method of claim 14 a nucleic acid fragment of interest under hybridisation conditions, and observing a pattern of hybridisation.
24. The method of claim 23, wherein a plurality of nucleic acid fragments of interest are separately presented to the set of arrays, and the resulting patterns of hybridisation are compared.
25. The method of claim 24, wherein the plurality of nucleic acid fragments of interest are drawn from the mixture of DNA fragments, enriched in fragments that are characteristic of a phenotype of interest, of claim 13.
27. A double-stranded DNA molecule having the sequence a-A-b-B...X-y-Y-z where A, B...X and Y are unique restriction sites for n different restriction endonuclease enzymes, and a, b...y, z denotes distances in base pairs, characterised in that each fragment, obtainable by cutting the DNA molecule by means of any one or more up to n of the restriction enzymes, has a different length from every other fragment.
28. (once amended) The double-stranded DNA molecule of claim 27, wherein:

- a) inter-fragment length differences are greater for larger fragments;
- b) all possible fragments are unambiguously resolvable by electrophoresis from one another;
- c) size gaps between bands comprising different numbers of inter-restriction-site units are larger than size gaps between bands comprising the same number of inter-restriction-site units;
- d) the size gaps and size spread from the largest to the smallest fragment are electrophoretically compatible.

29. (new) The method of claim 1 which further comprises steps a), b) and c) one or more times.
30. (new) The method of claim 14 wherein one of the restriction enzyme specific primers is tagged.

5/PTS

GENETIC ANALYSIS

5 INTRODUCTION

Limitations of the current approaches

There are a number of limitations to carrying out association studies using single nucleotide polymorphisms (SNPs) and linkage disequilibrium within human populations (see Science, Vol 278, p1580, (1997) for a review of such methods). We have no control over recombination frequency around a given locus or over past human genetic crossing. Some mutations will be closely correlated with nearby SNPs and others will not.

15

The need for whole genome analysis

With the SNP and linkage disequilibrium approach (and many others), markers are essentially used as a surrogate for sequencing - the more markers, the better. The logical endpoint of the above argument is to look at every base in the human genome - and carry out what could be termed a whole genome association study. In essence, the sequence at every base would be determined for the genome of each member of a phenotypically 'affected' and a phenotypically 'unaffected' population. Statistical correlations (associations) could then be drawn between sequence differences and phenotype. Such associations would have future predictive values for the phenotype of interest, knowing the genotype and could be of great value in medicine and pharmacogenetics.

The current invention selectively enriches for DNA fragments that determine phenotype in the 'affected' population and thus makes the prospect of carrying out whole genome association studies for humans and other species a very real possibility.

09/914604-0100

Definition of terms used with the current invention

Within the scope of the current invention, the individuals chosen for whole genome analysis may be human, animal or plant and they may be eukaryotic, prokaryotic or archaebacterial in origin.

5 The terms 'affected' and 'unaffected' are used without limitation in order to categorise individuals into two groups – those that possess a defined phenotype of interest ('affected' individuals) and those that do not possess the phenotype of interest ('unaffected' individuals). The phenotype common to the 'affected' individuals may be either
10 beneficial (e.g. for these individuals, a particular pharmaceutical entity might show high efficacy in a phase II clinical trial) or detrimental (e.g. for these individuals, a particular pharmaceutical entity might show adverse toxicology in a phase I clinical trial).

 The 'affected' population may comprise one or more
15 individuals and the 'unaffected' population may similarly comprise one or more individuals according to the particular embodiment of the invention (see below).

 The term DNA is used throughout for simplicity. Within the scope of the current invention, the term DNA may equally well apply to all
20 or part of the haploid, diploid or polyploid genomic DNA content of one or more germ line or somatic cell(s). The DNA may be extracted from cells taken directly from the individual(s), the DNA may be extracted from cells cultured or immortalised from the individual(s) or the DNA may be prepared from a library of clones - with inserts derived from the individual(s) and
25 propagated in some appropriate host and cloning vector system. For the particular case wherein the term DNA refers to the expressed part of the haploid, diploid or polyploid genomic DNA content of one or more somatic cells and the DNA is prepared from a library of clones - with inserts derived from the individual(s) and propagated in some appropriate host and cloning
30 vector system, a cDNA library (normalised or otherwise) may be used.

In the current invention, DNA is compared in fragmented form. Fragmentation can be performed after DNA extraction, prior to cloning and/or after cloning. Restriction enzyme digestion is the preferred method for such fragmentation - though other methods (e.g. shearing or
5 sonication) will be obvious to those skilled in the art.

For the particular case wherein the DNA is prepared from a library of clones (either genomic clones or cDNA clones) - with inserts derived from the individual(s) and propagated in some appropriate host and cloning vector system and wherein restriction enzyme fragmentation is
10 used prior to cloning, polymerase chain reaction amplification can be used to prepare the DNA for comparison in fragmented form. Priming sites within the vector sequence flanking the cloned restriction enzyme fragmented inserts may be usefully employed for one or more cycles of polymerase chain reaction amplification of the fragmented DNA of interest.
15 The primers used for polymerase chain reaction amplification of the fragmented DNA of interest could again be used after the phenotype-determining fragment enrichment process to 'rescue' and clone the enriched fragments.

Within the scope of the current invention, the terms
20 biotinylation and streptavidin capture are used both as an example and as the currently preferred embodiment for the invention. The streptavidin may be surface attached to inert particles (magnetic or otherwise) or to vessel walls (e.g. microtitre plate wells). The biotin may be introduced via a deoxynucleotide triphosphate analogue using a polymerase; by using a
25 biotin-conjugated primer and polymerase chain reaction amplification; chemically or photochemically. The use of biotin and streptavidin is not a limitation for the invention. The invention could equally be used with other high affinity capture systems well known to those skilled in the art (e.g. 'his tag' introduction and metal ion affinity capture).

30 Within the scope of the current invention, the term 'abnormal' - used with respect to the term 'normal' - is used without limitation in order

09914604-013102

to denote a somatic cell (or somatic cells) with a discernable phenotypic characteristic (or characteristics) arising from the acquisition of a different somatic mutation (or set of somatic mutations) from that (or those) seen in the 'normal' counterpart. Cells will most usually be considered 'abnormal' with respect to their 'normal' counterparts through the acquisition of a different somatic mutation (or set of somatic mutations) leading to one or more of the following phenotypic characteristics: altered marker gene expression, altered genomic organisation, growth under certain selective culture conditions, immortalised growth in culture, unrestrained growth *in vivo* or *in vitro*, failure of normal apoptotic control mechanisms *in vivo* or *in vitro*, induction of neovascularisation, escape of cells across epithelium, migratory cell survival or metastasis.

Within the scope of the current invention, the term mismatch recognition protein is used without limitation to denote a protein of eukaryotic, prokaryotic or archaebacterial origin capable of the selective recognition of (and binding to) a DNA duplex that is not perfectly matched along its entire length. Recognition of (and binding to) will be preferably for bases that are not engaged in correct Watson and Crick pairing and for small deletions or insertions. Many such proteins are known to those skilled in the art. Prokaryotic and eukaryotic mutS homologues, phage T4 endonuclease VII, phage T7 endonuclease I and the plant enzyme CEL-1 are just some examples.

Inter-population perfectly matched duplex depletion

In the inter-population perfectly matched duplex depletion approach, we compare (in fragmented form) the pooled DNA of 'affected' individuals with the pooled DNA of 'unaffected' individuals (both from populations as outbred and otherwise similar to each other as possible). We are only interested in those regions where differences occur between 'affected' and 'unaffected' DNA molecules. For populations as above, the only prevalent sequence differences within the 'affected' pool (compared to

the 'unaffected' pool) should be somewhere within the gene(s) (using the term gene in its widest sense to include exons, introns and all associated upstream and downstream regulatory sequences) that actually determine(s) their common phenotype. This means that we are no longer
5 tied into working with rare (and perhaps atypical) populations where there is high genetic homogeneity.

Pooling the DNA from entire phenotypically-defined populations massively reduces the amount of labour involved.

10 DNA sequence variation in populations

Nickerson *et al*, Nature Genetics, Vol 19, p233 (1998), sequenced 9.7 kb of the lipoprotein lipase gene from 71 individuals (24 African-Americans, 24 Europeans and 23 European-Americans). This gene is fairly typical (90 % intron and 10 % exon - total size 30 kb with 10 exons).

15 88 sequence variants were found (i.e. one per 110 bp on average). Most variations were found in non-coding sequence. 90 % of these were SNPs (60 % of which were transitions and 40 % were transversions). All of the SNPs were biallelic. 10 % of the sequence variants were insertions or deletions at repeat sequences.

20 58 % of the sequence variants were present in all three ethnic populations. Half of these were found in heterozygous form and half in homozygous form.

Nucleotide diversity (defined as the expected number of nucleotide differences per site between a random pair of chromosomes
25 drawn from the population) is 1/500 for DNA in general and 1/2,000 for coding sequence DNA. This means that, on average, any two DNA fragments annealed together from such a population will contain a mismatch every 500 bases.

DNA sequence variants are therefore very common. They
30 are not, however, totally random - the variants that occur every 500 bases or so are limited; they are generally biallelic at just that single base. It is

- 6 -

this fact that the inter-population perfectly matched duplex depletion approach selectively exploits.

Statistical analysis of the inter-population perfectly matched duplex depletion fragmentation process

If the length of the DNA fragments is F bases and the average length between sequence differences between any two DNA molecules is 500 bases, the probability that a hybrid duplex between any two random DNA fragments will contain no mismatches is given by

$$\Pr(0) = e^{-(F/500)}$$

and the probability that a hybrid duplex between any two random DNA molecules will contain any mismatches is given by

$$\Pr(\geq 1) = \{1 - (e^{-(F/500)})\}$$

Example values for $\Pr(0)$ and $\Pr(\geq 1)$ for different values of F are given below

F	$\Pr(0)$	$\Pr(\geq 1)$
10	0.98	0.02
20	0.96	0.04
50	0.90	0.10
100	0.82	0.18
200	0.67	0.33
300	0.55	0.45
400	0.45	0.55
500	0.37	0.63
600	0.30	0.70
700	0.25	0.75
800	0.20	0.80
900	0.17	0.83
1000	0.14	0.86

Example average restriction fragment sizes for DNA digestion with six 6 bp cutters and up to four 4 bp cutters

Each 6 bp cutter will cut DNA every 4.096 bp on average and each 4 bp cutter will cut DNA every 256 bp on average.

- 5 For a given set of A 6 bp cutters and B 4 bp cutters (with no duplication of restriction enzyme cutting sequence), the average fragment length F will be

$$\{(A/4.096)+(B/256)\}^{-1}$$

10

Example values for F as A and B are varied are given in the following table

A	B	average size (F)
6	0	683
6	1	186
6	2	108
6	3	76
6	4	59

- 15 We should note that the above situation assumes that none of the 4 bp cutter recognition sites lie within any of the 6 bp cutter recognition sites. If, for example, we have a 4 bp cutter recognition site nested within a 6 bp cutter recognition site (e.g. from the use of *MboI* and *BamHI* in the fragmentation), then we should reduce the value of A from 6 to 5 in the
- 20 above.

In general, if we have a given set of A enzymes that cut DNA every a bp on average, B enzymes that cut DNA every b bp on average, ... , Z enzymes that cut DNA every z bp on average (with no duplication of restriction enzyme cutting sequence), the average fragment length F will be

25

$$\{(A/a)+(B/b)+ \dots + (Z/z)\}^{-1}$$

207670-4094T660

Example sets of 6 bp cutters and 4 bp cutters for fragmentation

Example sets of 6 bp cutters and 4 bp cutters that contain panels of six 6 bp cutters that are compatible with terminal restriction site profiling array (TRSPA) analysis (see below) are given in the following:

5

Example enzyme set 1

For restriction in M buffer + BSA

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>Mbol</i>	GATC	100	37	25	NEB	65/20
2	<i>HaeIII</i>	GGCC	100	37	50	NEB	80/20
3	<i>MseI</i>	TTAA	100	37	20	NEB	65/20

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>BamHI</i>	GGATCC	100	37	100	NEB	80 / 20
2	<i>BsrGI</i>	TGTACA	100	37	10	NEB	80 / 20
3	<i>HindIII</i>	AAGCTT	100	37	>40	APB	65 / 20
4	<i>NcoI</i>	CCATGG	100	37	50	NEB	65 / 20
5	<i>SpeI</i>	ACTAGT	100	37	50	NEB	65 / 20
6	<i>AflIII</i>	CTTAAG	100	37	10	APB	60 / 15

10

Example enzyme set 2

For restriction in M buffer + BSA

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>Mbol</i>	GATC	100	37	25	NEB	65/20
2	<i>HaeIII</i>	GGCC	100	37	50	NEB	80/20
3	<i>MseI</i>	TTAA	100	37	20	NEB	65/20

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>EcoRI</i>	GAATTC	100	37	>40	APB	65 / 20
2	<i>BspHI</i>	TCATGA	100	37	10	APB	65 / 20
3	<i>BglII</i>	AGATCT	75	37	>40	APB	No
4	<i>XbaI</i>	TCTAGA	100	37	100	NEB	65 / 20
5	<i>Acc65I</i>	GGTACC	75	37	10	NEB	65 / 20
6	<i>ApaLI</i>	GTGCAC	100	37	10	NEB	No

Aspect (I)

In one aspect the invention provides a method of providing a mixture of DNA fragments enriched in fragments that are characteristic of a phenotype of interest, by providing affected DNA in fragmented form and providing unaffected DNA in fragmented form, which method comprises:

- a) mixing the fragments of the affected DNA and the fragments of the unaffected DNA under hybridising conditions;
 - b) recovering a mixture of hybrids that contain mismatches;
 - 10 c) recovering fragments of the affected DNA from the mixture of hybrids that contain mismatches;
- and optionally repeating steps a), b) and c) one or more times.

15 Inter-population mismatch-containing duplex selection

'Affected' versus 'unaffected' (i.e. inter-population) mismatch-containing duplex selection can be achieved by attaching a mismatch-binding protein to a solid support (or using the mismatch-binding protein in solution followed by subsequent solid-phase capture), taking fragmented and denatured 'affected' DNA and hybridising this to an excess of
20 fragmented, denatured and biotinylated 'unaffected' DNA with ensuing capture of mismatch-containing duplex molecules. Releasing the mismatch-containing duplex molecules without denaturation, streptavidin


2013-10-10 10:10:10

- 10 -


capture and then release of the non-biotinylated strands will give only the desired species as shown below.

Method

- 5 Fragment the 'affected' DNA. Fragment and derivatise the 'unaffected' DNA with biotin. Only DNA from this population will be streptavidin captured. Melt and anneal to give

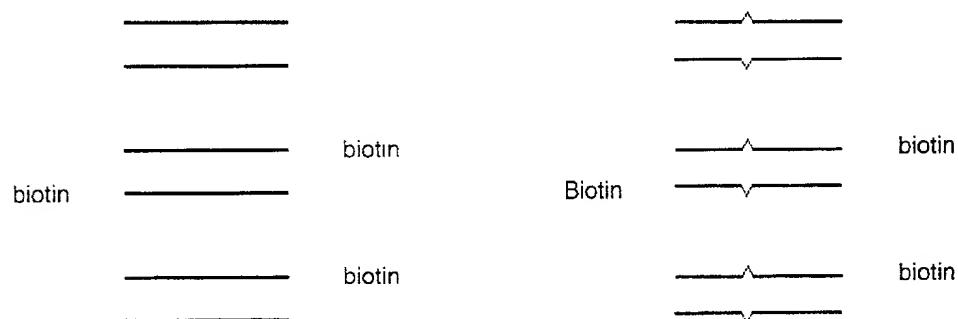


- 10 Mismatch-binding protein select. Capture only the mismatch-containing duplexes. Release without denaturation to give

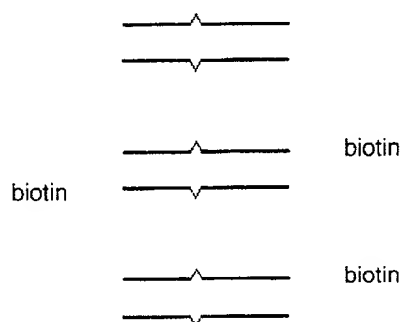


- 11 -

Streptavidin capture to give

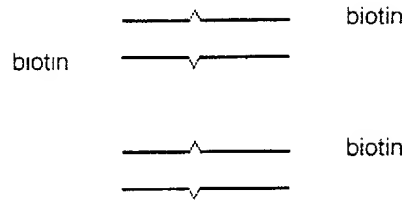

Release the non-biotinylated strands to give

5

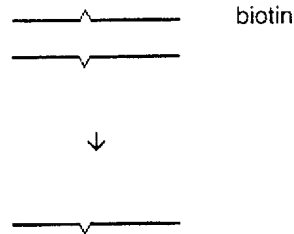

Repeat as necessary.

Repetition of the above sequence of reactions will lead to inter-population perfectly matched duplex depletion.

10

What will be purified ?

Inter-population mismatch-containing duplex selection as above ensures that all of the various phenotype-determining fragments (unique to the 'affected' population) are captured for subsequent analysis - but it also causes the co-purification of very many SNP-containing ('noise') fragments.

We now need to consider the fate of the various types of fragment (i.e. those that determine the phenotype and those that do not) as we carry out inter-population perfectly matched duplex depletion cycles as above.

20

0944634.013102

Recovery of 'affected' DNA molecules after streptavidin capture

For a particular fragment, if we have X molecules of 'affected' DNA and Y molecules of 'unaffected' biotinylated DNA, after complete hybridisation, there will be a ratio of $\{Y/(X+Y)\}$ streptavidin-capturable molecules to $\{X/(X+Y)\}$ streptavidin-non-capturable molecules. We can thus manipulate the yield of streptavidin-capturable hybrids by varying X and Y.

Example recovery and loss figures for various X and Y are shown in the following table

ratio (Y/X)	recovery = $\{Y/(X+Y)\}$	Loss = $\{X/(X+Y)\}$
1	50 %	50 %
2	67 %	33 %
3	75 %	25 %
4	80 %	20 %
9	90 %	10 %
19	95 %	5 %
99	99 %	1 %
999	99.9 %	0.1 %

After n cycles, the recovery for a phenotype-determining fragment will be given by $\{Y/(X+Y)\}^n$

Loss of general SNP-containing ('noise') fragments during inter-population perfectly matched duplex depletion cycles

If we anneal fragmented DNA molecules together and capture only the mismatch-containing duplexes, then n repetitions of such a process will reduce the original number of fragments to the following fraction

$$\{1 - (e^{-F/500})\}^n \cdot \{Y/(X+Y)\}^n$$

The enrichment for phenotype-determining fragments over SNP-containing ('noise') fragments during inter-population perfectly matched duplex depletion cycles will therefore be given by

$$\{1 - (e^{-(F/500)})\}^n$$

Example figures for enrichment are given below with $F = 50, 100, 200, 300, 400, 500, 600, 700, 800, 900$ and $1,000$ and $n = 1, 2, 3, 4, 5, 6, 7, 8, 9$ and 10 .

F	n =									
	1	2	3	4	5	6	7	8	9	10
50	11	110	1,160	12,194	128,135	1,346,489	1.E+07	1.E+08	2.E+09	2.E+10
100	6	30	168	926	5,110	28,187	155,500	857,841	4,732,415	3.E+07
200	3	9	28	85	257	779	2,362	7,166	21,735	65,929
300	2	5	11	24	53	119	263	582	1,291	2,860
400	2	3	6	11	20	36	65	118	215	390
500	2	3	4	6	10	16	25	39	62	98
600	1	2	3	4	6	9	12	18	25	36
700	1	2	2	3	4	5	7	10	13	17
800	1	2	2	2	3	4	5	6	8	10
900	1	1	2	2	2	3	4	4	5	6
1,000	1	1	2	2	2	2	3	3	4	4

Loss of specific SNP-containing ('noise') fragments during inter-population perfectly matched duplex depletion cycles

The non-polymorphic and SNP-containing ('noise') fragments will be depleted as described above.

Not all fragments will, however, be depleted at the same rate. An individual SNP-containing ('noise') fragment will be depleted with every cycle of inter-population perfectly matched duplex depletion as outlined below.

Let us assume that there are two alleles for the polymorphism within a particular fragment. Let these be called P and Q. Let p be the

- 14 -

fraction of the P allele in the (outbred) 'affected' and 'unaffected' populations and let q be the fraction of the Q allele. Let $(p+q)=1$, so that $q=(1-p)$.

After denaturation and annealing, the four possible events are
 5 PP, PQ, QP and QQ. PP and QQ will form perfectly matched duplexes and will therefore be lost - whereas PQ and QP will form mismatch-containing duplexes and will consequently be recovered.

After one cycle, the fraction of recovered molecules will be

$$\begin{aligned}
 & \{2pq/(p^2+2pq+q^2)\} \\
 & = \{2pq/((p+q)^2)\} \\
 & = 2pq \\
 & = 2p(1-p)
 \end{aligned}$$

Hence if we start out with M molecules of DNA from the population, there will be 2Mp q molecules remaining after the first round of
 20 inter-population mismatch-containing duplex selection. Let us denote the number of molecules entering the second round of inter-population mismatch-containing duplex selection by M', where

$$M' = 2Mp q$$

25 and the fraction of lost molecules will be

$$\begin{aligned}
 & \{(p^2+q^2)/(p^2+2pq+q^2)\} \\
 & = \{(p^2+q^2)/((p+q)^2)\}
 \end{aligned}$$

200914604-01310

- 15 -

$$= p^2 + q^2$$

$$= p^2 + \{1 - 2p + p^2\}$$

$$= 1 - 2p + 2p^2$$

$$= 1 - 2p(1 - p)$$

$$= 1 - 2pq$$

10

The fractional loss of P-allelic molecules will be p^2 and the fractional loss of Q-allelic molecules will be q^2 ($= (1-p)^2$).

If we start out with M molecules of DNA from the population, there will be Mp of the P-allelic molecules and $Mq = M(1-p)$ of the Q-allelic molecules before inter-population mismatch-containing duplex selection. After inter-population mismatch-containing duplex selection, there will therefore be $Mp - Mp^2 = Mp(1-p) = Mqp$ of the P-allelic molecules and $Mq - Mq^2 = Mq(1-q) = Mqp$ of the Q-allelic molecules. In other words, after the first round of mismatch-containing duplex selection, there will be the same number of P-allelic molecules as Q-allelic molecules.

20

We can define new allelic frequencies p' and q' as follows

$$p' = q' = 0.5$$

25

If we now perform a second round of inter-population mismatch-containing duplex selection, we start out with M' molecules of DNA from the first round. There will be $M'p'$ of the P-allelic molecules and $M'q'$ of the Q-allelic molecules before inter-population mismatch-containing duplex selection.

30

After inter-population mismatch-containing duplex selection, there will therefore be $M'p' - M'p'^2 = M'p'(1-p) = M'p'q$ of the P-allelic

094464-01318

- 16 -

molecules and $M'q' - Mqq' = M'q'(1-q) = M'q'p$ of the Q-allelic molecules.

The total number of molecules (M'') will be

$$M'p'q + M'q'p$$

$$= M'(p'q + q'p)$$

$$\text{but } p' = q' = 0.5$$

$$\text{hence } M'' = 0.5 * M'(p+q)$$

$$\text{but } (p+q) = 1$$

$$\text{so } M'' = (M'/2)$$

In other words, after the second round of inter-population mismatch-containing duplex selection, there will again be the same number of P-allelic molecules as Q-allelic molecules. Thus the new allelic frequencies p'' and q'' remain as previously $p'' = q'' = 0.5$.

A pattern now emerges. After the first cycle of inter-population mismatch-containing duplex selection, the allelic frequencies are equalled and the number of molecules is reduced to $2Mp q$ ($Mp q$ of each allelic molecule). Thereafter, every cycle halves the number of molecules and keeps both alleles at the same frequency.

Consequently, after n cycles, the number of P-allelic molecules will be

$$Mp q (0.5)^{n-1}$$

and the number of Q allelic molecules will be

00914604-013102

- 17 -

$$Mpq(0.5)^{n-1}$$

both are, of course, equal.

If we now take the capture yield (see above) into
 5 consideration, the SNP-containing ('noise') fragment yield will be given by

$$2Mpq(0.5)^{n-1} \cdot \{Y/(X+Y)\}^n$$

where both allelic variants are deemed to be captured 'noise'.

10

Polymorphisms that interfere with the pattern of restriction digestion

For both the loss of general and specific SNP-containing
 ('noise') fragments during inter-population perfectly matched duplex
 depletion cycles (described above) and where the SNP interferes with the
 15 pattern of restriction digestion, if the mismatch-binding protein also binds to
 duplex molecules with unequal lengths (e.g. from inter-population
 annealing around a site of restriction site polymorphism), then the above
 analysis still holds (with perfectly matched duplex being replaced by equal
 length duplex and mismatch-containing duplex being replaced by unequal
 20 partner-length duplex).

In the rare cases where a restriction site is lost due to a
 sequence change that actually determines the phenotype of interest, a
 double-length fragment will be obtained. This will give rise to a double
 terminal restriction site profiling array (TRSPA) signature (see below).
 25 Multiple isolates of the particular double signature will be indicative of an
 association between the fragment and the phenotype of interest.

00914604-013103

Further 'kinetic' enrichment to enhance the selective removal of SNP-containing 'noise' from the pool of phenotype determining fragments

After multiple cycles of enrichment by the above procedure, the enriched DNA pool should contain many copies of all phenotype-determining fragments but also low numbers of copies of many different phenotype non-determining fragments. The total number of 'noise' fragments may exceed the number of phenotype determining fragments, despite the number of each individual 'noise' species being very small. The 'noise' fragments would therefore increase the number of probes required for TRSPA analysis before a pattern emerges. To largely eliminate this problem, a further kinetic enrichment procedure is used. Either one or both of strategies A and B below can be employed to achieve 'kinetic' enrichment.

Strategy A – Subtraction of the enriched DNA from inter-population mismatch containing duplex depletion

The enriched fragment pool from inter-population mismatch containing duplex depletion is rapidly self-hybridised - enabling the common phenotype-determining fragments to form perfectly matched duplexes with greater efficiency than the rare 'noise' fragments. Selection for perfectly matched duplexes then yields a selectively further enriched pool of fragments. Multiple cycles of subtraction could be carried out if necessary.

Strategy B – Hybridisation of the enriched DNA from inter-population mismatch containing duplex depletion against the 'affected' DNA pool

The enriched fragment pool from inter-population mismatch containing duplex depletion is then hybridised to an excess of biotinylated DNA from the 'affected' pool. This allows the common phenotype-determining fragments to form perfectly matched duplexes with greater efficiency than the rare 'noise' fragments. Selection for perfectly matched

duplexes followed by streptavidin capture and denaturation to release the non-biotinylated strands then yields a further enriched pool of fragments. Multiple such 'affected' pool back-hybridisations could be carried out if necessary.

5

Extension of the invention to the case of single phenotypically 'affected' individuals within populations where the distinction between 'affected' and 'unaffected' is clear

The above has described inter-population perfectly matched duplex depletion between non-biotinylated DNA fragments from an 'affected' population and biotinylated DNA fragments from an 'unaffected' population. Provided the 'unaffected' population is sufficiently complex that it contains all the non-phenotype-determining sequence variants found in a single 'affected' individual, then inter-population perfectly matched duplex depletion should be possible for single phenotypically 'affected' individuals against an 'unaffected' population where the distinction between 'affected' and 'unaffected' is clear. The latter proviso is needed in order to ensure that a small number of misdiagnosed 'affected' individuals in the 'unaffected' population do not cause the removal of phenotype-determining fragments during inter-population perfectly matched duplex molecular depletion.

Extension of the invention to the case of disease gene identification in cases where novel phenotype-determining mutations arise spontaneously within a family

Except for a small number of sequence changes, each of us contains DNA sequence derived from our parents - our individuality resulting from precisely which parental alleles we receive. If one of the above small number of sequence changes results in a change in phenotype, then we can use inter-population perfectly matched duplex depletion to enrich for fragments encoding this change in phenotype.

If we take 'unaffected total ancestral' cells (by which we mean cells derived from a complete set of 'unaffected' ancestors - e.g. both parents, or mother plus two paternal grandparents, or father plus two maternal grandparents, or two maternal grandparents and two paternal grandparents etc.) as the source of our biotinylated fragments and cells from an 'affected' descendent as the source of our non-biotinylated fragments, any fragments that have acquired phenotype-determining sequence changes between 'unaffected' ancestral generations and the 'affected' descendent generation will be unable to form perfectly matched duplexes with the biotinylated 'unaffected total ancestral' fragments. Successive cycles of such inter-population perfectly matched duplex depletion will thus lead to the enrichment of fragments carrying all such sequence - the degree of enrichment per cycle being as described below.

15 Statistical analysis

Let us assume that equal numbers of fragments are used from each of the 'unaffected' ancestors. Let the number of such ancestors be A.

$\{1/A\}$ of the annealings in the inter-population perfectly matched duplex depletion will be self-against-ancestral-transmitted alleles - statistically equivalent to self-against-self inter-population perfectly matched duplex depletion (see below).

$\{(A-1)/A\}$ of the annealings in the inter-population perfectly matched duplex depletion will be self-against-ancestral-nontransmitted alleles - statistically equivalent to inter-population perfectly matched duplex depletion between unrelated individuals.

A	$\{1/A\}$	$\{(A-1)/A\}$
2	$1/2$	$1/2$
3	$1/3$	$2/3$
4	$1/4$	$3/4$

From the data of Nickerson *et al.* a DNA sequence variation between unrelated individuals should occur every 500 base pairs (it is this figure that we use for inter-population perfectly matched duplex depletion with self-against-ancestral-nontransmitted alleles and inter-population perfectly matched duplex depletion between unrelated individuals). In addition, a given individual should be heterozygous once every 573 base pairs (this figure is used for inter-population perfectly matched duplex depletion with self-against-ancestral-transmitted alleles and self-against-self inter-population perfectly matched duplex molecular depletion). Inter-population perfectly matched duplex depletion against the transmitted alleles and the nontransmitted alleles will now be considered separately.

Self-against-ancestral-transmitted alleles inter-population perfectly matched duplex depletion

If we anneal the two complementary strands for a fragment using DNA from an 'affected' descendent and DNA containing the transmitted alleles of 'unaffected' ancestors, then $\{1 - e^{(-F/573)}\}$ of fragments will contain one or more site of heterozygosity. $\{1/A\}$ of the annealings will be of this type. For such annealings, where a site of heterozygosity is present, the probability of obtaining a mismatch-containing duplex between a biotinylated fragment and a non-biotinylated fragment containing the site of heterozygosity is 0.5.

Self-against-ancestral-nontransmitted alleles inter-population perfectly matched duplex depletion

If we anneal the two complementary strands for a fragment using DNA from an 'affected' descendent and DNA containing the nontransmitted alleles of 'unaffected' ancestors, then $\{1 - e^{(-F/500)}\}$ of fragments will contain one or more site of DNA sequence variation. $\{(A-1)/A\}$ annealings will be of this type.

Phenotype-determining fragment enrichment

The fraction of fragments carried through the first cycle of inter-population perfectly matched duplex depletion will therefore be

$$5 \quad ([0.5 \cdot \{1/A\} \cdot \{1 - e^{(-F/573)}\}] + [\{(A-1)/A\} \cdot \{1 - e^{(-F/500)}\}]) \cdot \{Y/(X+Y)\}$$

where $\{Y/(X+Y)\}$ represents the yield of streptavidin-capturable molecules.

Hence n repetitions of such a process will reduce the original
10 number of fragments to the following fraction

$$([0.5 \cdot \{1/A\} \cdot \{1 - e^{(-F/573)}\}] + [\{(A-1)/A\} \cdot \{1 - e^{(-F/500)}\}])^n \cdot \{Y/(X+Y)\}^n$$

The enrichment for phenotype-determining fragments over
15 SNP-containing ('noise') fragments during inter-population perfectly matched duplex depletion cycles will therefore be given by

$$([0.5 \cdot \{1/A\} \cdot \{1 - e^{(-F/573)}\}] + [\{(A-1)/A\} \cdot \{1 - e^{(-F/500)}\}])^{-n}$$

20 Example figures for enrichment are given for $A = 2, 3$ and 4 below with $F = 50, 100, 200, 300, 400, 500, 600, 700, 800, 900$ and $1,000$ and $n = 1, 2, 3, 4, 5, 6, 7, 8, 9$ and 10 .

- 23 -

For A=2

F	n=									
	1	2	3	4	5	6	7	8	9	10
50	15	213	3,115	45,494	664,427	9,703,694	1.E+08	2.E+09	3.E+10	4.E+11
100	8	59	448	3,430	26,250	200,887	1,537,367	1.E+07	9.E+07	7.E+08
200	4	18	74	309	1,296	5,434	22,782	95,524	400,520	1,679,340
300	3	9	28	87	265	811	2,475	7,558	23,078	70,468
400	2	6	16	39	97	241	600	1,497	3,734	9,313
500	2	5	10	22	48	103	224	485	1,051	2,277
600	2	4	7	15	29	56	109	213	416	813
700	2	3	6	11	19	35	63	114	207	374
800	2	3	5	8	14	24	42	71	121	205
900	2	3	4	7	11	18	30	49	79	128
1,000	2	2	4	6	9	15	23	36	57	89

For A=3

5

F	n=									
	1	2	3	4	5	6	7	8	9	10
50	13	167	2,154	27,823	359,340	4,640,941	6.E+07	8.E+08	1.E+10	1.E+11
100	7	46	311	2,104	14,251	96,516	653,682	4,427,241	3.E+07	2.E+08
200	4	14	51	191	709	2,633	9,786	36,365	135,135	502,173
300	3	7	20	54	146	396	1,073	2,909	7,882	21,361
400	2	5	11	24	53	119	263	582	1,290	2,859
500	2	4	7	14	27	51	99	190	367	707
600	2	3	5	9	16	28	48	84	147	255
700	2	3	4	7	11	18	28	46	74	119
800	2	2	4	5	8	12	19	28	43	66
900	1	2	3	4	6	9	14	20	29	41
1,000	1	2	3	4	5	8	11	15	21	29

201604101650

For A=4

	n=									
F	1	2	3	4	5	6	7	8	9	10
50	12	149	1,826	22,316	272,759	3,333,771	4.E+07	5.E+08	6.E+09	7.E+10
100	6	41	264	1,690	10,834	69,465	445,373	2,855,512	2.E+07	1.E+08
200	4	12	44	154	540	1,902	6,696	23,570	82,967	292,045
300	3	7	17	44	112	287	737	1,894	4,864	12,494
400	2	4	9	20	41	86	181	381	800	1,681
500	2	3	6	11	20	37	68	125	229	418
600	2	3	5	7	12	20	34	56	92	152
700	2	2	4	5	8	13	20	30	46	71
800	1	2	3	4	6	9	13	19	27	39
900	1	2	3	4	5	7	9	13	18	25
1,000	1	2	2	3	4	6	7	10	13	17

Extension of the invention to the case of fully comprehensive

5 'abnormal' cell mutational profiling within an individual

If we now take 'normal' cells as the source of our biotinylated fragments and 'abnormal' cells from the same individual as the source of our non-biotinylated fragments, any fragments that have acquired sequence changes on the way to becoming 'abnormal' will be unable to
 10 form perfectly matched duplexes with the biotinylated 'normal' fragments. Successive cycles of such inter-population perfectly matched duplex depletion will thus lead to the enrichment of fragments carrying all those sequence differences between the 'normal' cells and the 'abnormal' cells - the degree of enrichment per cycle being as described below.

15 From the data of Nickerson *et al*, a given individual should be heterozygous about once every 573 base pairs.

If we anneal the two complementary strands for a fragment using DNA from 'abnormal' cells and DNA from 'normal' cells, then
 $\{1 - e^{(-F/573)}\}$ of fragments will contain one or more site of heterozygosity.

20 For a heterozygous site, p and q are both 0.5. The probability of obtaining a perfectly matched duplex between a biotinylated fragment

- 25 -

and a non-biotinylated fragment containing the site is 0.5. Similarly, the probability of obtaining a mismatch-containing duplex between a biotinylated fragment and a non-biotinylated fragment containing the site is also 0.5.

5 The fraction of fragments carried through the first cycle of inter-population perfectly matched duplex depletion will therefore be

$$0.5 \cdot \{1 - e^{(-F/573)}\} \cdot \{Y/(X+Y)\}$$

10 hence n repetitions of such a process will reduce the original number of fragments to the following fraction

$$[0.5 \cdot \{1 - e^{(-F/573)}\}]^n \cdot \{Y/(X+Y)\}^n$$

15 The enrichment for phenotype-determining fragments over SNP-containing ('noise') fragments during inter-population perfectly matched duplex depletion cycles will therefore be given by

$$[0.5 \cdot \{1 - e^{(-F/573)}\}]^n$$

20

Example figures for enrichment are given below with $F = 50$, 100, 200, 300, 400, 500, 600, 700, 800, 900 and 1,000 and $n = 1, 2, 3, 4, 5, 6, 7, 8, 9$ and 10.

09944604.043403
2012-10-19 14:55:00

F	n=									
	1	2	3	4	5	6	7	8	9	10
50	24	573	13.711	328.171	7.854.627	2.E+08	4.E+09	1.E+11	3.E+12	6.E+13
100	12	156	1.948	24.329	303.844	3.794.731	5.E+07	6.E+08	7.E+09	9.E+10
200	7	46	313	2.123	14.412	97.830	664.076	4.507.784	3.E+07	2.E+08
300	5	24	118	580	2.845	13.958	68.490	336.072	1.649.059	8.091.708
400	4	16	63	251	999	3.977	15.831	63.012	250.815	998.346
500	3	12	41	139	479	1.644	5.650	19.411	66.688	229.114
600	3	9	29	90	278	856	2.638	8.128	25.047	77.179
700	3	8	23	65	183	520	1.475	4.183	11.862	33.639
800	3	7	19	50	133	353	937	2.491	6.621	17.599
900	3	6	16	41	103	259	654	1.652	4.171	10.532
1,000	2	6	14	34	84	202	490	1.188	2.880	6.978

This approach will enrich for fragments containing *all* sequence differences within the 'abnormal' cells - no prior knowledge of the genes (oncogenes, tumour suppressor genes etc.) that may need to be investigated is required.

Having thus isolated a fragment (or fragments) determining differences between the 'affected' and 'unaffected' populations, we can now proceed to their analysis.

Aspect (II)

In another aspect the invention provides a method of making a set of arrays of fragments of DNA of interest, which method comprises:

- a) selecting, from a set of n restriction endonuclease enzymes, a subset of r restriction endonuclease enzymes;
- b) digesting genomic DNA with the subset of r enzymes;
- c) ligating to the resulting fragments restriction-enzyme-cutting-site-specific adapters with unique polymerase chain reaction amplifiable sequences;
- d) splitting the resulting fragments into r^2 aliquots;
- e) amplifying each aliquot with two-restriction enzyme-specific

- 27 -

primers;

f) forming an array of the r^2 aliquots of the amplimers;

and

g) repeating steps a) to f) using a different subset of r restriction

5 endonuclease enzymes. The invention also includes sets of arrays obtained or obtainable by the method.

The n restriction endonuclease enzymes may be selected from 4-cutters and 5-cutters and 6-cutters, and a set may include enzymes from one or two or three of these categories. The value of n is preferably 3
10 to 10, for reasons discussed below. The value of r is less than n and is preferably 2 to 4, chosen with reference to the frequency with which the chosen enzymes cut nucleic acids, and ease of fragment amplification by PCR.

15 **Terminal restriction site profiling arrays (TRSPAs)**

If we use a total of n 6 bp cutter restriction enzymes within the total set of enzymes used for fragmentation, let us use subsets of r 6 bp cutter enzymes (taken from the total set n) to make (rxr) TRSPA test matrices as follows.

20

Version 1 - (TRSPA-1)

For each (rxr) TRSPA test matrix, digest DNA to completion with all r restriction enzymes. Ligate on restriction enzyme cutting-site-specific adaptors with unique polymerase chain reaction amplifiable tags.

25 Split into r^2 aliquots and for each aliquot, amplify with biotinylated restriction enzyme _{j} and non-biotinylated restriction enzyme _{k} tag-specific primers and array the non-biotinylated strands for all values of j and k between 1 and r .

0914604.03103

Example

Consider the following dsDNA

--[A]--1--[B]--2--[B]--3--[C]--4--[A]--5--[B]-- +
5 --[A]--1--[B]--2--[B]--3--[C]--4--[A]--5--[B]-- -

where the restriction enzyme cutting sites are denoted A, B
and C and the fragments after restriction digestion are denoted 1, 2, 3, 4
and 5. + and - denote the sense of the strands.

10 Cut to completion with A, B and C to give

09944604-049403
207670-409460

- 29 -

A--1--B B--2--B B--3--C C--4--A A--5--B +

A--1--B B--2--B B--3--C C--4--A A--5--B -

Polymerase chain reaction amplify according to the following
 5 primer matrix and streptavidin capture to give

	A	B	C
bio-A		bio-A--1--B + A--1--B - bio-A--5--B + A--5--B -	C--4--A + C--4--A-bio -
bio-B	A--1--B + A--1--B-bio - A--5--B + A--5--B-bio -	bio-B--2--B + B--2--B - B--2--B + B--2--B-bio -	bio-B--3--C + B--3--C -
bio-C	bio-C--4--A + C--4--A -	B--3--C + B--3--C-bio -	

Keep only the non-biotinylated strands to give

	A	B	C
bio-A		A--1--B - A--5--B -	C--4--A +
bio-B	A--1--B + A--5--B +	B--2--B - B--2--B +	B--3--C -
bio-C	C--4--A -	B--3--C +	

- 30 -

In simplified form

	A	B	C
bio-A		1-, 5-	4+
bio-B	1+, 5+	2-, 2+	3-
bio-C	4-	3+	

Repeat for all nC_r combinations of restriction enzymes to
 5 generate the TRSPA.

TRSPA-1 hybridisation patterns

The two types of TRSPA-1 hybridisation pattern we should
 expect using a probe resulting from inter-population perfectly matched
 10 duplex depletion are

1. Hybridisation to an off-diagonal element (e.g. row x, column y
 - where x and y are different) and its complementary element reflected
 across the diagonal (i.e. row y, column x) and
2. Hybridisation to an on-diagonal element (e.g. the element at
 15 row z, column z).

TRSPA-1 analysis - a worked example for n=3 and r=2

Let us take the dsDNA

20 -C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

where

A, B and C denote restriction enzyme cutting sites. 1, 2, 3, 4,
 5, 6, 7, 8 and 9 denote the restriction fragments after digestion.

25

The TRSPA-1 test matrices

There are ${}^3C_2=3$ TRSPA-1 test matrices - AB, BC and AC.

- 31 -

The test matrix hybridisation patterns

There are $\{(2 \cdot (2+1))/2\} = 3$ possible hybridisation patterns for each of the test matrices



5

The combinatorial diversity

If there are three possible TRSPA-1 test matrix hybridisation patterns and three possible test matrices, then there will be $3^3 = 27$ possible TRSPA-1 signatures.

10

Fragment 4 analysis - the AB TRSPA-1 matrix

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

15

Cut with A and B to give

-C-1-B

B-2-B

20 B-3-A

A-4-A

A-5-C-6-C-7-B

B-8-A

A-9-A-

25

The fragment complementary to fragment 4 will be

A-4-A

- 32 -

Polymerase chain reaction amplify and hybridise with just
fragment 4 to give

AB matrix	A	B
bio-A		
bio-B		

5 Fragment 4 analysis - the BC TRSPA-1 matrix

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

10 Cut with B and C to give

-C-1-B

B-2-B

B-3-A-4-A-5-C

15 C-6-C

C-7-B

B-8-A-9-A-

The fragment complementary to fragment 4 will be

20

B-3-A-4-A-5-C

Polymerase chain reaction amplify and hybridise with just
fragment 4 to give

25

094404-0440

- 33 -

BC matrix	B	C
bio-B		
bio-C		

Fragment 4 analysis - the AC TRSPA-1 matrix

Take the dsDNA

5

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

10 -C

C-1-B-2-B-3-A

A-4-A

A-5-C

C-6-C

15 C-7-B-8-A

A-9-A

A-

The fragment complementary to fragment 4 will be

20

A-4-A

Polymerase chain reaction amplify and hybridise with just
fragment 4 to give

AC matrix	A	C
bio-A		
bio-C		

25

207670-40947650

- 34 -

Fragment 5 analysis - the AB TRSPA-1 matrix

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

5

Cut with A and B to give

-C-1-B

B-2-B

10 B-3-A

A-4-A

A-5-C-6-C-7-B

B-8-A

A-9-A-

15

The fragment complementary to fragment 5 will be

A-5-C-6-C-7-B

20

Polymerase chain reaction amplify and hybridise with just
fragment 5 to give

AB matrix	A	B
bio-A		
bio-B		

Fragment 5 analysis - the BC TRSPA-1 matrix

25

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

- 35 -

Cut with B and C to give

-C-1-B

B-2-B

5 B-3-A-4-A-5-C

C-6-C

C-7-B

B-8-A-9-A-

10 The fragment complementary to fragment 5 will be

B-3-A-4-A-5-C

Polymerase chain reaction amplify and hybridise with just

15 fragment 5 to give

BC matrix	B	C
bio-B		
bio-C		

Fragment 5 analysis - the AC TRSPA-1 matrix

Take the dsDNA

20

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

25 -C

C-1-B-2-B-3-A

A-4-A

A-5-C

0944604-039402

- 36 -

C-6-C

C-7-B-8-A

A-9-A

A-

5

The fragment complementary to fragment 5 will be

A-5-C

10

Polymerase chain reaction amplify and hybridise with just
fragment 5 to give

AC matrix	A	C
bio-A		
bio-C		

Fragment 6 analysis - the AB TRSPA-1 matrix

15

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and B to give

20

-C-1-B

B-2-B

B-3-A

A-4-A

25

A-5-C-6-C-7-B

B-8-A

A-9-A-

- 37 -

The fragment complementary to fragment 6 will be

A-5-C-6-C-7-B

- 5 Polymerase chain reaction amplify and hybridise with just
fragment 6 to give

AB matrix	A	B
bio-A		
bio-B		

Fragment 6 analysis - the BC TRSPA-1 matrix

- 10 Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with B and C to give

15

-C-1-B

B-2-B

B-3-A-4-A-5-C

C-6-C

20

C-7-B

B-8-A-9-A-

The fragment complementary to fragment 6 will be

- 25 C-6-C

Polymerase chain reaction amplify and hybridise with just
fragment 6 to give

- 38 -

BC matrix	B	C
bio-B		
bio-C		

Fragment 6 analysis - the AC TRSPA-1 matrix

Take the dsDNA

5

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

10

-C

C-1-B-2-B-3-A

A-4-A

A-5-C

C-6-C

15

C-7-B-8-A

A-9-A

A-

The fragment complementary to fragment 6 will be

20

C-6-C

Polymerase chain reaction amplify and hybridise with just
fragment 6 to give


AC matrix	A	C
bio-A		
bio-C		

25


0944504-04340

Overall results**Fragment 4 TRSPA-1 analysis**

5

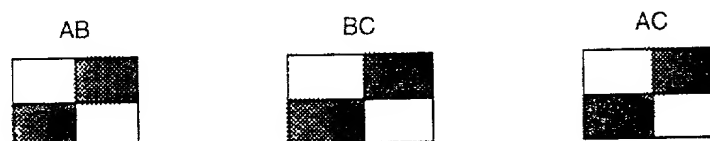Fragment 5 TRSPA-1 analysis**Fragment 6 TRSPA-1 analysis**

10

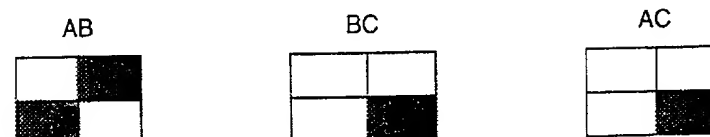**Version 2 - (TRSPA-2)**

For each (rxr) TRSPA test matrix, digest DNA to completion with all r restriction enzymes. Ligate on restriction enzyme cutting-site-specific adaptors with unique polymerase chain reaction amplifiable tags. Split into r^2 aliquots and for each aliquot, amplify with non-biotinylated restriction enzyme_j and non-biotinylated restriction enzyme_k tag-specific primers and array the denatured strands for all values of j and k between 1 and r .

20

- 40 -

Example

Consider the following dsDNA

--[A]--1--[B]--2--[B]--3--[C]--4--[A]--5--[B]-- +
 5 --[A]--1--[B]--2--[B]--3--[C]--4--[A]--5--[B]-- -

where the restriction enzyme cutting sites are denoted A, B and C and the fragments after restriction digestion are denoted 1, 2, 3, 4 and 5. + and - denote the sense of the strands.

10 Cut to completion with A, B and C to give

A--1--B B--2--B B--3--C C--4--A A--5--B +
 A--1--B B--2--B B--3--C C--4--A A--5--B -

15 Polymerase chain reaction amplify according to the following primer matrix to give

09944604-043402

- 41 -

	A	B	C
A		A--1--B + A--1--B - B--2--B + B--2--B - A--5--B + A--5--B -	C--4--A + C--4--A -
B	A--1--B + A--1--B - B--2--B + B--2--B - A--5--B + A--5--B -	B--2--B + B--2--B -	B--2--B + B--2--B - B--3--C + B--3--C -
C	C--4--A + C--4--A -	B--2--B + B--2--B - B--3--C + B--3--C -	

In simplified form

	A	B	C
A		1, 2, 5	4
B	1, 2, 5	2	2, 3
C	4	2, 3	

5

Repeat for all $^n C_r$ combinations of restriction enzymes to generate the TRSPA.

TRSPA-2 hybridisation patterns

The two types of TRSPA-2 hybridisation pattern we should expect using a probe resulting from inter-population perfectly matched duplex depletion are

1. Hybridisation to an off-diagonal element (e.g. row x, column y - where x and y are different) and its complementary element reflected across the diagonal (i.e. row y, column x) and
2. Hybridisation to a whole row and column intersecting at an on-diagonal element (e.g. all of row z and all of column z).

TRSPA-2 analysis - a worked example for $n=3$ and $r=2$

Let us take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

where

A, B and C denote restriction enzyme cutting sites. 1, 2, 3, 4, 5, 6, 7, 8 and 9 denote the restriction fragments after digestion.

The TRSPA-2 test matrices

There are ${}^3C_2=3$ TRSPA-2 test matrices - AB, BC and AC.

The test matrix hybridisation patterns

There are $\{(2 \cdot (2+1))/2\}=3$ possible hybridisation patterns for each of the test matrices



The combinatorial diversity

If there are three possible TRSPA-2 test matrix hybridisation

- 43 -

patterns and three possible test matrices, then there will be $3^3=27$ possible TRSPA-2 signatures.

Fragment 4 analysis - the AB TRSPA-2 matrix

5 Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and B to give

10

-C-1-B

B-2-B

B-3-A

A-4-A

15 A-5-C-6-C-7-B

B-8-A

A-9-A-

The fragment complementary to fragment 4 will be

20

A-4-A

Polymerase chain reaction amplify and hybridise with just fragment 4 to give

25

AB matrix	A	B
A		
B		

Fragment 4 analysis - the BC TRSPA-2 matrix

Take the dsDNA

- 44 -

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with B and C to give

-C-1-B

5 B-2-B

B-3-A-4-A-5-C

C-6-C

C-7-B

B-8-A-9-A-

10

The fragment complementary to fragment 4 will be

B-3-A-4-A-5-C

15

Polymerase chain reaction amplify and hybridise with just
fragment 4 to give

BC matrix	B	C
B		
C		

Fragment 4 analysis - the AC TRSPA-2 matrix

20

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

25

-C

C-1-B-2-B-3-A

A-4-A

- 45 -

A-5-C

C-6-C

C-7-B-8-A

A-9-A

5 A-

The fragment complementary to fragment 4 will be

A-4-A

10

Polymerase chain reaction amplify and hybridise with just
fragment 4 to give

AC matrix	A	C
A		
C		

15 **Fragment 5 analysis - the AB TRSPA-2 matrix**

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

20

Cut with A and B to give

-C-1-B

B-2-B

B-3-A

25 A-4-A

A-5-C-6-C-7-B

B-8-A

A-9-A-

05914604-01340

- 46 -

The fragment complementary to fragment 5 will be

A-5-C-6-C-7-B

- 5 Polymerase chain reaction amplify and hybridise with just fragment 5 to give

AB matrix	A	B
A		
B		

Fragment 5 analysis - the BC TRSPA-2 matrix

- 10 Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with B and C to give

15

-C-1-B

B-2-B

B-3-A-4-A-5-C

C-6-C

20 C-7-B

B-8-A-9-A-

The fragment complementary to fragment 5 will be

25 B-3-A-4-A-5-C

Polymerase chain reaction amplify and hybridise with just fragment 5 to give

- 47 -

BC matrix	B	C
B		
C		

Fragment 5 analysis - the AC TRSPA-2 matrix

Take the dsDNA

5

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

10 -C

C-1-B-2-B-3-A

A-4-A

A-5-C

C-6-C

15 C-7-B-8-A

A-9-A

A-

The fragment complementary to fragment 5 will be

20

A-5-C

Polymerase chain reaction amplify and hybridise with just
fragment 5 to give

AC matrix	A	C
A		
C		

25

09914604.0340

- 48 -

Fragment 6 analysis - the AB TRSPA-2 matrix

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

5

Cut with A and B to give

-C-1-B

B-2-B

10 B-3-A

A-4-A

A-5-C-6-C-7-B

B-8-A

A-9-A-

15

The fragment complementary to fragment 6 will be

A-5-C-6-C-7-B

20

Polymerase chain reaction amplify and hybridise with just
fragment 6 to give

AB matrix	A	B
A		
B		

Fragment 6 analysis - the BC TRSPA-2 matrix

25

Take the dsDNA

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

- 49 -

Cut with B and C to give

-C-1-B

B-2-B

5 B-3-A-4-A-5-C

C-6-C

C-7-B

B-8-A-9-A-

10 The fragment complementary to fragment 6 will be

C-6-C

15 Polymerase chain reaction amplify and hybridise with just
fragment 6 to give

BC matrix	B	C
B		
C		

Fragment 6 analysis - the AC TRSPA-2 matrix

Take the dsDNA

20

-C-1-B-2-B-3-A-4-A-5-C-6-C-7-B-8-A-9-A-

Cut with A and C to give

25 -C

C-1-B-2-B-3-A

A-4-A

A-5-C

094464-0440

- 50 -

C-6-C

C-7-B-8-A

A-9-A

A-

5

The fragment complementary to fragment 6 will be

C-6-C


10

Polymerase chain reaction amplify and hybridise with just fragment 6 to give


AC matrix	A	C
A		
C		

Overall results

15

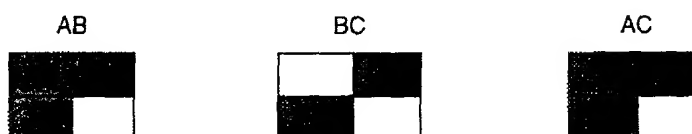Fragment 4 TRSPA-2 analysis**Fragment 5 TRSPA-2 analysis**

20

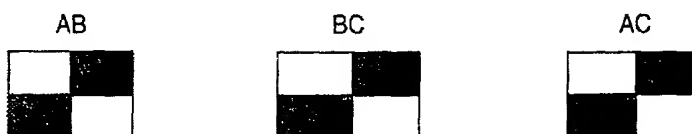

09944604-043402

$r^2 \cdot \{^nC_r\}$ for various n and r is given in the following table

R	n =							
	3	4	5	6	7	8	9	10
1	3	4	5	6	7	8	9	10
2	12	24	40	60	84	112	144	180
3	9	36	90	180	315	504	756	1,080
4	-	16	80	240	560	1,120	2,016	3,360
5	-	-	25	150	525	1,400	3,150	6,300
6	-	-	-	36	252	1,008	3,024	7,560
7	-	-	-	-	49	392	1,764	5,880
8	-	-	-	-	-	64	576	2,880
9	-	-	-	-	-	-	81	810
10	-	-	-	-	-	-	-	100

Each test as above will give rise to a particular hybridisation


5 'signature'.

TRSPA-1

(for $n=6$ and $r=3$)

TRSPA-2

(for $n=6$ and $r=3$)

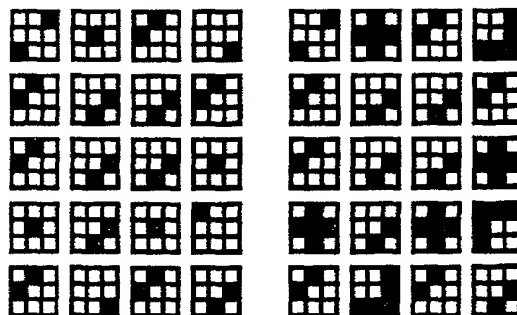

10

There will be $\{1+2+\dots+r\} = \{r(r+1)/2\}$ patterns per $(r \times r)$ test matrix. If there are nC_r $(r \times r)$ test matrices, the total possible number of signatures will be given by $\{r(r+1)/2\}$ raised to the power nC_r . We require

15 that the number of such different signatures is greatly in excess over the

number of fragments generated upon fragmentation.

If we choose $n=6$ and $r=3$, there will be 6 possible hybridisation patterns per test and 20 possible tests - giving a combinatorial diversity of $6^{20} \approx 3.7 \cdot 10^{15}$. The preferred scheme therefore employs six enzymes, tested in groups of three - giving 180 spots per experiment. In order to avoid restriction fragment length polymorphism problems, a duplicate analysis could be performed with a different set of enzymes sharing no cutting sites in common with the first set.

10 **Example restriction enzyme sets for the preparation of test matrices for TRSPAs**

Criteria for enzyme choice

In order to make TRSPAs, the selection of suitable enzymes is an important factor. Ideally, two sets of different enzymes are required to eliminate the small possibility that a phenotype-determining polymorphism might fall within a chosen restriction site and therefore compromise the specificity of the resulting signature. The selection of enzymes can be based upon a number of criteria

- 20 a) The enzymes should be 6 bp cutters.
- b) Cleavage by any selected enzyme should leave a 4 bp overhang at the 5' end.
- c) The selected enzymes in each set should all work efficiently under the same buffer conditions.
- 25 d) The selected enzymes in each set should ideally work efficiently at a single incubation temperature.
- e) The chosen enzymes should be commercially available - ideally at concentrations of 10 U / μ l or more.
- f) The 5' overhangs left by any two enzymes in the same set
- 30 should not be identical.
- g) No enzyme should appear in both sets for TRSPA fabrication.

h) Enzymes should be selected to avoid or minimise the effects of mammalian methylation patterns. In particular, enzymes with CG dinucleotides in their recognition sites should be avoided unless the enzyme is known to be able to restrict ^{m5}CpG sites.

5

DNA methylation

In vertebrates, DNA is often methylated at the 5th position of cytosine in the sequence of CpG and this is the only chemical modification that DNA of vertebrates contains under physiological conditions. By the careful selection of enzymes which do not contain CpG sequences within the recognition site, or the selection of enzymes which freely restrict ^{m5}CpG methylated sites, it is possible to remove the potentially adverse effects of DNA methylation from the TRSPA analysis. 6 bp cutters which are known to restrict ^{m5}CpG modified DNA efficiently to leave a 5' four base overhang are *BspEI* and *XmaI*. These enzymes are therefore the only enzymes with restriction sites containing CpG dinucleotides that are potentially useful in a TRSPA analysis.

10

15

Enzyme selection method

Sixteen possible four base pair overhangs exist (excluding unusual enzymes with asymmetrical recognition sequences such as *BssSI* or *BsiI*), five of which contain CG in the sequence. A further four overhangs could potentially contain CG sequences within the restriction recognition site if preceded by a C and followed by a G. Enzymes are therefore preferentially selected from the remaining six groups.

20

Excluding isoschizomers, there are up to four possible enzymes which would leave a particular 5' four bp overhang. For example, enzymes leaving a CTAG overhang are

25

30

- 55 -

AvrII - CCTAGG
NheI - GCTAGC
SpeI - ACTAGT
XbaI - TCTAGA

5

Enzymes to cleave sites with all the combinations of flanking bases are not available for all overhangs - hence the enzyme choice is more limited for some overhang groups than others.

As a primary step towards enzyme selection, the enzymes supplied by Amersham Pharmacia Biotech, New England Biolabs and Promega are ordered below by overhang sequence. Supplementary details such as the percentage activity in common buffers, the reaction temperature, concentration and supplier are also recorded. For three overhang sequences, there are no available enzymes, the remaining 13 are described below. Enzymes considered unsuitable due to methylation sensitivity are shaded (darker shade). Enzymes considered unfavourable due to the presence of CpG sites even though they do restrict methylated DNA to some degree are also shaded (lighter shaded).

20 **Candidate restriction enzymes leaving a 5' overhang of GTAC**

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>Acc65I</i>	GGTACC	10/50/100	37	10	NEB	65 / 20
<i>BsiWI</i>	CGTACG	100(50)/100(50)/25(-)	55 (37)	12	NEB	80 / 20
<i>BsrGI</i>	TGTACA	10/100/100	37	10	NEB	80 / 20
<i>SpII</i> (<i>BsiWI</i>)	CGTACG	100/20/60	55	10	APB	NI

00914504-013400

Candidate restriction enzymes leaving a 5' overhang of TTAA

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>AflIII</i>	CTTAAG	25/100/100	37	10	APB	60 / 15

Candidate restriction enzymes leaving a 5' overhang of AATT

5

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>EcoRI</i>	GAATTC	100/100/100	37	>40	APB/NEB	65 / 20
<i>MfeI</i>	CAATTG	10/50/100	37	10	NEB	65 / 20

Candidate restriction enzymes leaving a 5' overhang of CATG

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>BspHI</i>	TCATGA	50/100/100	37	10	NEB/APB	65 / 20
<i>NcoI</i>	CCATGG	100/100/100	37	50	NEB/APB	65 / 20

10 Candidate restriction enzymes leaving a 5' overhang of GATC

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>BamHI</i>	GGATCC	50/100/75	37	>40	APB/NEB	80 / 20
<i>BclI</i>	TGATCA	100/100/75	50 (37)	10	NEB	No
<i>BglII</i>	AGATCT	100/75/10	37	>40	APB/NEB	No

Candidate restriction enzymes leaving a 5' overhang of CCGG

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>BspEI</i>	TCCGGA	100/10/0	37	10	NEB	80 / 20
<i>NgoM IV</i>	GCCGGC	10/50/100	37	10	NEB	80 / 20
<i>XmaI</i>	CCCGGG	0/100/100	37	10	NEB	65 / 20

Candidate restriction enzymes leaving a 5' overhang of GCGC

5

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>KasI</i>	GGCGCC	75/100/75	37	5	NEB	65 / 20

Candidate restriction enzymes leaving a 5' overhang of TCGA

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>PaeR7I</i>	CTCGAG	10/100/100	37	20	NEB	80 / 20
<i>Sall</i>	GTCGAC	50/0/0	37	>40	APB/NEB	60 / 15
<i>XhoI</i>	CTCGAG	100/100/100	37	>40	APB/NEB	65 / 20

10 Candidate restriction enzymes leaving a 5' overhang of AGCT

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>HindIII</i>	AAGCTT	10/100/50	37	>40	NEB/APB	65 / 20

Candidate restriction enzymes leaving a 5' overhang of CGCG

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>BssHII</i>	GCGCGC	100/100/100	50	20	NEB/APB	80 / 20
<i>MluI</i>	ACGCGT	100/75/50	37	10	APB/NEB	65 / 20

Candidate restriction enzymes leaving a 5' overhang of GGCC

5

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>EagI</i>	CGGCCG	100/25/10	37	50	NEB	65 / 20
<i>Eco52I</i> (<i>EagI</i>)	CGGCCG	20/20/20	37	10	APB	60 / 15

Candidate restriction enzymes leaving a 5' overhang of TGCA

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>ApaLI</i>	GTGCAC	10/100/100	37	10	NEB/APB	70 / 15
<i>Alw44I</i>	GTGCAC	-/50/-	37	10	Promega	70 / 15

10 Candidate restriction enzymes leaving a 5' overhang of CTAG

Enzyme	Site	H/M/T buffer + BSA % activity	Optimu m °C	U / µl	Supplier	Inactivate °C/min
<i>AvrII</i>	CCTAGG	50/100/100	37	4	NEB	No
<i>NheI</i>	GCTAGC	10/100/100	37	10	APB/NEB	65 / 20
<i>SpeI</i>	ACTAGT	25/100/75	37	50	NEB	65 / 20
<i>XbaI</i>	TCTAGA	75/100/75	37	100	NEB	65 / 20
<i>BlnI</i> (<i>AvrII</i>)	CCTAGG	40/20/20	37	10	APB	No

- 59 -

From the above tables, a short-list of the most useful enzymes is given below for each of the buffer conditions shown

Overhang	Buffer H + BSA		Buffer M + BSA		Buffer T + BSA	
	Enzyme	Efficiency	Enzyme	Efficiency	Enzyme	Efficiency
AATT	<i>EcoRI</i>	100	<i>EcoRI</i>	100	<i>MfeI</i>	100
			<i>MfeI</i>	50	-	-
CATG	<i>NcoI</i>	100	<i>NcoI</i>	100	<i>NcoI</i>	100
	<i>BspHI</i>	50	<i>BspHI</i>	100	<i>BspHI</i>	100
GATC	<i>BglII</i>	100	<i>BamHI</i>	100	<i>BamHI</i>	75
	<i>BamHI</i>	50	<i>BglII</i>	75	-	-
TATA	-	-	-	-	-	-
ACGT	-	-	-	-	-	-
CCGG	<i>BspEI</i>	100	<i>XmaI</i>	100	<i>XmaI</i>	100
GCGC	-	-	-	-	-	-
TCGA	-	-	-	-	-	-
AGCT	-	-	<i>HindIII</i>	100	<i>HindIII</i>	50
CGCG	-	-	-	-	-	-
GGCC	-	-	-	-	-	-
TGCA	-	-	<i>ApaLI</i>	100	<i>ApaLI</i>	100
ATAT	-	-	-	-	-	-
CTAG	<i>XbaI</i>	75	<i>XbaI</i>	100	<i>NheI</i>	100
	<i>AvrII</i>	50	<i>SpeI</i>	100	<i>AvrII</i>	100
	<i>BlnI</i>	40	<i>NheI</i>	100	<i>XbaI</i>	75
	-	-	<i>AvrII</i>	100	<i>SpeI</i>	75
GTAC	<i>Acc65I</i>	100	<i>BsrGI</i>	100	<i>BsrGI</i>	100
			<i>Acc65I</i>	75	-	-
TTAA	-	-	<i>AflII</i>	100	<i>AflII</i>	100

- 60 -

- example selection of two six-enzyme sets is described below. Reserve enzymes, which could also be used, are shown as well. These reserve enzymes can be substituted (provided this does not cause overhang duplication) if practical problems regarding enzyme availability or performance should occur.

Example enzyme set 1 - for restriction in buffer M + BSA

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>BamHI</i>	GGATCC	100	37	100	NEB	80 / 20
2	<i>BsrGI</i>	TGTACA	100	37	10	NEB	80 / 20
3	<i>HindIII</i>	AAGCTT	100	37	>40	APB	65 / 20
4	<i>NcoI</i>	CCATGG	100	37	50	NEB	65 / 20
5	<i>SpeI</i>	ACTAGT	100	37	50	NEB	65 / 20
6	<i>AflII</i>	CTTAAG	100	37	10	APB	60 / 15

10 Example enzyme set 2 - for restriction in buffer M + BSA

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
1	<i>EcoRI</i>	GAATTC	100	37	>40	APB	65 / 20
2	<i>BspHI</i>	TCATGA	100	37	10	APB	65 / 20
3	<i>BglII</i>	AGATCT	75	37	>40	APB	No
4	<i>XbaI</i>	TCTAGA	100	37	100	NEB	65 / 20
5	<i>Acc65I</i>	GGTACC	75	37	10	NEB	65 / 20
6	<i>ApaI</i>	GTGCAC	100	37	10	NEB	No

1. Example reserve enzymes

Number	Enzyme	Site	M+BSA % activity	Optimum °C	U / μ l	Supplier	Inactivate °C/min
r1	<i>XmaI</i>	CCCGGG	100	37	10	NEB	65 / 20
r2	<i>BspEI</i>	TCCGGA	100	37	10	NEB	80 / 20
r3	<i>NheI</i>	GCTAGC	100	37	10	APB	65/20
r4	<i>AvrII</i>	CCTAGG	100	37	4	NEB	No

Example specific (3x3) terminal restriction site profiling array (TRSPA)

5 test matrices for above two sets of six enzymes

Triplet combinations from *BamHI*, *BsrGI*, *HindIII*, *NcoI*, *SpeI* and *AflIII*

For example set 1 - *BamHI*, *BsrGI*, *HindIII*, *NcoI*, *SpeI* and *AflIII* - the 20 ($= {}^6C_3$) triplet combinations are

10

<i>BamHI BsrGI HindIII</i>	<i>BamHI HindIII SpeI</i>	<i>BsrGI HindIII NcoI</i>	<i>BsrGI SpeI AflIII</i>
<i>BamHI BsrGI NcoI</i>	<i>BamHI HindIII AflIII</i>	<i>BsrGI HindIII SpeI</i>	<i>HindIII NcoI SpeI</i>
<i>BamHI BsrGI SpeI</i>	<i>BamHI NcoI SpeI</i>	<i>BsrGI HindIII AflIII</i>	<i>HindIII NcoI AflIII</i>
<i>BamHI BsrGI AflIII</i>	<i>BamHI NcoI AflIII</i>	<i>BsrGI NcoI SpeI</i>	<i>HindIII SpeI AflIII</i>
<i>BamHI HindIII NcoI</i>	<i>BamHI SpeI AflIII</i>	<i>BsrGI NcoI AflIII</i>	<i>NcoI SpeI AflIII</i>

Triplet combinations from *EcoRI*, *BspHI*, *BglII*, *XbaI*, *Acc65I* and *ApaLI*

For example set 2 - *EcoRI*, *BspHI*, *BglII*, *XbaI*, *Acc65I* and *ApaLI* - the 20 ($= {}^6C_3$) triplet combinations are

15

<i>EcoRI BspHI BglII</i>	<i>EcoRI BglII Acc65I</i>	<i>BspHI BglII XbaI</i>	<i>BspHI Acc65I ApaLI</i>
<i>EcoRI BspHI XbaI</i>	<i>EcoRI BglII ApaLI</i>	<i>BspHI BglII Acc65I</i>	<i>BglII XbaI Acc65I</i>
<i>EcoRI BspHI Acc65I</i>	<i>EcoRI XbaI Acc65I</i>	<i>BspHI BglII ApaLI</i>	<i>BglII XbaI ApaLI</i>
<i>EcoRI BspHI ApaLI</i>	<i>EcoRI XbaI ApaLI</i>	<i>BspHI XbaI Acc65I</i>	<i>BglII Acc65I ApaLI</i>
<i>EcoRI BglII XbaI</i>	<i>EcoRI Acc65I ApaLI</i>	<i>BspHI XbaI ApaLI</i>	<i>XbaI Acc65I ApaLI</i>

Example TRSPA-1 test matrices for set 1 - *Bam*HI, *Bsr*GI, *Hind*III, *Nco*I, *Spe*I and *Afl*II

matrix 1	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Hind</i> III-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Bsr</i> GI-primer	4	5	6
5'-biotin- <i>Hind</i> III-primer	7	8	9

matrix 2	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Bsr</i> GI-primer	4	5	6
5'-biotin- <i>Nco</i> I-primer	7	8	9

5

matrix 3	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Spe</i> I-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Bsr</i> GI-primer	4	5	6
5'-biotin- <i>Spe</i> I-primer	7	8	9

matrix 4	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Bsr</i> GI-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

matrix 5	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Nco</i> I-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Nco</i> I-primer	7	8	9

matrix 6	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Spe</i> I-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Spe</i> I-primer	7	8	9

matrix 7	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

10

2000041094660

matrix 8	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Spe</i> I-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Nco</i> I-primer	4	5	6
5'-biotin- <i>Spe</i> I-primer	7	8	9

matrix 9	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Nco</i> I-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

matrix 10	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Spe</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bam</i> HI-primer	1	2	3
5'-biotin- <i>Spe</i> I-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

matrix 11	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Nco</i> I-primer
5'-biotin- <i>Bsr</i> GI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Nco</i> I-primer	7	8	9

matrix 12	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Spe</i> I-primer
5'-biotin- <i>Bsr</i> GI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Spe</i> I-primer	7	8	9

5

matrix 13	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bsr</i> GI-primer	1	2	3
5'-biotin- <i>Hind</i> III-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

matrix 14	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Spe</i> I-primer
5'-biotin- <i>Bsr</i> GI-primer	1	2	3
5'-biotin- <i>Nco</i> I-primer	4	5	6
5'-biotin- <i>Spe</i> I-primer	7	8	9

matrix 15	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-biotin- <i>Bsr</i> GI-primer	1	2	3
5'-biotin- <i>Nco</i> I-primer	4	5	6
5'-biotin- <i>Afl</i> II-primer	7	8	9

- 64 -

matrix 16	5'-HO- <i>BsrGI</i> -primer	5'-HO- <i>SpeI</i> -primer	5'-HO- <i>AflII</i> -primer
5'-biotin- <i>BsrGI</i> -primer	1	2	3
5'-biotin- <i>SpeI</i> -primer	4	5	6
5'-biotin- <i>AflII</i> -primer	7	8	9

matrix 17	5'-HO- <i>HindIII</i> -primer	5'-HO- <i>NcoI</i> -primer	5'-HO- <i>SpeI</i> -primer
5'-biotin- <i>HindIII</i> -primer	1	2	3
5'-biotin- <i>NcoI</i> -primer	4	5	6
5'-biotin- <i>SpeI</i> -primer	7	8	9

matrix 18	5'-HO- <i>HindIII</i> -primer	5'-HO- <i>NcoI</i> -primer	5'-HO- <i>AflII</i> -primer
5'-biotin- <i>HindIII</i> -primer	1	2	3
5'-biotin- <i>NcoI</i> -primer	4	5	6
5'-biotin- <i>AflII</i> -primer	7	8	9

matrix 19	5'-HO- <i>HindIII</i> -primer	5'-HO- <i>SpeI</i> -primer	5'-HO- <i>AflII</i> -primer
5'-biotin- <i>HindIII</i> -primer	1	2	3
5'-biotin- <i>SpeI</i> -primer	4	5	6
5'-biotin- <i>AflII</i> -primer	7	8	9

matrix 20	5'-HO- <i>NcoI</i> -primer	5'-HO- <i>SpeI</i> -primer	5'-HO- <i>AflII</i> -primer
5'-biotin- <i>NcoI</i> -primer	1	2	3
5'-biotin- <i>SpeI</i> -primer	4	5	6
5'-biotin- <i>AflII</i> -primer	7	8	9

5

Example TRSPA-1 test matrices for set 2 - *EcoRI*, *BspHI*, *BglII*, *XbaI*, *Acc65I* and *ApaLI*

matrix 21	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>BglII</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BspHI</i> -primer	4	5	6
5'-biotin- <i>BglII</i> -primer	7	8	9

matrix 22	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>XbaI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BspHI</i> -primer	4	5	6
5'-biotin- <i>XbaI</i> -primer	7	8	9

10

matrix 23	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BspHI</i> -primer	4	5	6
5'-biotin- <i>Acc65I</i> -primer	7	8	9

matrix 24	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BspHI</i> -primer	4	5	6
5'-biotin- <i>ApaLI</i> -primer	7	8	9

matrix 25	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>XbaI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BglII</i> -primer	4	5	6
5'-biotin- <i>XbaI</i> -primer	7	8	9

matrix 26	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BglII</i> -primer	4	5	6
5'-biotin- <i>Acc65I</i> -primer	7	8	9

matrix 27	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>BglII</i> -primer	4	5	6
5'-biotin- <i>ApaLI</i> -primer	7	8	9

5

matrix 28	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>XbaI</i> -primer	4	5	6
5'-biotin- <i>Acc65I</i> -primer	7	8	9
matrix 29	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>XbaI</i> -primer	4	5	6
5'-biotin- <i>ApaLI</i> -primer	7	8	9

matrix 30	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>Acc65I</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-biotin- <i>EcoRI</i> -primer	1	2	3
5'-biotin- <i>Acc65I</i> -primer	4	5	6
5'-biotin- <i>ApaLI</i> -primer	7	8	9

matrix 31	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Xba</i> I-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Bgl</i> II-primer	4	5	6
5'-biotin- <i>Xba</i> I-primer	7	8	9

matrix 32	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Acc</i> 65I-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Bgl</i> II-primer	4	5	6
5'-biotin- <i>Acc</i> 65I-primer	7	8	9

matrix 33	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Apa</i> LI-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Bgl</i> II-primer	4	5	6
5'-biotin- <i>Apa</i> LI-primer	7	8	9

matrix 34	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Acc</i> 65I-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Xba</i> I-primer	4	5	6
5'-biotin- <i>Acc</i> 65I-primer	7	8	9

matrix 35	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Apa</i> LI-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Xba</i> I-primer	4	5	6
5'-biotin- <i>Apa</i> LI-primer	7	8	9

5

matrix 36	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Acc</i> 65I-primer	5'-HO- <i>Apa</i> LI-primer
5'-biotin- <i>Bsp</i> HI-primer	1	2	3
5'-biotin- <i>Acc</i> 65I-primer	4	5	6
5'-biotin- <i>Apa</i> LI-primer	7	8	9

matrix 37	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Acc</i> 65I-primer
5'-biotin- <i>Bgl</i> II-primer	1	2	3
5'-biotin- <i>Xba</i> I-primer	4	5	6
5'-biotin- <i>Acc</i> 65I-primer	7	8	9

matrix 38	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Apa</i> LI-primer
5'-biotin- <i>Bgl</i> II-primer	1	2	3
5'-biotin- <i>Xba</i> I-primer	4	5	6
5'-biotin- <i>Apa</i> LI-primer	7	8	9

matrix 39	5'-HO- <i>BglII</i> -primer	5'-HO- <i>Acc65I</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-biotin- <i>BglII</i> -primer	1	2	3
5'-biotin- <i>Acc65I</i> -primer	4	5	6
5'-biotin- <i>ApaLI</i> -primer	7	8	9

matrix 40	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Acc65</i> I-primer	5'-HO- <i>Apa</i> LI-primer
5'-biotin- <i>Xba</i> I-primer	1	2	3
5'-biotin- <i>Acc65</i> I-primer	4	5	6
5'-biotin- <i>Apa</i> LI-primer	7	8	9

Example TRSPA-2 test matrices for set 1 - *Bam*HI, *Bsr*GI, *Hind*III, *Nco*I, *Spe*I and *Afl*II

5

matrix 1	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Hind</i> III-primer
5'-HO- <i>Bam</i> HI-primer	1	2	3
5'-HO- <i>Bsr</i> GI-primer	4	5	6
5'-HO- <i>Hind</i> III-primer	7	8	9

matrix 2	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer
5'-HO- <i>Bam</i> HI-primer	1	2	3
5'-HO- <i>Bsr</i> GI-primer	4	5	6
5'-HO- <i>Nco</i> I-primer	7	8	9

matrix 3	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Spe</i> I-primer
5'-HO- <i>Bam</i> HI-primer	1	2	3
5'-HO- <i>Bsr</i> GI-primer	4	5	6
5'-HO- <i>Spe</i> I-primer	7	8	9

matrix 4	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Bam</i> HI-primer	1	2	3
5'-HO- <i>Bsr</i> GI-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

matrix 5	5'-HO- <i>Bam</i> HI-primer	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Nco</i> I-primer
5'-HO- <i>Bam</i> HI-primer	1	2	3
5'-HO- <i>Hind</i> III-primer	4	5	6
5'-HO- <i>Nco</i> I-primer	7	8	9

matrix 6	5'-HO-BamHI-primer	5'-HO-HindIII-primer	5'-HO-SpeI-primer
5'-HO-BamHI-primer	1	2	3
5'-HO-HindIII-primer	4	5	6
5'-HO-SpeI-primer	7	8	9

matrix 7	5'-HO-BamHI-primer	5'-HO-HindIII-primer	5'-HO-AflII-primer
5'-HO-BamHI-primer	1	2	3
5'-HO-HindIII-primer	4	5	6
5'-HO-AflII-primer	7	8	9

matrix 8	5'-HO-BamHI-primer	5'-HO-NcoI-primer	5'-HO-SpeI-primer
5'-HO-BamHI-primer	1	2	3
5'-HO-NcoI-primer	4	5	6
5'-HO-SpeI-primer	7	8	9

matrix 9	5'-HO-BamHI-primer	5'-HO-NcoI-primer	5'-HO-AflII-primer
5'-HO-BamHI-primer	1	2	3
5'-HO-NcoI-primer	4	5	6
5'-HO-AflII-primer	7	8	9

matrix 10	5'-HO-BamHI-primer	5'-HO-SpeI-primer	5'-HO-AflII-primer
5'-HO-BamHI-primer	1	2	3
5'-HO-SpeI-primer	4	5	6
5'-HO-AflII-primer	7	8	9

5

matrix 11	5'-HO-BsrGI-primer	5'-HO-HindIII-primer	5'-HO-NcoI-primer
5'-HO-BsrGI-primer	1	2	3
5'-HO-HindIII-primer	4	5	6
5'-HO-NcoI-primer	7	8	9

matrix 12	5'-HO-BsrGI-primer	5'-HO-HindIII-primer	5'-HO-SpeI-primer
5'-HO-BsrGI-primer	1	2	3
5'-HO-HindIII-primer	4	5	6
5'-HO-SpeI-primer	7	8	9

matrix 13	5'-HO-BsrGI-primer	5'-HO-HindIII-primer	5'-HO-AflII-primer
5'-HO-BsrGI-primer	1	2	3
5'-HO-HindIII-primer	4	5	6
5'-HO-AflII-primer	7	8	9

matrix 14	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Spe</i> I-primer
5'-HO- <i>Bsr</i> GI-primer	1	2	3
5'-HO- <i>Nco</i> I-primer	4	5	6
5'-HO- <i>Spe</i> I-primer	7	8	9

matrix 15	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Bsr</i> GI-primer	1	2	3
5'-HO- <i>Nco</i> I-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

matrix 16	5'-HO- <i>Bsr</i> GI-primer	5'-HO- <i>Spe</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Bsr</i> GI-primer	1	2	3
5'-HO- <i>Spe</i> I-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

matrix 17	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Spe</i> I-primer
5'-HO- <i>Hind</i> III-primer	1	2	3
5'-HO- <i>Nco</i> I-primer	4	5	6
5'-HO- <i>Spe</i> I-primer	7	8	9

matrix 18	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Hind</i> III-primer	1	2	3
5'-HO- <i>Nco</i> I-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

5

matrix 19	5'-HO- <i>Hind</i> III-primer	5'-HO- <i>Spe</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Hind</i> III-primer	1	2	3
5'-HO- <i>Spe</i> I-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

matrix 20	5'-HO- <i>Nco</i> I-primer	5'-HO- <i>Spe</i> I-primer	5'-HO- <i>Afl</i> II-primer
5'-HO- <i>Nco</i> I-primer	1	2	3
5'-HO- <i>Spe</i> I-primer	4	5	6
5'-HO- <i>Afl</i> II-primer	7	8	9

- 70 -

Example TRSPA-2 test matrices for set 2 - *EcoRI*, *BspHI*, *BglII*, *XbaI*,
Acc65I and *ApaLI*

matrix 21	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>BglII</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BspHI</i> -primer	4	5	6
5'-HO- <i>BglII</i> -primer	7	8	9

matrix 22	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>XbaI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BspHI</i> -primer	4	5	6
5'-HO- <i>XbaI</i> -primer	7	8	9

matrix 23	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BspHI</i> -primer	4	5	6
5'-HO- <i>Acc65I</i> -primer	7	8	9

matrix 24	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BspHI</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

matrix 25	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>XbaI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>XbaI</i> -primer	7	8	9

matrix 26	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>Acc65I</i> -primer	7	8	9

matrix 27	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

matrix 28	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>XbaI</i> -primer	4	5	6
5'-HO- <i>Acc65I</i> -primer	7	8	9

matrix 29	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>XbaI</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

matrix 30	5'-HO- <i>EcoRI</i> -primer	5'-HO- <i>Acc65I</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>EcoRI</i> -primer	1	2	3
5'-HO- <i>Acc65I</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

matrix 31	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>XbaI</i> -primer
5'-HO- <i>BspHI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>XbaI</i> -primer	7	8	9

matrix 32	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-HO- <i>BspHI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>Acc65I</i> -primer	7	8	9

5

matrix 33	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>BglII</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>BspHI</i> -primer	1	2	3
5'-HO- <i>BglII</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

matrix 34	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>Acc65I</i> -primer
5'-HO- <i>BspHI</i> -primer	1	2	3
5'-HO- <i>XbaI</i> -primer	4	5	6
5'-HO- <i>Acc65I</i> -primer	7	8	9

matrix 35	5'-HO- <i>BspHI</i> -primer	5'-HO- <i>XbaI</i> -primer	5'-HO- <i>ApaLI</i> -primer
5'-HO- <i>BspHI</i> -primer	1	2	3
5'-HO- <i>XbaI</i> -primer	4	5	6
5'-HO- <i>ApaLI</i> -primer	7	8	9

- 72 -

matrix 36	5'-HO- <i>Bsp</i> HI-primer	5'-HO- <i>Acc</i> 65I-primer	5'-HO- <i>Apa</i> LI-primer
5'-HO- <i>Bsp</i> HI-primer	1	2	3
5'-HO- <i>Acc</i> 65I-primer	4	5	6
5'-HO- <i>Apa</i> LI-primer	7	8	9

matrix 37	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Acc</i> 65I-primer
5'-HO- <i>Bgl</i> II-primer	1	2	3
5'-HO- <i>Xba</i> I-primer	4	5	6
5'-HO- <i>Acc</i> 65I-primer	7	8	9

matrix 38	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Apa</i> LI-primer
5'-HO- <i>Bgl</i> II-primer	1	2	3
5'-HO- <i>Xba</i> I-primer	4	5	6
5'-HO- <i>Apa</i> LI-primer	7	8	9

matrix 39	5'-HO- <i>Bgl</i> II-primer	5'-HO- <i>Acc</i> 65I-primer	5'-HO- <i>Apa</i> LI-primer
5'-HO- <i>Bgl</i> II-primer	1	2	3
5'-HO- <i>Acc</i> 65I-primer	4	5	6
5'-HO- <i>Apa</i> LI-primer	7	8	9

matrix 40	5'-HO- <i>Xba</i> I-primer	5'-HO- <i>Acc</i> 65I-primer	5'-HO- <i>Apa</i> LI-primer
5'-HO- <i>Xba</i> I-primer	1	2	3
5'-HO- <i>Acc</i> 65I-primer	4	5	6
5'-HO- <i>Apa</i> LI-primer	7	8	9

5

Hybridisation patterns

There are six possible hybridisation patterns for a given probe fragment from inter-population perfectly matched duplex depletion enrichment with a (3x3) test matrix in TRSPA-1 analysis.

10

There are also six possible hybridisation patterns for a given probe fragment from inter-population perfectly matched duplex depletion enrichment with a (3x3) test matrix in TRSPA-2 analysis.

We can denote these patterns as follows

15

- 73 -

TRSPA-1 analysis**Pattern a (1 only)**

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer		2	3
5'-biotin-Y-primer	4	5	6
5'-biotin-Z-primer	7	8	9

5

Pattern b (2 and 4)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer	1		3
5'-biotin-Y-primer		5	6
5'-biotin-Z-primer	7	8	9

Pattern c (3 and 7)

10

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer	1	2	
5'-biotin-Y-primer	4	5	6
5'-biotin-Z-primer		8	9

Pattern d (5 only)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer	1	2	3
5'-biotin-Y-primer	4		6
5'-biotin-Z-primer	7	8	9

15

20 T E O 4 0 9 1 6 0

- 74 -

Pattern e (6 and 8)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer	1	2	3
5'-biotin-Y-primer	4	5	
5'-biotin-Z-primer	7		9

Pattern f (9 only)

5

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-biotin-X-primer	1	2	3
5'-biotin-Y-primer	4	5	6
5'-biotin-Z-primer	7	8	

TRSPA-2 analysis**Pattern a (row 1 and column 1)**

10

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-HO-X-primer			
5'-HO-Y-primer		5	6
5'-HO-Z-primer		8	9

Pattern b (2 and 4)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'-HO-X-primer	1		3
5'-HO-Y-primer		5	6
5'-HO-Z-primer	7	8	9

15

2025-04-10 10:54:00

- 75 -

Pattern c (3 and 7)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'- HO -X-primer	1	2	
5'- HO -Y-primer	4	5	6
5'- HO -Z-primer		8	9

Pattern d (row 2 and column 2)

5

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'- HO -X-primer	1		3
5'- HO -Y-primer			
5'- HO -Z-primer	7		9

Pattern e (6 and 8)

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'- HO -X-primer	1	2	3
5'- HO -Y-primer	4	5	
5'- HO -Z-primer	7		9

10 **Pattern f (row 3 and column 3)**

matrix #	5'-HO-X-primer	5'-HO-Y-primer	5'-HO-Z-primer
5'- HO -X-primer	1	2	
5'- HO -Y-primer	4	5	
5'- HO -Z-primer			

Signatures

If we have 40 such (3x3) test matrices and we denote the possible hybridisation patterns a, b, c, d, e and f, then we can write the overall TRSPA hybridisation signature as follows

15

Matrix	Pattern	matrix	pattern
1	a, b, c, d, e or f	21	a, b, c, d, e or f
2	a, b, c, d, e or f	22	a, b, c, d, e or f
3	a, b, c, d, e or f	23	a, b, c, d, e or f
4	a, b, c, d, e or f	24	a, b, c, d, e or f
5	a, b, c, d, e or f	25	a, b, c, d, e or f
6	a, b, c, d, e or f	26	a, b, c, d, e or f
7	a, b, c, d, e or f	27	a, b, c, d, e or f
8	a, b, c, d, e or f	28	a, b, c, d, e or f
9	a, b, c, d, e or f	29	a, b, c, d, e or f
10	a, b, c, d, e or f	30	a, b, c, d, e or f
11	a, b, c, d, e or f	31	a, b, c, d, e or f
12	a, b, c, d, e or f	32	a, b, c, d, e or f
13	a, b, c, d, e or f	33	a, b, c, d, e or f
14	a, b, c, d, e or f	34	a, b, c, d, e or f
15	a, b, c, d, e or f	35	a, b, c, d, e or f
16	a, b, c, d, e or f	36	a, b, c, d, e or f
17	a, b, c, d, e or f	37	a, b, c, d, e or f
18	a, b, c, d, e or f	38	a, b, c, d, e or f
19	a, b, c, d, e or f	39	a, b, c, d, e or f
20	a, b, c, d, e or f	40	a, b, c, d, e or f

Aspect (III)

- In another aspect the invention provides a nucleic acid
- 5 characterisation method which comprises presenting to the set of arrays as defined above a nucleic acid fragment of interest under hybridisation conditions, and observing a pattern of hybridisation. Preferably, a plurality of nucleic acid fragments of interest are separately presented to the set of arrays, and the resulting patterns of hybridisation are compared.
 - 10 Preferably, the plurality of nucleic acid fragments of interest are drawn from the mixture of DNA fragments, enriched in fragments that are characteristic of a phenotype of interest, as described under the invention (1) above.

Thus in this aspect the invention provides a method of identifying fragments of DNA that are characteristic of a phenotype of

interest, which method comprises recovering, cloning and amplifying individual DNA fragments from the mixture of DNA fragments obtained under invention (1) above, presenting the individual DNA fragments to the set of arrays as defined under hybridisation conditions, observing a pattern of hybridisation generated by each individual DNA fragment. and subjecting to further investigation any two or more individual DNA fragments whose hybridisation patterns are similar or identical.

TRSPA signatures and whole genome association studies

After a given number of cycles of inter-population perfectly matched duplex depletion, phenotype determining fragments will be enriched but will not be entirely free from 'noise' fragments. Noise may result from unequal allelic frequencies for certain SNPs between the two populations. Noise will also result from the presence of somatic mutations in the cells used to prepare DNA fragments and from the use of polymerase chain reaction in some of the embodiments of the current invention.

For the preferred embodiment, DNA is prepared from a library of clones (either genomic clones or cDNA clones) - with inserts derived from the individual(s) and propagated in some appropriate host and cloning vector system. Restriction enzyme fragmentation is used prior to cloning and polymerase chain reaction amplification is used to prepare the DNA for comparison in fragmented form. Priming sites within the vector sequence flanking the cloned restriction enzyme fragmented inserts are employed for one or more cycles of polymerase chain reaction amplification of the fragmented DNA of interest. The primers used for polymerase chain reaction amplification of the fragmented DNA of interest are again used after the phenotype-determining fragment enrichment process to 'rescue' and clone the enriched fragments. Cloned enriched fragments are colony purified, picked into appropriate storage containers, catalogued and archived, DNA probes are prepared from these single clones and are

individually hybridised against 180 spot TRSPA arrays as above. TRSPA signatures are determined for many colonies. Most noise fragments will be random in nature and will thus be randomly distributed amongst the $6^{20} = 3.7 \cdot 10^{15}$ possible types of TRSPA signature. Phenotype-determining signal fragments, however, will be those where repeat TRSPA signatures are obtained as more and more colonies are sampled. The more frequently the repetition of a given TRSPA signature occurs per unit colonies sampled, the greater the signal to noise ratio and the more successful has been the enrichment.

10 A great many of the steps in this process are amenable to high throughput automation - enabling very large numbers of single colony TRSPA signatures to be determined with ease and extending the power of the current invention to cases where signal to noise ratios are beyond current approaches.

15 Statistical correlations (associations) can initially be drawn between TRSPA signatures and phenotype. The clones giving rise to a particular TRSPA signature showing a useful association with a phenotype of interest can then be sequenced in order to determine at a DNA sequence level the association(s) with the phenotype of interest. Such associations have future predictive values for the phenotype of interest, knowing the genotype and will be of great value in medicine and pharmacogenetics.

25 If the genome of interest is wholly or partially sequenced, we can also *in silico* restrict the DNA with all n enzymes, calculate the expected signature for each fragment and pattern match these expected signatures with the observed signature (taking into account any loss or gain of restriction sites due to polymorphic variation compared to the reference sequence) to immediately identify the fragment of interest within a gene, genomic region, chromosome or whole genome. This latter method will be of great value in those cases where a great many phenotype-determining fragments are obtained and repeat signatures are rare or unobtained. The

30

clustering of phenotype-determining fragments to adjacent DNA regions thus gives an association between those genomic regions and the phenotype of interest.

5 **Aspect (IV)**

In yet another aspect the invention provides a double-stranded DNA molecule having the sequence a-A-b-B...X-y-Y-z where A, B...X and Y are unique restriction sites for n different restriction endonuclease enzymes, and a, b...y, z denotes distances in base pairs, characterised in that each fragment, obtainable by cutting the DNA molecule by means of any one or more up to n of the restriction enzymes, has a different length from every other fragment.

15 **An example totally diagnostic internal control DNA which allows both the extent and exact nature of any example set 1 (or example set 2) 6 bp cutter partial digestion to be unambiguously determined for inter-population perfectly matched duplex depletion or TRSPA restriction**

In both of the above schemes, it is important that limit digestion products are obtained. Monitoring the extent of partial digestion resulting from multi-enzymatic restriction and determining precisely which enzymes have failed to cut is a task of great importance.

20 If we have up to six enzymes for DNA digestion - let us label these A, B, C, D, E and F. We need to somehow determine that these have all cut to completion during the fragmentation stage for inter-population perfectly matched duplex depletion and also for the digestion step prior to adaptor ligation in TRSPA fabrication. If any of the enzymes have failed to cut to completion, we need to know which ones and to what degree in order to effectively rectify the problem.

30 **The structure of the internal control DNA**

If we construct a double stranded DNA molecule with the

- 80 -

following structure

end---t---A---u---B---v---C---w---D---x---E---y---F---z---end

5 where the A, B, C, D, E and F denote the sites for restriction enzyme cutting and t, u, v, w, x, y and z denote distances in base pairs.

This internal control DNA is either uniformly pre-labelled and added to the DNA of interest at an appropriate concentration prior to restriction or is Southern blot probed with a complementary sequence not
10 found in the DNA of interest after restriction.

All six enzymes can cut in only one way.

One enzyme can fail to cut in ${}^6C_1 = 6$ ways, these are: A, B, C, D, E or F failing to cut.

Two enzymes can fail to cut in ${}^6C_2 = 15$ ways, these are: AB, AC, AD, AE, AF, BC, BD, BE, BF, CD, CE, CF, DE, DF or EF failing to cut.
15

Three enzymes can fail to cut in ${}^6C_3 = 20$ ways, these are: ABC, ABD, ABE, ABF, ACD, ACE, ACF, ADE, ADF, AEF, BCD, BCE, BCF, BDE, BDF, BEF, CDE, CDF, CEF or DEF failing to cut.

Four enzymes can fail to cut in ${}^6C_4 = 15$ ways, these are:
20 ABCD, ABCE, ABCF, ABDE, ABDF, ABEF, ACDE, ACDF, ACEF, ADEF, BCDE, BCDF, BCEF, BDEF or CDEF failing to cut.

Five enzymes can fail to cut in ${}^6C_5 = 6$ ways, these are: ABCDE, ABCDF, ABCEF, ABDEF, ACDEF or BCDEF failing to cut.

All six enzymes can fail to cut in only one way.

25 Each of the above possibilities will generate one or more fragments from the internal control DNA. If each possible fragment has a discernible size from any other, then we can determine exactly which enzymes have cut and which have not from the size distribution of the fragments generated. The task is therefore to design such a DNA
30 molecule.

0914504-03400

- 81 -

Example simulations

Seven simulations are given below - varying the size of inter-site fragments. Criteria for a successful outcome include the following

1. The inter-fragment spacing should be greater for larger fragments (so as to aid electrophoretic resolution).
2. All possible fragments should be unambiguously resolvable in size from each-other.
3. Size gaps between bands comprising different numbers of inter-site units should be greater than the size gaps between bands comprising the same number of inter-site units.
4. The size gaps and size spread from largest to smallest fragment should be electrophoretically compatible.

Simulation 1**Inter-site fragment sizes (in bp)**

End---A	A---B	B---C	C---D	D---E	E---F	F---end
80	90	100	110	120	130	140

Possible digestion products obtained (in bp)

	Fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	80	
A---B	90	10
B---C	100	10
C---D	110	10
D---E	120	10
E---F	130	10
F---end	140	10
two unit fragments		
end---B	170	30
A---C	190	20
B---D	210	20
C---E	230	20
D---F	250	20
E---end	270	20
three unit fragments		
end---C	270	0
A---D	300	30
B---E	330	30
C---F	360	30
D---end	390	30
four unit fragments		
end---D	380	-10
A---E	420	40
B---F	460	40
C---end	500	40
five unit fragments		
end---E	500	0
A---F	550	50
B---end	600	50
six unit fragments		
end---F	630	30
B---end	690	60
seven unit fragments		
end---end	770	80

05914604.013102

- 83 -

Spread = 690 bp

Simulation 2**5 Inter-site fragment sizes (in bp)**

end---A	A---B	B---C	C---D	D---E	E---F	F---end
90	100	110	120	130	140	150

Possible digestion products obtained (in bp)0944604-043403
2012-10-24 10:54:50

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	90	
A---B	100	10
B---C	110	10
C---D	120	10
D---E	130	10
E---F	140	10
F---end	150	10
two unit fragments		
end---B	190	40
A---C	210	20
B---D	230	20
C---E	250	20
D---F	270	20
E---end	290	20
three unit fragments		
end---C	300	10
A---D	330	30
B---E	360	30
C---F	390	30
D---end	420	30
four unit fragments		
end---D	420	0
A---E	460	40
B---F	500	40
C---end	540	40
five unit fragments		
end---E	550	10
A---F	600	50
B---end	650	50
six unit fragments		
end---F	690	40
B---end	750	60
seven unit fragments		
end---end	840	90

Spread = 750 bp

0944604-0409T650

Simulation 3

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---D	D---E	E---F	F---end
100	110	120	130	140	150	160

5

Possible digestion products obtained (in bp)

00914604-013400

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	100	
A---B	110	10
B---C	120	10
C---D	130	10
D---E	140	10
E---F	150	10
F---end	160	10
two unit fragments		
end---B	210	50
A---C	230	20
B---D	250	20
C---E	270	20
D---F	290	20
E---end	310	20
three unit fragments		
end---C	330	20
A---D	360	30
B---E	390	30
C---F	420	30
D---end	450	30
four unit fragments		
end---D	460	10
A---E	500	40
B---F	540	40
C---end	580	40
five unit fragments		
end---E	600	20
A---F	650	50
B---end	700	50
six unit fragments		
end---F	750	50
B---end	810	60
seven unit fragments		
end---end	910	100

Spread = 810 bp

0944304 033

Simulation 4

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---D	D---E	E---F	F---end
110	120	130	140	150	160	170

5

Possible digestion products obtained (in bp)

0944504.04469

- 88 -

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	110	
A---B	120	10
B---C	130	10
C---D	140	10
D---E	150	10
E---F	160	10
F---end	170	10
two unit fragments		
end---B	230	60
A---C	250	20
B---D	270	20
C---E	290	20
D---F	310	20
E---end	330	20
three unit fragments		
end---C	360	30
A---D	390	30
B---E	420	30
C---F	450	30
D---end	480	30
four unit fragments		
end---D	500	20
A---E	540	40
B---F	580	40
C---end	620	40
five unit fragments		
end---E	650	30
A---F	700	50
B---end	750	50
six unit fragments		
end---F	810	60
B---end	870	60
seven unit fragments		
end---end	980	110

Spread = 870 bp

05944604-013102

- 89 -

Simulation 5**Inter-site fragment sizes (in bp)**

end---A	A---B	B---C	C---D	D---E	E---F	F---end
120	130	140	150	160	170	180

5

Possible digestion products obtained (in bp)

03944604-019102

- 90 -

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	120	
A---B	130	10
B---C	140	10
C---D	150	10
D---E	160	10
E---F	170	10
F---end	180	10
two unit fragments		
end---B	250	70
A---C	270	20
B---D	290	20
C---E	310	20
D---F	330	20
E---end	350	20
three unit fragments		
end---C	390	40
A---D	420	30
B---E	450	30
C---F	480	30
D---end	510	30
four unit fragments		
end---D	540	30
A---E	580	40
B---F	620	40
C---end	660	40
five unit fragments		
end---E	700	40
A---F	750	50
B---end	800	50
six unit fragments		
end---F	870	70
B---end	930	60
seven unit fragments		
end---end	1050	120

Spread = 930 bp

09044666 044666

- 91 -

Simulation 6

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---D	D---E	E---F	F---end
130	140	150	160	170	180	190

5

Possible digestion products obtained (in bp)

09914604-04910
207670-70947560

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	130	
A---B	140	10
B---C	150	10
C---D	160	10
D---E	170	10
E---F	180	10
F---end	190	10
two unit fragments		
end---B	270	80
A---C	290	20
B---D	310	20
C---E	330	20
D---F	350	20
E---end	370	20
three unit fragments		
end---C	420	50
A---D	450	30
B---E	480	30
C---F	510	30
D---end	540	30
four unit fragments		
end---D	580	40
A---E	620	40
B---F	660	40
C---end	700	40
five unit fragments		
end---E	750	50
A---F	800	50
B---end	850	50
six unit fragments		
end---F	930	80
B---end	990	60
seven unit fragments		
end---end	1120	130

Spread = 990 bp

0044604-094760

Simulation 7

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---D	D---E	E---F	F---end
140	150	160	170	180	190	200

5

Possible digestion products obtained (in bp)

0914604-0440

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	140	
A---B	150	10
B---C	160	10
C---D	170	10
D---E	180	10
E---F	190	10
F---end	200	10
two unit fragments		
end---B	290	90
A---C	310	20
B---D	330	20
C---E	350	20
D---F	370	20
E---end	390	20
three unit fragments		
end---C	450	60
A---D	480	30
B---E	510	30
C---F	540	30
D---end	570	30
four unit fragments		
end---D	620	50
A---E	660	40
B---F	700	40
C---end	740	40
five unit fragments		
end---E	800	60
A---F	850	50
B---end	900	50
six unit fragments		
end---F	990	90
B---end	1050	60
seven unit fragments		
end---end	1190	140

Spread = 1050 bp

0994604-0400

According to the above criteria for success, simulation 7 above clearly fulfils all of the requirements.

An example totally diagnostic internal control DNA which allows both the extent and exact nature of any example set 1 (or example set 2) 4 bp cutter partial digestion to be unambiguously determined for inter-population perfectly matched duplex depletion

If we have up to three enzymes for DNA digestion - let us label these A, B and C. We need to somehow determine that these have all cut to completion during the fragmentation stage for inter-population perfectly matched duplex depletion. If any of the enzymes have failed to cut to completion, we need to know which ones and to what degree in order to effectively rectify the problem.

The structure of the internal control DNA

If we construct a double stranded DNA molecule with the following structure

end---t---A---u---B---v---C---w---end

where the A, B and C denote the sites for restriction enzyme cutting and t, u, v and w denote distances in base pairs.

This internal control DNA is uniformly pre-labelled and added to the DNA of interest at an appropriate concentration prior to restriction or is Southern blot probed with a complementary sequence not found in the DNA of interest after restriction.

All three enzymes can cut in only one way.

One enzyme can fail to cut in ${}^3C_1 = 3$ ways, these are: A, B or C failing to cut.

Two enzymes can fail to cut in ${}^3C_2 = 3$ ways, these are: AB, AC or BC failing to cut.

03914504-013102

- 96 -

All three enzymes can fail to cut in only one way.

Each of the above possibilities will generate one or more fragments from the internal control DNA. If each possible fragment has a discernible size from any other (and from any of the fragments in simulation 5 7 above for the panel of up to 6 enzymes), then we can determine exactly which enzymes have cut and which have not from the size distribution of the fragments generated. The task is therefore to design such a DNA molecule.

10 Example simulations

Six simulations are given below - varying the size of inter-site fragments. Criteria for a successful outcome included the following:

1. The inter-fragment spacing should be greater for larger fragments (so as to aid electrophoretic resolution).
- 15 2. All possible fragments should be unambiguously resolvable in size from each-other.
3. Size gaps between bands comprising different numbers of inter-site units should ideally be greater than the size gaps between bands comprising the same number of inter-site units.
- 20 4. The size gaps and size spread from largest to smallest fragment should be electrophoretically compatible.
5. The largest fragment obtained should ideally be smaller than the smallest fragment obtained in simulation 7 above for the panel of up to six enzymes.

25

- 97 -

Simulation 1**Inter-site fragment sizes (in bp)**

end---A	A---B	B---C	C---end
10	20	30	40

5

Possible digestion products obtained (in bp)

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	10	
A---B	20	10
B---C	30	10
C---end	40	10
two unit fragments		
end---B	30	-10
A---C	50	20
B---end	70	20
three unit fragments		
end---C	60	-10
A---end	90	30
four unit fragments		
end---end	100	10

Spread = 90 bp

10

Simulation 2**Inter-site fragment sizes (in bp)**

end---A	A---B	B---C	C---end
15	25	35	45

15

Possible digestion products obtained (in bp)

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	15	
A---B	25	10
B---C	35	10
C---end	45	10
two unit fragments		
end---B	40	-5
A---C	60	20
B---end	80	20
three unit fragments		
end---C	75	-5
A---end	105	30
four unit fragments		
end---end	120	15

Spread = 105 bp

5

Simulation 3

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---end
20	30	40	50

10

- 99 -

Possible digestion products obtained (in bp)

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	20	
A---B	30	10
B---C	40	10
C---end	50	10
two unit fragments		
end---B	50	0
A---C	70	20
B---end	90	20
three unit fragments		
end---C	90	0
A---end	120	30
four unit fragments		
end---end	140	20

Spread = 120 bp

5

Simulation 4

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---end
25	35	45	55

10

09944604 0109

5

Simulation 5

Inter-site fragment sizes (in bp)

10

- 101 -

Possible digestion products obtained (in bp)

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	20	
A---B	25	5
B---C	30	5
C---end	35	5
two unit fragments		
end---B	45	10
A---C	55	10
B---end	65	10
three unit fragments		
end---C	75	10
A---end	90	15
four unit fragments		
end---end	110	20

Spread = 90 bp

5

Simulation 6

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---end
25	30	35	40

10

0944604-0340

- 102 -

Possible digestion products obtained (in bp)

	fragment length (in bp)	size gap to the next smallest fragment (in bp)
one unit fragments		
end---A	25	
A---B	30	5
B---C	35	5
C---end	40	5
two unit fragments		
end---B	55	15
A---C	65	10
B---end	75	10
three unit fragments		
end---C	90	15
A---end	105	15
four unit fragments		
end---end	130	25

Spread = 105 bp

5

According to the above criteria for success, simulation 6 clearly fulfils all of the requirements.

Example determination of the entire set of internal control DNA limit and partial digestion patterns for a panel of up to six restriction enzymes

For the example simulation 7, the entire set of internal control DNA limit and partial digestion patterns for a panel of up to six restriction enzymes can be determined as below.

15

- 103 -

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---D	D---E	E---F	F---end
140	150	160	170	180	190	200

Possible digestion products obtained

5 All six enzymes can cut in only one way.

failed digests	bp	bp	bp	bp	bp	bp	bp
none	140	150	160	170	180	190	200

One enzyme can fail to cut in ${}^6C_1 = 6$ ways, these are: A, B, C, D, E or F failing to cut.

10

failed digest	bp	bp	bp	bp	bp	bp
A	290	160	170	180	190	200
B	140	310	170	180	190	200
C	140	150	330	180	190	200
D	140	150	160	350	190	200
E	140	150	160	170	370	200
F	140	150	160	170	180	390

20140404T0600

- 104 -

Two enzymes can fail to cut in ${}^6C_2 = 15$ ways, these are: AB, AC, AD, AE, AF, BC, BD, BE, BF, CD, CE, CF, DE, DF or EF failing to cut.

failed digests	bp	bp	bp	bp	bp
AB	450	170	180	190	200
AC	290	330	180	190	200
AD	290	160	350	190	200
AE	290	160	170	370	200
AF	290	160	170	180	390
BC	140	480	180	190	200
BD	140	310	350	190	200
BE	140	310	170	370	200
BF	140	310	170	180	390
CD	140	150	510	190	200
CE	140	150	330	370	200
CF	140	150	330	180	390
DE	140	150	160	540	200
DF	140	150	160	350	390
EF	140	150	160	170	570

- 105 -

Three enzymes can fail to cut in ${}^6C_3 = 20$ ways, these are:
 ABC, ABD, ABE, ABF, ACD, ACE, ACF, ADE, ADF, AEF, BCD, BCE, BCF,
 BDE, BDF, BEF, CDE, CDF, CEF or DEF failing to cut.

failed digests	bp	bp	bp	bp
ABC	620	180	190	200
ABD	450	350	190	200
ABE	450	170	370	200
ABF	450	170	180	390
ACD	290	510	190	200
ACE	290	330	370	200
ACF	290	330	180	390
ADE	290	160	540	200
ADF	290	160	350	390
AEF	290	160	170	570
BCD	140	660	190	200
BCE	140	480	370	200
BCF	140	480	180	390
BDE	140	310	540	200
BDF	140	310	350	390
BEF	140	310	170	570
CDE	140	150	700	200
CDF	140	150	510	390
CEF	140	150	330	570
DEF	140	150	160	740

- 106 -

Four enzymes can fail to cut in ${}^6C_4 = 15$ ways, these are:
ABCD, ABCE, ABCF, ABDE, ABDF, ABEF, ACDE, ACDF, ACEF, ADEF,
BCDE, BCDF, BCEF, BDEF or CDEF failing to cut.

failed digests	bp	bp	bp
ACCD	800	190	200
ABCE	620	370	200
ABCF	620	180	390
ABDE	450	540	200
ABDF	450	350	390
ABEF	450	170	570
ACDE	290	700	200
ACDF	290	510	390
ACEF	290	330	570
ADEF	290	160	740
BCDE	140	850	200
BCDF	140	660	390
BCEF	140	480	570
BDEF	140	310	740
CDEF	140	150	900

- 107 -

Five enzymes can fail to cut in ${}^6C_5 = 6$ ways, these are:
 ABCDE, ABCDF, ABCEF, ABDEF, ACDEF or BCDEF failing to cut.

failed digests	bp	bp
ABCDE	990	200
ABCDF	800	390
ABCEF	620	570
ABDEF	450	740
ACDEF	290	900
BCDEF	140	1050

5

All six enzymes can fail to cut in only one way.

failed digests	bp
all	1190

**Example determination of the entire set of internal control DNA limit
 and partial digestion patterns for a panel of up to three restriction
 enzymes**

10

For the example simulation 6, the entire set of internal control
 DNA limit and partial digestion patterns for a panel of up to three restriction
 enzymes can be determined as below.

15

Inter-site fragment sizes (in bp)

end---A	A---B	B---C	C---end
25	30	35	40

Possible digestion products obtained

All three enzymes can cut in only one way.

20

failed digests	bp	bp	bp	bp
none	25	30	35	40

- 108 -

One enzyme can fail to cut in ${}^3C_1 = 3$ ways, these are: A, B or C failing to cut.

failed digests	bp	bp	bp
A	55	35	40
B	25	65	40
C	25	30	75

5 Two enzymes can fail to cut in ${}^3C_2 = 3$ ways, these are: AB, AC or BC failing to cut.

failed digests	bp	bp
AB	90	40
AC	55	75
BC	25	105

All three enzymes can fail to cut in only one way.


10

failed digests	bp
all	130

09914604.013102

Example 1a – The digestion internal control plasmid for the 6 bp
cutter set 1 TRSPA enzymes *Bam*HI, *Bsr*GI, *Hind*III, *Nco*I, *Spe*I, and
*Afl*III

The plasmid pNW33 (shown below) was constructed to
5 contain an insert with all of the 6 bp cutter TRSPA enzyme sites.

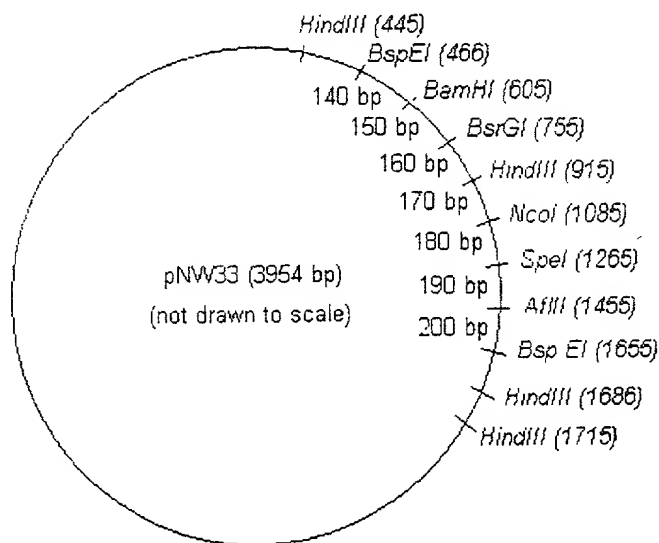

BspEI sites define the outer ends of the 140 bp and the 200
10 bp fragments. The full sequence for pNW33 is shown below:

```

tcgcgcgttcggtgatgacggtgaaaacctctgacacatgcagctcccgagacggtcacagctgtctgtaagcggatgcc
gggagcagacaagcccgtagggcgctcagcgggtgtggcgggtgtcggggctggcttaactatcgggcatcagagca
gattgtactgagatgcaccatatcggtgtgaaataccgcacagatgcgtaaggagaaaataccgcatcaggcgccattc
15 gccattcaggctgcgcaactgttgggaagggcgatcggtgcgggcctcttcgctattacgccagctggcgaaaggggatgt
gctgcaaggcgattaagtgtggaacgccagggtttccagtcacgacgttgtaaacgacggccagtgaattcgagctcgg
taccgggccccccctcgaggtcgacggtatcgataagcttgatcgacgttgtaatccggacgcccgcgtgaagatgttg
ggtgttgtaacaatatcgattccaattcagcggggccacctgatctcttgtatttaataaagactcaagcgggtcaactatga
agaagtgtctcttcgtcccagtaaggatccgcactttgaattttgtaatcctgaagggatcgtaaaaacagctcttctcaaatc
20 atacattaagacgactcgaaatccacatatcaaatatccgagtgtagtaaacattccaaaaccgtgatggaatggaacaaca
cttaaatgtacaccctggaatccggttttagaatccatgataataattttctggatttggtaattttttgcacgttcaaaattttgca
accccttttggaaacaaacactacggtaggctcggaatgtcactactgttgagcaattcacgttcattataagctttcactgcac
acgacgattctgtgattgtattcagcccatatcgttcatagctctcgccaaccgaacggacattcgaagtattccgcgtacgtg

```

atgttcacctcgatatgcatcgtgtaaaagcaattgtccagggaaccagggcgatctctcatagccatggaatacgccttttc
agtggtgcgatgctaataccggttacaaatattccgagcaccaagaatggctgcgcgcttgccctggacttgacgtcgatttgacg
gggtccttgagaaagtatttaaactggaacacaaatctgaggaatgatcaaaagcaaccaacgccaacgcataataactagt
caataaccaagacctcccaataatagcaccagactgtgtaataacctctggctctgatattgtccagatggaattggagcat
5 atggctcattaattgcgtcgataicicacataaccagtcgttgattgtctgtatagccagtaagacaaggaccagacatcatca
tgcaaaagaatcgcttaagccctcttgccctttataggatctctctgattttctgctcgagtttccggtaagacctttcggactt
cgtccacaaacacaaactcctcgcaacttttcgcggttggtacttgactggccacgtaatccacgatctcttttccgctacgt
ctttccgtgctccaaaacaacacggcgggcggtccggtaccagctgcgatcaagcttaicgataccgtcgacctcgacct
gcaggcatgcaagctggcgtaatcatggtcatagctgttctgtgaaattgtatccgctcacaattccacacacatacga
10 gccggaagcataaagtgtaaagcctggggtgctaatgagtgaactacacattatgcttgctcactgcccgtttc
cagtcgggaaacctgtcgtgccagctgcattaatgaatcgcccaacgcggggagagggcgtttgctgattggcgctcttc
cgcttctcgtcactgactcgtcgtcgtcggtcgttcggctgcggcgagcgggtatcagctcactcaaggcggttaacgggt
atccacagaatcaggggataacgcaggaagaacatgtgagcaaaaggccagcaaaaggccaggaaaccgtaaaaag
ggcggttgctggcgttttccataggtcgcgcctccctgacgagcatcacaataatcgacgctcaagtcagaggtggcgaa
15 acccgacaggactataagataccaggcgttccccctggaagctccctcgtcgtcctcttccgacctgcccgttaccg
gatacctgtccgcttttcccttcgggaagcgtggcgcttctcatagctcacgctgtaggtatcagttcggtgtaggtcgttcg
ctccaagctggcgtgtgacgaacccccgttcagcccgaaccgtcgccttatccggtaactatcgtcttgagccaacc
ggtaagacacgactatcgccactggcagcagccactggttaacaggattagcagagcaggtatgtagggcgtgctacag
agttctgaagtgtggcctaactacggctacactagaaggacagtaattggatctcgcgtcgtgtaagccagttaccttcgga
20 aaaagagttgtagctcttgatccggcaacaaaccaccgctggtagcgggtgtttttgttgcaagcagcagattacgcgc
agaaaaaaaggatcgaagaatccttgatcttttaccggggtgacgctcagtggaacgaaacacagtttaagggtt
ttggtcatgagattacaaaaagatcttcacatagatcttttaataaaaaatgaagtttaatacaatcaaatatataatgagt
aaacttggtcgtacagttaccaatgcttaacagtgaggcacctatctcagcgaatctgtctatttcgttcatcattgctgact
ccccgtcgtgtagataactacgatacgggaggggttaccatctggccccagtgctgcaatgataccgcgagaccacgctca
25 ccggctccagatttatcagcaataaaccagccagccggaaggcgagcgcagaagtggtcctgcaactttatccgcctcc
atccagctctattaattgttccgggaagctagagtaagtagtccagttatagtttgcgcaacgttgttgcattgctacaggc
atcgttggtgcacgctcgtcgttttgtaiggttcattcagctccgggttcccaacgatcaaggcgagttacatgatccccatgttgt
gcaaaaaagcggtagctccttcggtcctccgatcgtgtcagaagtaagttggccgcagtggtatcactcatggttatggcagc
actgcataaattcttactgtcatgccatccgaagaigtcttctgtgactggtgagtactcaaccaagtcattctgagaafagtga
30 tgcggcgaccgagttgctcttcccggtcaatacgggataataccgcgccacatagcagaactttaaagtgtcatcatt
ggaaaacgttctcggggcgaaaactcgaagatcttaccgctgttgagatccagttcgaatgaaccactcgtgcaccaaac
tgatcttcagcatctttactttaccagcgtttctgggtgagcaaaaacaggaaggcaaaatgccgcaaaaaagggaataag
ggcgacacggaaatgttgaatactcactcttcttcaatattatgaagcattatcaggggtattgtctcatgagcggtataca
tatttgatgtatttagaaaaataaacaatagggttccgcgcacattccccgaaaagtccacctgacgtctaagaaacc
35 attattatcatgacattaacctataaaaaataggcgatcacaggccctttcgtc

0914504-013100

- 111 -

The insert was introduced as a four fragment ligation of 140 (KpnI-BamHI), 150 / 160 (BamHI-HindIII), 170 / 180 (HindIII-SpeI) and 190 / 200 (SpeI-XhoI) into KpnI / Sall digested pUC19c DNA (Genbank X02514). The Sall and XhoI sites were lost as a result of the joining of their compatible sticky ends. The insert can be removed with BspEI or PvuII.

The insert region of 1190 bp and the short flanking regions to the vector junctions were sequenced twice in each direction in order to establish the plasmid sequence. In addition, a total of 63 analytical restriction digests and one minus-enzyme control were performed as detailed in the following table:

Digest	Well	Expected band sizes						
		bp	bp	bp	bp	bp	bp	bp
Uncut	1	1190						
BamHI	2	140	1050					
BsrGI	3	290	900					
HindIII	4	450	740					
NcoI	5	620	570					
SpeI	6	800	390					
AflIII	7	990	200					
BamHI BsrGI	8	140	150	900				
BamHI HindIII	9	140	310	740				
BamHI NcoI	10	140	480	570				
BamHI SpeI	11	140	660	390				
BamHI AflIII	12	140	850	200				
BsrGI HindIII	13	290	160	740				
BsrGI NcoI	14	290	330	570				
BsrGI SpeI	15	290	510	390				
BsrGI AflIII	16	290	700	200				
HindIII NcoI	17	450	170	570				
HindIII SpeI	18	450	350	390				
HindIII AflIII	19	450	540	200				
NcoI SpeI	20	620	180	390				
NcoI AflIII	21	620	370	200				
SpeI AflIII	22	800	190	200				
BamHI BsrGI HindIII	23	140	150	160	740			
BamHI BsrGI NcoI	24	140	150	330	570			
BamHI BsrGI SpeI	25	140	150	510	390			
BamHI BsrGI AflIII	26	140	150	700	200			

<i>BamHI HindIII NcoI</i>	27	140	310	170	570			
<i>BamHI HindIII SpeI</i>	28	140	310	350	390			
<i>BamHI HindIII AflII</i>	29	140	310	540	200			
<i>BamHI NcoI SpeI</i>	30	140	480	180	390			
<i>BamHI NcoI AflII</i>	31	140	480	370	200			
<i>BamHI SpeI AflII</i>	32	140	660	190	200			
<i>BsrGI HindIII NcoI</i>	33	290	160	170	570			
<i>BsrGI HindIII SpeI</i>	34	290	160	350	390			
<i>BsrGI HindIII AflII</i>	35	290	160	540	200			
<i>BsrGI NcoI SpeI</i>	36	290	330	180	390			
<i>BsrGI NcoI AflII</i>	37	290	330	370	200			
<i>BsrGI SpeI AflII</i>	38	290	510	190	200			
<i>HindIII NcoI SpeI</i>	39	450	170	180	390			
<i>HindIII NcoI AflII</i>	40	450	170	370	200			
<i>HindIII SpeI AflII</i>	41	450	350	190	200			
<i>NcoI SpeI AflII</i>	42	620	180	190	200			
<i>BamHI BsrGI HindIII NcoI</i>	43	140	150	160	170	570		
<i>BamHI BsrGI HindIII SpeI</i>	44	140	150	160	350	390		
<i>BamHI BsrGI HindIII AflII</i>	45	140	150	160	540	200		
<i>BamHI BsrGI NcoI SpeI</i>	46	140	150	330	180	390		
<i>BamHI BsrGI NcoI AflII</i>	47	140	150	330	370	200		
<i>BamHI BsrGI SpeI AflII</i>	48	140	150	510	190	200		
<i>BamHI HindIII NcoI SpeI</i>	49	140	310	170	180	390		
<i>BamHI HindIII NcoI AflII</i>	50	140	310	170	370	200		
<i>BamHI HindIII SpeI AflII</i>	51	140	310	350	190	200		
<i>BamHI NcoI SpeI AflII</i>	52	140	480	180	190	200		
<i>BsrGI HindIII NcoI SpeI</i>	53	290	160	170	180	390		
<i>BsrGI HindIII NcoI AflII</i>	54	290	160	170	370	200		
<i>BsrGI HindIII SpeI AflII</i>	55	290	160	350	190	200		
<i>BsrGI NcoI SpeI AflII</i>	56	290	330	180	190	200		
<i>HindIII NcoI SpeI AflII</i>	57	450	170	180	190	200		
<i>BamHI BsrGI HindIII NcoI SpeI</i>	58	140	150	160	170	180	390	
<i>BamHI BsrGI HindIII NcoI AflII</i>	59	140	150	160	170	370	200	
<i>BamHI BsrGI HindIII SpeI AflII</i>	60	140	150	160	350	190	200	
<i>BamHI BsrGI NcoI SpeI AflII</i>	61	140	150	330	180	190	200	
<i>BamHI HindIII NcoI SpeI AflII</i>	62	140	310	170	180	190	200	
<i>BsrGI HindIII NcoI SpeI AflII</i>	63	290	160	170	180	190	200	
<i>BamHI BsrGI HindIII NcoI SpeI AflII</i>	64	140	150	160	170	180	190	200

All of these digests produced the expected fragment patterns on agarose gel electrophoresis and a number of these, shaded in the table above, are shown in figure 1. The restriction digests illustrated in figure 1

- 113 -

were carried out using the following conditions:

Digest 1 (Minus enzyme control)

10 μ l	100 μ g/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 μ g/ml BSA
90 μ l	1x NEB buffer #2 + 100 μ g/ml BSA

5

Digest 5

15 μ l	0.33 U/ μ l <i>Nco</i> I in 1x NEB buffer #2 + 100 μ g/ml BSA
10 μ l	100 μ g/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 μ g/ml BSA
75 μ l	1x NEB buffer #2 + 100 μ g/ml BSA

Digest 9

10

15 μ l	0.33 U/ μ l <i>Bam</i> HI in 1x NEB buffer #2 + 100 μ g/ml BSA
15 μ l	0.33 U/ μ l <i>Hind</i> III in 1x NEB buffer #2 + 100 μ g/ml BSA
10 μ l	100 μ g/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 μ g/ml BSA
60 μ l	1x NEB buffer #2 + 100 μ g/ml BSA

Digest 37

15 μ l	0.33 U/ μ l <i>Bsr</i> GI in 1x NEB buffer #2 + 100 μ g/ml BSA
15 μ l	0.33 U/ μ l <i>Nco</i> I in 1x NEB buffer #2 + 100 μ g/ml BSA
15 μ l	0.33 U/ μ l <i>A</i> III in 1x NEB buffer #2 + 100 μ g/ml BSA
10 μ l	100 μ g/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 μ g/ml BSA
45 μ l	1x NEB buffer #2 + 100 μ g/ml BSA

15

2004-04-09 10:44:50

Digest 49

15 µl	0.33 U/µl <i>Bam</i> HI in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Hind</i> III in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Nco</i> I in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Spe</i> I in 1x NEB buffer #2 + 100 µg/ml BSA
10 µl	100 µg/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 µg/ml BSA
30 µl	1x NEB buffer #2 + 100 µg/ml BSA

Digest 61


15 µl	0.33 U/µl <i>Bam</i> HI in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Bsr</i> GI in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Nco</i> I in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Spe</i> I in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	0.33 U/µl <i>Afl</i> III in 1x NEB buffer #2 + 100 µg/ml BSA
10 µl	100 µg/ml internal control plasmid DNA in 1x NEB buffer #2 + 100 µg/ml BSA
15 µl	1x NEB buffer #2 + 100 µg/ml BSA

All 64 restriction digests were incubated at 37°C for 6 hours and samples were then electrophoresed on 2.5 % FMC MetaPhor agarose gels as illustrated in figure 1. The lanes marked M contain mixed samples from all of the digests 1-64 and these lanes were used as size markers after confirmation of the fragment sizes against Stratagene Kb ladder markers.

When spiked into a genomic DNA digest the internal control restriction fragment pattern produced is indicative of both the degree of digestion and the nature of any partial restriction at less than limit digestion. It is possible to deduce from the bands present which, if any, of the enzymes have failed to cut and therefore to take action to correct this before the DNA is used in a subsequent analysis or enrichment procedure.

Preparation of *BspEI*-released internal control DNA on a large scale

The method used for generating internal control PCR product DNA (free from contaminating dNTPs) is depicted in the following figure:

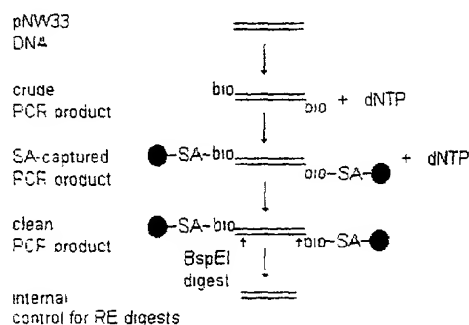

5

Primers and PCR

20 μ M BIO140UP

10 5' biotin-CGCAGCTGGTAATCCGGACGCCGCGTCGAAGATGTT 3'

20 μ M BIO200DOWN

5' biotin-CGCAGCTGGTAATCCGGACGCCGCGCGTTGTTGTT 3'

15

Bulk PCR amplification (192x 100 μ l reactions) was carried out according to the conditions described below:

	per 20 ml	final conditions
100 mM dATP	20 μ l	100 μ M
100 mM dCTP	20 μ l	100 μ M
100 mM dGTP	20 μ l	100 μ M
100 mM dTTP	20 μ l	100 μ M
10x PCR buffer	2 ml	1x
25 mM MgCl ₂	400 μ l	1/50th volume
1 ng/ μ l pNW33	400 μ l	2 ng per PCR
20 μ M BIO140UP	200 μ l	20 pmol per PCR
20 μ M BIO200DOWN	200 μ l	20 pmol per PCR
Taq DNA polymerase	100 μ l	2.5 U per PCR

- 116 -

The master mix was rapidly dispensed into PCR tubes. Thermal cycling was initiated using the following parameters: 30 cycles of 97°C for 1 min, 50°C for 2 min, and 72°C for 3 min; 72°C for 5 min; and
5 then 4°C.

After PCR amplification, samples were pooled and subjected to capture of biotinylated PCR product termini and BspEI release of internal control DNA.

10 **Capture of biotinylated PCR product termini and *BspEI* release of internal control DNA**

All separations were carried out using a Dynal MPC-1 separator (Dynal, product #12001).

20 ml of pooled PCR reaction were mixed with 20 ml of
15 Dynabeads M-280 in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M NaCl. The tube was incubated at room temperature for 1 hour with rolling.

The Dynabeads were then washed four times with 20 ml of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M NaCl.

BspEI digestion was carried out for 1 hour at 37°C in 2 ml of
20 0.5 U/μl BspEI in 1x NEB buffer #3. The digest was then ethanol precipitated by the addition of 0.1 volume of 3 M sodium acetate (pH 5.2) and 2.5 volumes of ethanol, chilling to -20°C and then centrifugation. Pellets were rinsed with 70 % ethanol prior to redissolution in 500 μl of 1x TE buffer.

25 1 μl of the BspEI-released internal control DNA was mixed with 10 μl of 50 % glycerol AGE loading dye and electrophoresed on a 1.5 % agarose gel to confirm that the size of the purified DNA was in accordance with that expected.

30

0914504.043402

- 117 -

High molecular weight genomic DNA digestion in the presence of internal control DNA with a dilution series of a mixture of the six set 1 6 bp cutters

High molecular weight canine genomic DNA was first mixed
 5 with a dilution series of a mixture of the six set 1 6 bp cutters – each at the same number of units. Aliquots were then removed and mixed with the BspEI-released internal control DNA described above. 20 µg of canine genomic DNA was digested with 0.25, 0.025, 0.0025, 0.00025, 0.000025, 0.0000025 and 0 U/µl BamHI / BsrGI / HindIII / NcoI / SpeI / AflII in 200 µl.
 10 0.4 µg of canine genomic DNA and 1 µl of BspEI-released internal control DNA were digested with 0.25, 0.025, 0.0025, 0.00025, 0.000025, 0.0000025 and 0 U/µl BamHI / BsrGI / HindIII / NcoI / SpeI / AflII in 4 µl. Vistra Green staining was used to monitor the extent of internal control DNA digestion and high molecular weight DNA digestion. The fragment
 15 pattern from the internal control DNA is diagnostic of both the degree of digestion and the nature of any partial restriction at less than limit digestion.

A premix of restriction enzymes, buffer and BSA was prepared as detailed below:

component	µl
100 U/µl <i>BamHI</i>	1 µl
10 U/µl <i>BsrGI</i>	10 µl
40 U/µl <i>HindIII</i>	2.5 µl
50 U/µl <i>NcoI</i>	2 µl
50 U/µl <i>SpeI</i>	2 µl
20 U/µl <i>AflII</i>	5 µl
10x NEB buffer #2	4 µl
100x BSA	0.4 µl
water	13.1 µl

20

Total volume = 40 µl

The premix therefore contained each restriction enzyme at

- 118 -

2.5 U/ μ l in 1x NEB buffer #2 and 1x BSA. Serial 10-fold dilutions of this premix were then prepared in 1x NEB buffer #2 and 1x BSA.

Premixes of canine genomic DNA, buffer and BSA were also prepared as detailed below:

5

Component	μ l
1 mg/ml canine genomic DNA	20 μ l
10x NEB buffer #2	8 μ l
100x BSA	0.8 μ l
Water	51.2 μ l

These premixes therefore contain canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA.

Canine genomic DNA and restriction enzyme mixes were
10 then set up as follows:

tube	component #1	component #2
1	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 2.5 U/ μ l in 1x NEB buffer #2 and 1x BSA
2	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 0.25 U/ μ l in 1x NEB buffer #2 and 1x BSA
3	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 0.025 U/ μ l in 1x NEB buffer #2 and 1x BSA
4	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 0.0025 U/ μ l in 1x NEB buffer #2 and 1x BSA
5	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 0.00025 U/ μ l in 1x NEB buffer #2 and 1x BSA
6	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of each restriction enzyme at 0.000025 U/ μ l in 1x NEB buffer #2 and 1x BSA
7	80 μ l of canine genomic DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	20 μ l of 1x NEB buffer #2 and 1x BSA

20070404 013102

- 119 -

2 μ l aliquots were then removed from tubes 1-7 and added to
1 μ l of BspEI-released internal control DNA and 1 μ l of 2x NEB buffer #2
and 2x BSA to give samples 1-7ic. Samples 1-7ic were then overlaid with
5 50 μ l of mineral oil in order to prevent evaporation.

The remaining volume from tubes 1-7 was then added to
100 μ l of 1x NEB buffer #2 and 1x BSA.

All samples were finally incubated at 37°C overnight.

After digestion, 20 μ l of digests 1-7 were mixed with 10 μ l of
10 50 % glycerol AGE loading dye and 4 μ l of digests 1-7ic were mixed with
2 μ l of 50 % glycerol AGE loading dye. Digests in loading dye were then
electrophoresed on a 2.5 % MetaPhor agarose gel in 1x TBE. The gel was
stained for 60 min in 500 ml of 1x TBE containing 50 μ l of Vistra Green.
The stained gel was finally imaged on a Fluorimager with the following
15 settings: a 488 nm laser; a 570 DF 30 filter; a PMT setting of 700 V;
200 μ m resolution; and low sensitivity.

Restriction digests - final conditions

Restriction digests 1-7 were therefore performed under the
20 following conditions:

tube	μ g canine genomic DNA	U <i>Bam</i> H1	U <i>Bsr</i> GI	U <i>Hind</i> III	U <i>Nco</i> I	U <i>Spe</i> I	U <i>Afl</i> III	NEB buffer #2	μ g/ml BSA
1	20	50	50	50	50	50	50	1x	100
2	20	5	5	5	5	5	5	1x	100
3	20	0.5	0.5	0.5	0.5	0.5	0.5	1x	100
4	20	0.05	0.05	0.05	0.05	0.05	0.05	1x	100
5	20	0.005	0.005	0.005	0.005	0.005	0.005	1x	100
6	20	0.0005	0.0005	0.0005	0.0005	0.0005	0.0005	1x	100
7	20	0	0	0	0	0	0	1x	100

Total volume = 200 μ l

0944604-043460

- 120 -

Likewise, restriction digests 1-7ic were performed under the following conditions:

tube	µg canine genomic DNA	µl IC DNA	U <i>Bam</i> H1	U <i>Bsr</i> GI	U <i>Hind</i> III	U <i>Nco</i> I	U <i>Spe</i> I	U <i>Afl</i> II	NEB buffer #2	µg/ml BSA
1ic	0.4	1	1	1	1	1	1	1	1x	100
2ic	0.4	1	0.1	0.1	0.1	0.1	0.1	0.1	1x	100
3ic	0.4	1	0.01	0.01	0.01	0.01	0.01	0.01	1x	100
4ic	0.4	1	0.001	0.001	0.001	0.001	0.001	0.001	1x	100
5ic	0.4	1	0.0001	0.0001	0.0001	0.0001	0.0001	0.0001	1x	100
6ic	0.4	1	0.00001	0.00001	0.00001	0.00001	0.00001	0.00001	1x	100
7ic	0.4	1	0	0	0	0	0	0	1x	100

5

IC = BspEI-released internal control DNA


Total volume = 4 µl

The results are shown in figure 2.

2000-04-04 10:04:04

Example 1b – The digestion internal control plasmid for the 4 bp cutter set 1 TRSPA enzymes *HaeIII*, *MboI*, and *MseI*

The plasmid pNW35 (shown below) was constructed to contain an insert with all of the 4 bp cutter TRSPA enzyme sites.

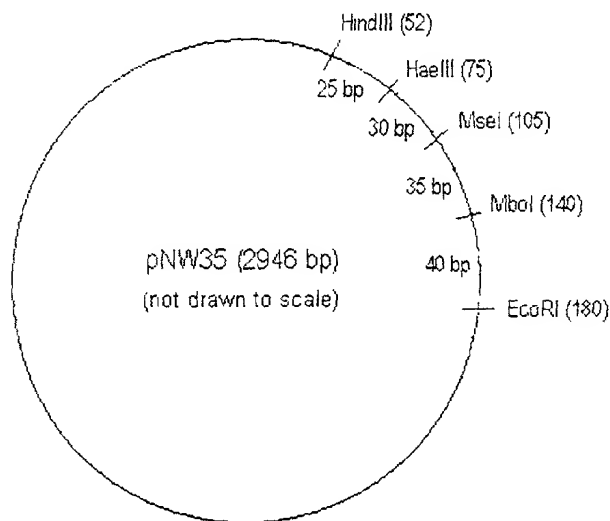

5

HindIII and EcoRI sites define the outer ends of the 25 bp and the 40 bp fragments. The sequence of pNW35 is shown below with the inserted region shown in bold type:

- 10 atgaccatgattacgccaagctctaatacgaactcactatagggaa**gcttccggacgtctcaggcta**atgttggcccacc
gacgttccacgatggggcgctctaagggttagaccctcgctcgggagtatttctgtgatctggcgacactcacgcg
agaagtcattaccggcgatatgaattcactggccgtcgtttacaacgtcgtgactgggaaaaccctggcgttacccaactta
atcgccctgcagcacatcccccttcgccagctggcgtaatagcgaagaggcccgaccgatcgccctcccaacagttgcg
cagcctgaatggcgaatgggaaattgtaaacgttaatatgtttaaattcgcgtaaatgtttaaatacagctcatttttaacca
15 ataggccgaaatcggcgaatcccttataaatcaaaagaatagaccgagataggggtgagtggttccagtttgaacaag
agtcactattaaagaacgtggactccaacgtcaaggcgaaaaaccgtctatcaggcgatggccactacgtgaacc
atcacctaatacaagttttggggtcgaggtgccgtaaagcactaaatcggaaccctaaggagccccgatttagagcttg
acggggaaagccggcgaacgtggcgagaaggaagggaagaaagcgaaaggagcgggctaggggcgctggcaa
gtgtagcgggtcacgtgcgcgttaaccaccacccgcgcgttaatgcgcgcgtacaggcgcgctcaggtggcacttttgcg
20 gggaaatgtgcggaacccctattgttttttctaatacattcaaatatgtatccgctcatgagacaataaccctgataaat
gcttcaataatattgaaaaaggaagagtatgagtattcaacattccgtgtcgccttattcccttttgcggcattttgccttctgttt
ttgctcaccagaaaacgtggtgaaagtaaaagatgctgaagatcagttgggtgcacgagtggttacatcgaactggatctc
aacagcggtaagatccttgagagtttgcggccgaagaacgtttccaatgatgagcacttttaaagttctgctatgtggcgcggt

attatcccgattgacgccgggcaagagcaactcggcgccgcatacactattctcagaatgacttggttgagtactcaccagtc
acagaaaagcatcttacggatggcatgacagtaagagaattatgcagtgctgccataaccatgagtgataaacactgcggcc
aacttactctgacaacgatcggaggaccgaaggagctaaccgctttttgcacaacatgggggatcatgtaactgccttgat
cgttgggaaccggagctgaatgaagccataccaaacgacgagcgtgacaccacgatgcctgtagcaatggcaacaacggt
5 ggcgaaactattaactggcgaactacttactctagcttcccggaacaattaatagactggatggagcgggataaagtgcag
gaccacttctcgctcggcccttccggctggctggtttattgctgataaatctggagccggtgagcgtgggtctcgcggtatcatt
gcagcactggggccagatggaagccctcccgatcgtagtattctacacgacggggagtcaggcaactatggatgaacga
aatagacagatcgctgagataggtgcctcactgattaagcattggtaactgtcagaccaagtttactcatatatactttgattgat
ttaaacttcattttaatttaaaggatctaggtaagatccttttgataaactcatgacaaaatcccttaacgtgagtttctgctcc
10 actgagcgcagaccccgtagaaaagatcaaaggatcttcttgagatcctttttctgcgcgtaactctgctgcttgcacaacaaa
aaaccaccgctaccagcgggtggtttgtttgccggatcaagagctaccaactcttttccgaaggtaactggcttcagcagagcg
cagataccaaatactgcttctagtgtagccgtagttagggccaccactcaagaactctgtagcaccgcctacatacctcgctc
tgctaactctgttaccagtggtcgtgccagtgccgataagtcgtgtcttaccgggttgactcaagacgatagttaccggataa
ggcgcagcggctgggctgaacgggggttcgtgcacacagccagcttgagcgaacgacctacaccgaactgagatac
15 ctacacgctgagctatgagaaagccacgcttcccgaggagaaaggcggacaggtatccggttaagcggcagggctcg
gaacaggagagcgcacgaggagcttccagggggaacgcctggtatctttagtctcgtcgggttcgccacctctgactt
gagcgtcgattttgtgatgctcgtcagggggcgagcctatggaaaaacgccagcaacgcggccttttaccggttctggcc
tttctgctgctttgtcgcacatgttcttctcgttatccctgattctgttgataaccgtattaccgctttgagtgagctgataccgc
tcgccgcagccgaacgaccgagcgcagcagtgagcaggaagcgggaagagcggccaataacgcaaaccgcctc
20 tccccgcgctgtggccgattcattaatgcagctggcacgacaggttcccgactggaaagcgggcagtgagcgaacgcaat
taatgtgagtgactcactcattaggcaccacaggctttacactttatgcttccggctcgtatgtgtgtggaattgtgagcgggataa
caatttcacacaggaaacagct

The 4 bp internal control plasmid for restriction enzymes


25 HaeIII, MboI, and MseI was prepared by the insertion of a synthetic 130 bp
fragment into HindIII / EcoRI-digested pMOSblue (Amersham Pharmacia
Biotech).

The insert region of 130 bp was sequenced twice in each
direction in order to establish the plasmid sequence. The presence of the
30 restriction sites and the mobility of the fragments released were also
checked by restriction digestion and polyacrylamide gel electrophoresis.

Preparation of the internal control spike DNA from the plasmid

The following illustration describes the strategy used for the
35 preparation of EcoRI-released internal control DNA:

- 123 -

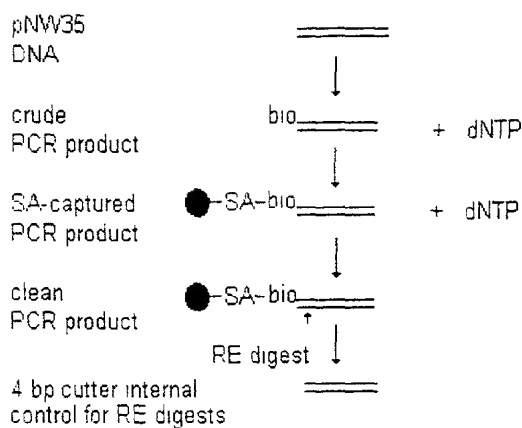

Primers and PCR

- 5 U-19 mer bio primer 5' bio-GTTTTCCCAGTCACGACGT 3'
- ICPCR(F) primer 5' TCCGGACGTCTCAGGCTAATGTT 3'

Bulk PCR amplification (96x 100 µl reactions, repeated to
 10 give 192 reactions in total) was carried out according to the conditions
 described below (all volumes are in µl):

	x1 (µl)	x100 (µl)	Final conditions
Water	84.4	8440	-
10x PCR buffer	10	1000	1x
25 mM dNTP mix	0.8	80	200 µM
20 µM T7 promoter-bio primer	1	100	20 pmol per PCR
20 µM U19 mer-bio primer	1	100	20 pmol per PCR
1 ng/µl pNW35	2	200	2 ng per PCR
5 U/µl <i>Taq</i> polymerase	0.8	80	4U per PCR
Total	100	10000	

The PCR master mix was rapidly dispensed into 96 PCR tubes. Thermal cycling was initiated using the following parameters: 94°C for 2 min; 50 °C for 2 min; 29 cycles of 72°C for 2 min, 94°C for 45 sec, and 50°C for 1 min; 72°C for 8 min; and then 4°C.

5

Capture of biotinylated PCR product termini and *EcoRI*- release of internal control DNA

All separations were carried out using a Dynal MPC-1 separator (Dynal, product #12001).

10 20 ml of pooled PCR reaction were mixed with 20 ml of Dynabeads M-280 in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M NaCl. The tube was incubated at room temperature for 1 hour with mixing on a Denley Spiramix 5.

15 The Dynabeads were then washed four times with 20 ml of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M NaCl. A fifth wash was performed in 20 ml of 1x buffer M.

20 EcoRI-digestion was carried out for 1 hour at 37°C in 5 ml of 0.25 U/μl EcoRI in 1x buffer M. The digest was then divided into ten 500 μl aliquot parts. Each aliquot was ethanol precipitated by the addition of 1 μl of See DNA, 0.1 volume of 3 M sodium acetate (pH 5.2), and 2.5 volumes of ethanol. The precipitations were mixed and chilled to 0°C on ice for 30 minutes and then centrifugation at 20,000 maxRCF for 10 minutes. The pellets were rinsed with 70 % ethanol before dissolving in a total of 500 μl of 1x TE buffer.

25 A 1 μl aliquot of the EcoRI-released internal control DNA was electrophoresed on an 8 % polyacrylamide gel to confirm that the size of the purified DNA was in accordance with that expected.

05914504-013102

- 125 -

High molecular weight genomic DNA digestion in the presence of the 4 bp cutter internal control DNA with a dilution series of the enzymes *HaeIII*, *MseI*, and *MboI*

High molecular weight human placental DNA was mixed with a dilution series of a mixture of the three set 1 4 bp cutter restriction endonucleases – each at the same number of units. Aliquots were then removed and mixed with the EcoRI-released internal control DNA. 3.6 µg of Human placental DNA (Sigma) was digested with 0.5 U/µl, 0.1 U/µl, 0.02 U/µl, 0.004 U/µl, and 0 U/µl each of *HaeIII*, *MboI*, and *MseI* in a total volume of 36 µl. 0.4 µg of placental DNA and 1 µl of EcoRI-released internal control DNA were digested with 0.5 U/µl, 0.1 U/µl, 0.02 U/µl, 0.004 U/µl, and 0 U/µl each of *HaeIII*, *MboI*, and *MseI* in a total volume of 4 µl. Vistra Green (Amersham Pharmacia Biotech) staining was used to monitor the extent of internal control DNA digestion and high molecular weight DNA digestion.

A premix of restriction enzymes (5 U/µl each enzyme), buffer, and BSA was prepared as described below:

Component	µl
50 U/µl <i>HaeIII</i>	1
20 U/µl <i>MseI</i>	2.5
25 U/µl <i>MboI</i>	2
10x NEB buffer #2	1
10x BSA	1
Water	2.5
Total	10

Serial 5-fold dilutions of the premix were prepared in 1x NEB buffer #2 and 1x BSA.

A 6x premix of human placental DNA, buffer, and BSA were also prepared as described below:

- 126 -

Component	6x mix (μ l)	per reaction (μ l)
1mg/ml human placental DNA	24	4
10x NEB buffer #2	9.6	1.6
10x BSA	9.6	1.6
Water	52.8	8.8

This premix contained placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA.

Human placental DNA and restriction mixes were then set up
5 as follows:

Tube	Component #1	Component #2
1	16 μ l of placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	4 μ l of each restriction enzyme at 5 U/ μ l in 1x NEB buffer #2 and 1x BSA
2	16 μ l of placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	4 μ l of each restriction enzyme at 1 U/ μ l in 1x NEB buffer #2 and 1x BSA
3	16 μ l of placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	4 μ l of each restriction enzyme at 0.2 U/ μ l in 1x NEB buffer #2 and 1x BSA
4	16 μ l of placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	4 μ l of each restriction enzyme at 0.04 U/ μ l in 1x NEB buffer #2 and 1x BSA
5	16 μ l of placental DNA at 0.25 mg/ml in 1x NEB buffer #2 and 1x BSA	4 μ l of 1x NEB buffer #2 and 1x BSA

2 μ l aliquots were then removed from tubes 1-5 and added to 1 μ l of EcoRI-released internal control DNA and 1 μ l of 2x NEB buffer #2
10 and 2x BSA to give samples 1-5ic. Samples 1-5ic were then overlaid with 50 μ l of mineral oil in order to prevent evaporation.

The remaining volume from tubes 1-5 was then added to 18 μ l of 1x NEB buffer #2 and 1x BSA.

All samples were finally incubated at 37°C overnight.

15

- 127 -

Restriction digests - final conditions

Restriction digests 1-5 were therefore performed under the following conditions:

tube	µg human placental DNA	U <i>HaeIII</i>	U <i>MseI</i>	U <i>MboI</i>	NEB buffer #2	µg/ml BSA
1	3.6	20	20	20	1x	100
2	3.6	4	4	4	1x	100
3	3.6	0.8	0.8	0.8	1x	100
4	3.6	0.16	0.16	0.16	1x	100
5	3.6	0	0	0	1x	100

Total volume = 36 µl

Likewise, restriction digests 1-5ic were performed under the following conditions:

tube	µg human placental DNA	µl IC DNA	U <i>HaeIII</i>	U <i>MseI</i>	U <i>MboI</i>	NEB buffer #2	µg/ml BSA
1	0.4	1	2	2	2	1x	100
2	0.4	1	0.4	0.4	0.4	1x	100
3	0.4	1	0.08	0.08	0.08	1x	100
4	0.4	1	0.016	0.016	0.016	1x	100
5	0.4	1	0	0	0	1x	100

IC = *EcoRI*-released internal control DNA

Total volume = 4 µl

10 µl of digests 1-5 were each mixed with 3 µl of loading dye and 4 µl of digests 1-5ic were mixed with 1 µl of loading dye. To sample number 5ic, 180 ng of PCR molecular weight markers were added to serve as size standards. The band sizes for these markers are 50 bp, 100 bp, 200 bp, 300 bp, 400 bp, 500 bp, 525 bp, 700 bp, and 1000 bp. The digests

- 128 -

were then electrophoresed on an 8 % polyacrylamide gel in 1x TBE as described below:

10 ml of 10x TBE

5 19.5 ml of 40 % acrylamide

19 ml of 2 % methylene bisacrylamide

50.4 ml of water

1 ml of freshly prepared 10 % (w/v) ammonium persulphate

100 µl of TEMED

10

A Cambridge Electrophoresis Ltd. vertical protein electrophoresis unit was used with 1 mm plate spacing. The samples were electrophoresed at 30 mA for 2 hours in 1x TBE. The gel was then stained for 30 minutes in 500 ml of 1x TBE containing 50 µl of Vistra Green. The
15 stained gel was imaged on a Fluorimager with the following settings: a 488 nm laser; a 570 DF 30 filter; a PMT setting of 700 V; 200 µm resolution; and low sensitivity.

The results are shown in figure 3.


20 **Example 2**

TRSPA-2 analysis with the pNW33 *Bam*HI, *Hind*III, and *Afl*II matrix and the pNW33 *Hind*III, *Nco*I, and *Spe*I matrix probed with the 140 bp *Bsp*EI to *Bam*HI fragment from pNW33

25

As an example of the case where the probe hybridizes to an arrayed PCR fragment with a different restriction site at each end, the pNW33 *Bam*HI, *Hind*III, and *Afl*II matrix (matrix #7) was probed with a PCR product from the 140 bp *Bsp*EI (466) to *Bam*HI (605) restriction fragment within pNW33. In this example, the probe binds to a 204 bp *Hind*III to
30 *Bam*HI PCR product derived from the 160 bp *Hind*III (445) to *Bam*HI (605) restriction fragment within pNW33 (see below).

- 129 -

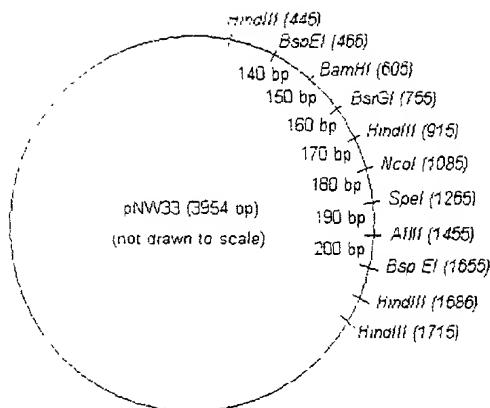

As an example of the case where the probe hybridizes to an
 5 arrayed PCR fragment with the same restriction site at each end, the
 pNW33 HindIII, NcoI, and SpeI matrix (matrix #17) was probed with a PCR
 product from the 140 bp BspEI (466) to BamHI (605) restriction fragment
 within pNW33. In this example, the probe binds to a 514 bp HindIII to
 HindIII PCR product derived from the 470 bp HindIII (445) to HindIII (915)
 10 restriction fragment within pNW33 (see above).

Oligonucleotides

<i>BamHI</i> short PCR primer	5' TGTAACGACACATTGCTGGATACC 3'
15 <i>HindIII</i> short PCR primer	5' ATATAACTCTCGCTCCTTGATAAC 3'
<i>NcoI</i> short PCR primer	5' AGGCGTCTGAGGCTGCGGCTATGG 3'
<i>SpeI</i> short PCR primer	5' AACCCGTCGCGACGAGAGTCTAAG 3'
<i>AflIII</i> short PCR primer	5' GATATACGTGATATATTTTGATTG 3'
20 <i>BamHI</i> adaptor	5' pGATCGGTATCCAGCAATGTGTCGTTACA 3'
<i>HindIII</i> adaptor	5' pAGCTGTTATCAAGGAGCGAGAGTTATAT 3'
<i>NcoI</i> adaptor	5' pCATGCCATAGCCGCAGCCTCAGACGCCT 3'
<i>SpeI</i> adaptor	5' pCTAGCTTAGACTCTCGTCGCGACGGGTT 3'

- 130 -

*Afl*III adaptor 5' pTTAACAATCAAAATATATCACGTATATC 3'

*Bam*HI long PCR primer 5' TGTAACGACACATTGCTGGATACCGATCC 3'

*Hind*III long PCR primer 5' ATATAACTCTCGCTCCTTGATAACAGCTT 3'

5 *Nco*I long PCR primer 5' AGGCGTCTGAGGCTGCGGCTATGGCATGG 3'

*Spe*I long PCR primer 5' AACCCGTCGCGACGAGAGTCTAAGCTAGT 3'

*Afl*III long PCR primer 5' GATATACGTGATATATTTTGATTGTTAAG 3'

Luc140down primer

10 5' GCGCTAGGGATCCTTACTGGGACGAAGACGAA 3'

Luc140up-bio primer

5' biotin-CGCAGCTGGTAATCCGGACGCCCCGCGTCGAAGATGTT3'

15 Restriction digestion of pNW33

Restriction digests of pNW33 were set up by combining the following:

Matrix #7

20

450 µl	22 µg/ml pNW33 in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Bam</i> HI in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Hind</i> III in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Afl</i> III in 1x NEB buffer #2 + 100 µg/ml BSA

Matrix #17

450 µl	22 µg/ml pNW33 in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Hind</i> III in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Nco</i> I in 1x NEB buffer #2 + 100 µg/ml BSA
50 µl	0.6 U/µl <i>Spe</i> I in 1x NEB buffer #2 + 100 µg/ml BSA

- 131 -

The samples were then vortex mixed, and incubated at 37°C overnight.

Calf intestinal alkaline phosphatase (CIAP) digestion of pNW33

5 restriction digests

400 µl fractions of the 20 restriction digests (each containing 6.6 µg of digested pNW33) were ethanol precipitated by the addition of 1 µl of See DNA, 0.1 volume of 3 M sodium acetate (pH 5.2), and 2.5 volumes of ethanol, chilling to -20°C, and then centrifugation. Pellets were rinsed
10 with 70 % ethanol prior to re-dissolution in 20 µl of 1x CIAP buffer containing 40 U of CIAP. The CIAP digests were carried out for 5 hours at 37°C and were then made up to 400 µl with 1x TE buffer.

Phenol extraction of pNW33 CIAP digests

15 The diluted digests were extracted with 400 µl of phenol and then ethanol precipitated as described above but with a 100 % ethanol wash after the 70 % ethanol wash. Samples were finally re-dissolved in 20 µl of TE buffer.

20 Annealing of short PCR primers to cognate adaptors

Short PCR primers and their cognate adaptors were annealed by adding 1 µl of 200 µM short PCR primer to 1 µl of 200 µM cognate adaptor in 20 µl of 50 mM NaCl, 1x TE buffer. The mixed oligonucleotides were overlaid with 30 µl of light mineral oil and were then heated to 90°C
25 for 5 minutes followed by slow cooling to room temperature. The annealed short PCR primer / cognate adaptor complexes were then diluted with 1 ml of 1x TE buffer and stored frozen at -20°C.

Ligation to annealed short PCR primers and cognate adaptors

30 1 µl of phenol extracted pNW33 digest was used per ligation. Ligations were performed in 100 µl of 1x ligase buffer containing 1 mM ATP

- 132 -

and 10 µl (100 U) of T4 DNA ligase. 20 µl aliquots from the 1 ml of annealed short PCR primer / cognate adaptor complexes were added according to the following table.

#7	<i>Bam</i> HI / <i>Hind</i> III / <i>Afl</i> II
#17	<i>Hind</i> III / <i>Nco</i> I / <i>Spe</i> I

5

Ligation reactions were carried out for 24 hours at 16°C.

Samples were then diluted to 1 ml with TE buffer and stored at -20°C.

The diluted ligation reactions were then further diluted 1 in 10 and 10 µl was used as PCR template per 100 µl reaction.

10

Restriction digestion of human placental DNA

1 U (~50 µg) of human placental DNA was digested overnight at 37°C with 100 U each of *Bam*HI, *Bsr*GI, *Hind*III, *Nco*I, *Spe*I, and *Afl*II in 400 µl of 1x NEB buffer #2 containing 100 µg/ml BSA.

15

CIAP digestion of human placental DNA restriction digest

The digest was ethanol precipitated by the addition of 1 µl of See DNA, 0.1 volume of 3 M sodium acetate (pH 5.2), and 2.5 volumes of ethanol, chilling to -20°C, and then centrifugation. The pellet was rinsed with 70 % ethanol prior to re-dissolution in 50 µl of 1x CIAP buffer containing 100 U of CIAP. The CIAP digest was carried out for 5 hours at 37°C and was then made up to 400 µl with 1x TE buffer.

25

Phenol extraction of human placental DNA CIAP digest

The diluted CIAP digest was extracted with 400 µl of phenol and then ethanol precipitated as described above – again with a 100 % ethanol wash after the 70 % ethanol wash. The sample was finally re-dissolved in 10 µl of TE buffer.

Ligation to annealed short PCR primers and cognate adaptors

Short PCR primers were annealed to their cognate adaptors as described above.

The ligation to annealed short PCR primers and cognate adaptors was carried out in 100 μ l of 1x ligase buffer with 1 mM ATP, 100 U of T4 DNA ligase and 10 μ l of each of the six short PCR primer / cognate adaptor complexes as above.

The ligation reaction was carried out for 24 hours at 16°C. The sample was then diluted to 1 ml with 1x TE buffer and stored at -20°C. 0.2 μ l was used as PCR template per 100 μ l reaction.

PCR amplification conditions

An initial touch-down reaction was carried out in 50 μ l of 1x PCR buffer with all four dNTPs at 200 μ M and *Taq* DNA polymerase at 0.05 U/ μ l. Long PCR primers were used at 400 nM. 10 μ l of pNW33 PCR template was used per reaction and 0.2 μ l of human placental DNA PCR template was used per reaction. The samples were overlaid with 40 μ l of light mineral oil and were touch-down thermocycled as described below:

98°C for 1 min, 72°C for 2 min, 72°C for 5 min
98°C for 1 min, 69°C for 2 min, 72°C for 5 min
98°C for 1 min, 66°C for 2 min, 72°C for 5 min
98°C for 1 min, 63°C for 2 min, 72°C for 5 min
98°C for 1 min, 60°C for 2 min, 72°C for 5 min

25

The main thermal cycling was then carried out in 100 μ l of 1x PCR buffer with all four dNTPs at 200 μ M and *Taq* DNA polymerase at 0.05 U/ μ l. Long PCR primers were used at 400 nM. The samples were subjected to 20 cycles of 98°C for 1 min, 60°C for 2 min, 72°C for 5 min and then 72°C for 10 min followed by chilling at 10°C and recovery from under oil.

30

Arraying onto nylon membranes

1 μ l aliquots from the PCR amplified samples were spotted onto Hybond N+ nylon membranes. The membranes were then transferred to a stack of three sheets of 3MM paper, saturated with 0.4 M NaOH, and incubated for 10 minutes. The NaOH was rinsed away in 2x SSC and the membranes were used directly for the pre-hybridization.

Probe synthesis

10 A PCR master mix was prepared as follows (all volumes are in μ l):

Master Mix Preparation	140
Water	183.6
10x PCR buffer	30
2 mM dATP	15
2 mM dGTP	15
2 mM dTTP	15
2 mM dCTP	3
140 probe template	1
Luc140up-bio	2.5
Luc140down	2.5
Taq DNA polymerase (5 U/ μ l)	2.4
Total	270

15 Five PCR reactions were performed. For each reaction, 5 μ l of ³³P-labelled dCTP was added to the PCR tube. The reactions were then made up to 50 μ l by the addition of 45 μ l of the master mix to each tube. Each reaction was gently mixed.

PCR cycling parameters

94°C	2 min		
50 °C	2 min		
5	72°C	2 min	
	94°C	45 sec	} 25 cycles
	50 °C	1 min	
	72°C	8 min	
	4°C	hold	

10 After thermal cycling, the five PCR amplifications were pooled. The 250 µl of labelled DNA was then used for the following magnetic bead purification procedure.

15 **Capture of biotinylated PCR product and release of non-biotinylated strands**

 All separations were carried out using a Dynal MPC-4 separator (Dynal, product #12004).

 250 µl of the pooled PCR was mixed with an equal volume of 10 mg/ml streptavidin-coated colloidal Fe₃O₄ particles in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M NaCl. The tube was incubated at room
20 temperature for 1 hour with mixing on a Denley Orbital Mixer.

 The streptavidin-coated colloidal Fe₃O₄ particles were then washed with 1 ml of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M NaCl. Three more identical washes were performed.

25 The washed streptavidin-coated colloidal Fe₃O₄ particles were incubated in 500 µl of 0.1 M NaOH for 10 minutes at room temperature. The supernatant was removed and added to 500 µl of 0.5 M

- 136 -

HEPES. The samples were then ethanol precipitated by the addition of 0.1 volume of 3 M sodium acetate (pH 5.2) and 2.5 volumes of ethanol, chilling to 0°C (on ice for 30 minutes), and then centrifugation at 20,000 maxRCF for 10 minutes. The pellets were rinsed with 70 % ethanol before
5 dissolving in 100 µl of TE buffer.

Hybridization

Each membrane was placed in a 55 mm x 35 mm x 21 mm plastic box and 1.25 ml of pre-hybridization solution (5x SSC; Denhardt's
10 solution; 1 % SDS; 10 % dextran sulphate [Mw 500,000]; 0.3 % tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA – pre-warmed to 65°C) was added. Each box was closed and incubated at 65°C for 50 minutes on a rocking platform. The pre-hybridization solution was removed and replaced with hybridization solution (5x SSC; Denhardt's
15 solution; 1 % SDS; 10 % dextran sulphate [Mw 500,000]; 0.3 % tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA – containing 5 µl of the appropriate 33P-labelled probe) and the box was incubated at 65°C for 3 hours on a rocking platform.

20 Washing

The membranes were drained and transferred to 200 ml of 2x SSC, 0.1 % SDS at 68°C for 30 minutes. A further wash was carried out in 0.2x SSC, 0.1 % SDS at 71°C for 30 minutes. The membranes were rinsed in 2x SSC at room temperature and laid out on blotting paper to
25 remove excess liquid. Once dry, the membrane was covered in Saran Wrap and exposed to a Kodak Phosphor Screen for 1 hour. The phosphor screen was subsequently imaged using a Molecular Dynamics Storm 860 Phosphorimager.

Results are shown in figures 4 and 5.

30

Example #3a

A single cycle of inter-population perfectly matched duplex depletion wherein *E.coli* MutS protein is used to capture an A/A mismatch-
5 containing duplex

'Affected' versus 'unaffected' (i.e. inter-population) mismatch-containing duplex selection can be achieved by: attaching a mismatch-binding protein to a solid support (or using the mismatch-binding protein in solution followed by subsequent solid-phase capture); taking denatured
10 'affected' DNA fragments and hybridizing these to denatured and biotinylated 'unaffected' DNA fragments; and capture of mismatch-containing duplex molecules with the mismatch-binding protein. Releasing the mismatch-containing duplex molecules (without strand denaturation), streptavidin capture and then release of the non-biotinylated strands will
15 give only the desired species.

In this example, PCR fragments are prepared and used to demonstrate each of the individual steps for a single cycle of inter-population perfectly matched duplex depletion using *E.coli* MutS protein.

20 Clone insert design

The clone inserts were constructed using standard cloning methodology well known to those skilled in the art and were inserted between the *AvaI* site and *EcoRI* of pMOSBlue (Amersham Pharmacia Biotech).

25 The clone inserts contain a common 9 base pair internal core sequence in which a single nucleotide change or an insertion can be located. The internal core sequence is derived from codons 272-274 of human p53. These codons (GTG CGT GGT) correspond to a mutational hotspot found in lung and other types of cancer (R273L). For the design of
30 the clone inserts containing a mismatch (only one from a complete series of which is shown below), the nucleotide in position 5 of this core sequence is

- 138 -

modified (GTGCXTGGT).

The internal core sequence is flanked by a random sequence
- allowing the independent detection of the clone #1 and the clone #7
insert sequence in a mixed population of clone inserts.

5

Clone #1**Mutant sequence (#1M)**

10 5' CCCGGGGGATCCTCGTTTTATTGGGCCGAGTTTGGTCCGTAGTGCTTGGTTAGATATGCTTAT
3' GGGCCCCCTAGGAGCAAAATAACCCGGCTCAAACCAGGCATCACGAACCATCTATACGAATA

GTTCACAAAATCATCCTTGACAGAATTC3'
CAAGTGTTTTAGTAGGAACATGCTTAAG5'

15

Control sequence (#1C)

5' CCCGGGGGATCCTCGTTTTATTGGGCCGAGTTTGGTCCGTAGTGCATGGTTAGATATGCTTAT
20 3' GGGCCCCCTAGGAGCAAAATAACCCGGCTCAAACCAGGCATCACGTACCAATCTATACGAATA

GTTCACAAAATCATCCTTGACAGAATTC3'
CAAGTGTTTTAGTAGGAACATGCTTAAG5'

Clone #7

25

Control sequence (#7C)

5' CCCGGGTGTACACAAAAGTTTACCTGAAGAACGTGGGGGTCGTGCCTGGTCTTGCGTCACCTG
30 3' GGGCCACATGTGTTTTCAAATGGACTTCTTGACCCCCCAGCACGGACCAGAACGCAGTGGAC

GTCTCAGGAGAGGGTCCCCATGGGAATTC3'
CAGAGTCCTCTCCAGGGGTACCCCTTAAG5'

35

Preparation of upper strand 5'-biotinylated and lower strand 5'-biotinylated #1C double-stranded PCR product, and upper strand 5'-biotinylated and lower strand 5'-biotinylated #7C double-stranded PCR product

5

Oligonucleotides

BIOUPST2 5' bio-CTACTGATCGGATCCCCG 3'

BIODOWN3 5' bio-AAACGACGGCCAGTGAAT 3'

10 PCR reaction set-up (#1C)

#1C in pMOSBlue at 2.50 µg/ml	100 µl
Water	8400 µl
10x PCR buffer	1000 µl
15 50 µM BIOUPST2	100 µl
50 µM BIODOWN3	100 µl
10 mM dNTPs	200µl
Taq DNA polymerase (5 U/µl)	100 µl

20 PCR reaction set-up (#7C)

#7C in pMOSBlue at 3.32 µg/ml	100 µl
Water	8400 µl
10x PCR buffer	1000 µl
25 50 µM BIOUPST2	100 µl
50 µM BIODOWN3	100 µl
10 mM dNTPs	200µl
Taq DNA polymerase (5 U/µl)	100 µl

30

- 140 -

96x 200 µl reactions were carried out for template #1C and
96x 200 µl reactions were carried out for template #7C on a 96-well Perkin
Elmer Cetus GeneAmp PCR System 9600 machine as described below:

95 °C, 5 minutes	1 cycle
95 °C, 1 minute	30 cycles
50 °C, 1 minute	
72 °C, 1 minute	
72 °C, 5 minutes	1 cycle
4°C, hold	

5

The PCR products were pooled together and precipitated by
adding 0.1 volumes of 3 M sodium acetate and 1 volume of isopropanol
followed by centrifugation at 16,000 rpm for 30 minutes at 4°C (in a
Centrikon T-2070 ultracentrifuge; swinging bucket Kontron rotor TST
41.14). Pellets were washed with 14 ml of ethanol and centrifuged at
20,000 rpm for 30 minutes. Finally, the pellets were air-dried and
resuspended in a total volume of 0.6 ml of TE buffer.

The 0.6 ml PCR sample was purified in twelve Microspin S-
300 HR columns (50 µl per column) following the manufacturer's protocol.
Briefly, the resin in the columns was resuspended by vortexing. Columns
were centrifuged at 735 x g (3000 rpm in a microfuge) for 1 minute. The
sample was then applied to the centre of the resin, being careful not to
disturb the bed. The columns were centrifuged at 735 x g for 2 minutes
and the flow-through containing the PCR product was collected. The
twelve eluted 50 µl volumes were pooled together (pool 1). Columns were
washed with 50 µl of TE buffer and the eluted fractions were loaded onto a
fresh S-300 HR column. Product yield and removal efficiency of the PCR
primers were analysed on a 1.5 % agarose gel.

25

Preparation of non-biotinylated #1M single-stranded DNA and non-biotinylated #7C single-stranded DNA

Oligonucleotides

5

BIOUPST2 5' bio-CTACTGATCGGATCCCCG 3'
 DOWN3 5' AAACGACGGCCAGTGAAT 3'

PCR reaction set-up (#1M)

10

#1M in pMOS <i>Blue</i> at 3.56 µg/ml	100 µl
Water	8400 µl
10x PCR buffer	1000 µl
50 µM BIOUPST2	100 µl
15 50 µM DOWN3	100 µl
10 mM dNTPs	200µl
Taq DNA polymerase (5 U/µl)	100 µl

PCR reaction set-up (#7C)

20

#7C in pMOS <i>Blue</i> at 3.32 µg/ml	100 µl
Water	8400 µl
10x PCR buffer	1000 µl
50 µM BIOUPST2	100 µl
25 50 µM DOWN3	100 µl
10 mM dNTPs	200µl
Taq DNA polymerase (5 U/µl)	100 µl

48x 200 µl reactions were carried out for template #1M and

30 48x 200 µl reactions were carried out for template #7C on a 96-well Perkin

- 142 -

Elmer Cetus GeneAmp PCR System 9600 machine as described below:

95 °C, 5 minutes	1 cycle
95 °C, 1 minute	30 cycles
50 °C, 1 minute	
72 °C, 1 minute	
72 °C, 5 minutes	1 cycle
4°C, hold	

Capture of biotinylated PCR product strands and release of non-

- 5 **biotinylated strands for the preparation of non-biotinylated #1M
single-stranded DNA and non-biotinylated #7C single-stranded DNA**

10 10 ml of pooled #1M and #7C PCR amplifications were each
mixed with an equal volume of 4 mg/ml streptavidin-coated colloidal Fe_3O_4
particles taken up in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M
NaCl. The tubes were incubated at room temperature for 60 minutes with
mixing.

The streptavidin-coated colloidal Fe_3O_4 particles were then
washed with 20 ml of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M
NaCl. Two more identical washes were performed.

- 15 The washed streptavidin-coated colloidal Fe_3O_4 particles
were finally incubated in 800 μl of 0.1 M NaOH for 10 minutes at room
temperature. The supernatants were removed and added to 200 μl of 2 M
HEPES (free acid). Samples were quantified by absorbance at 260 nm.

- 20 **Denaturation and annealing**

Denaturation and annealing reactions were prepared by
mixing:

A = 400 ng of upper strand 5'-biotinylated and lower strand 5'-biotinylated #1C double-stranded PCR product
B = 400 ng of upper strand 5'-biotinylated and lower strand 5'-biotinylated #7C double-stranded PCR product
C = 200 ng of non-biotinylated #1M single-stranded DNA (the lower of the two #1M strands shown above)
D = 200 ng of non-biotinylated #7C single-stranded DNA (the lower of the two #7C strands shown above)

Reannealing between the upper strand of A and the single-stranded C will therefore give rise to an A/A mismatch-containing duplex, 5'-biotinylated on the upper strand, for clone insert #1.

5 Reannealing between the upper strand of B and the single-stranded D will therefore give rise to a perfectly matched duplex, 5'-biotinylated on the upper strand, for clone insert #7.

0.1 volumes of 1 M NaOH were then added, followed by incubation at room temperature for 10 minutes. 0.25 volumes of 2 M
10 HEPES (free acid) were finally added followed by incubation at 42°C for 1 hour.

Samples were adjusted to 50 µl and were made 1x in PBS and 1 mg/ml in BSA ready for reaction with MutS protein-coated magnetic beads.

15 One 50 µl sample (the pre-enrichment control, sample 6) was used directly for capture of biotinylated PCR product strands and release of non-biotinylated strands.

Mismatch capture

20 20 µl of M₂B₂ MutS protein-coated magnetic particles (Genecheck, lot #20) were added to the annealed DNA above. Samples were incubated for 1 hour at room temperature with shaking.

Samples were then washed twice with 200 µl of ice-cold PBS.

- 144 -

Samples were finally eluted from the magnetic beads for 10 minutes at room temperature in 50 μ l of the following:

Sample 1	1 M NaCl
Sample 2	PBS
Sample 3	1 M urea
Sample 4	1% (w/v) SDS
Sample 5	10 mM NaOH

5 **Capture of biotinylated PCR product strands and release of non-biotinylated strands**

All separations were carried out using an Amersham magnetic separator (RPN1682, batch #1).

50 μ l of the eluates from the MutS protein-coated magnetic
10 beads and the pre-enrichment control (sample 6) were each mixed with an equal volume of 4 mg/ml streptavidin-coated colloidal Fe_3O_4 particles in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M NaCl. The tubes were incubated at room temperature for 30 minutes with regular mixing.

The streptavidin-coated colloidal Fe_3O_4 particles were then
15 washed twice with 500 μ l of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M NaCl at room temperature.

The washed streptavidin-coated colloidal Fe_3O_4 particles were incubated in 10 μ l of 0.1 M NaOH for 10 minutes at room temperature. The supernatant was removed and added to 2.5 μ l of 2 M
20 HEPES (free acid).

5 μ l fractions were finally spotted onto Hybond N+ nylon membranes along with 200 ng, 20 ng, 2 ng, and 200 pg amounts of the following:

25

- 145 -

C = non-biotinylated #1M single-stranded DNA
(the lower of the two #1M strands shown previously)

D = non-biotinylated #7C single-stranded DNA
(the lower of the two #7C strands shown previously)

Probe labelling

Oligonucleotides

5

#1 probe oligo 5' GGCCGAGTTTTGGTCCGTAG 3'

#7 probe oligo 5' GTCTTGCGTCACCTGGTCTCAG 3'

10 Preparation of probes

#1 probe oligo and #7 probe oligo were radioactively 5' end-labelled using T4 polynucleotide kinase as described below (all volumes are in μ l):

	#1 probe	#7 probe
#1 probe oligo (25 pmol)	1.25	-
#7 probe oligo (25 pmol)	-	1.25
10x PNK buffer	2.5	2.5
[γ - ³³ P] ATP >92 TBq/mmol, >2500 Ci/mmol, 370 mBq/ml, 10 mCi/ml, (Amersham Pharmacia Biotech., AH9968, lot #B0006)	12.5	12.5
T4 PNK (10 U/ μ l) (Amersham International, E70031Y, lot #201226)	2.5	2.5
Water	6.25	6.25

15

The reactions were incubated at 37°C for 30 minutes and

09914604-013402

- 146 -

then heated to 70°C for 5 minutes to denature the enzyme.

MicroSpin G-25 column purification

Two G-25 columns (APB 27-5325-01, lot #9015325011) were
5 resuspended by vortexing, and the bottom closures were snapped off as
described in the manufacturer's instructions. A pre-spin of the columns
was carried out for 1 minute at 730 maxRCF – 2670 rpm in a Hettich
Zeutrifugen EBA 12 benchtop centrifuge – and the eluates were discarded.
The 25 µl reactions were applied to each column, and the columns were
10 centrifuged for 2 minutes at 730 maxRCF. The eluates from the second
spin were stored and used as probes.

Hybridization

Each membrane was placed in a 55 mm x 35 mm x 21 mm
15 plastic box and 2.5 ml of pre-hybridization solution (5x SSC; Denhardt's
solution; 1 % SDS; 10 % dextran sulphate [MW 500,000]; 0.3 %
tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA – pre-
warmed to 42°C) was added. Each box was closed and incubated at 42°C
for 1 hour on a rocking platform. The pre-hybridization solution was
20 removed and replaced with hybridization solution (5x SSC; Denhardt's
solution; 1 % SDS; 10 % dextran sulphate [MW 500,000]; 0.3 %
tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA –
containing 2.5 µl of the appropriate ³³P-labelled probe) and the box was
incubated at 42°C overnight on a rocking platform.

25

Washing



The membranes were drained and transferred to 200 ml of 2x
SSC, 0.1 % SDS at 42°C for 10 minutes. A further wash was carried out in
0.2x SSC, 0.1 % SDS at 42°C for 10 minutes. The membranes were
30 rinsed in 2x SSC at room temperature and laid out on blotting paper to
remove excess liquid. Once dry, the membrane was covered in Saran

5 Dot blot layout

Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Sample 6
5 µl fraction	5 µl fraction	5 µl fraction	5 µl fraction	5 µl fraction	5 µl fraction
200 ng #1M ssDNA	20 ng #1M ssDNA	2 ng #1M ssDNA	200 pg #1M ssDNA	-	-
200 ng #7C ssDNA	20 ng #7C ssDNA	2 ng #7C ssDNA	200 pg #7C ssDNA	-	-

ssDNA = single-stranded DNA

10 Probe hybridisation results

#1 probe	#7 probe
	

Signal intensities from the spots above were reported using ImageQuant 5.0 software (Molecular Dynamics), with the SumAboveBG

15 figures being used after drawing a 6x3 grid over the array of spots.

- 148 -

#1 probe SumAboveBG signals

Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Sample 6
142257	24356	180359	394296	465242	658323
200 ng	20 ng	2 ng	200 pg	-	-
#1M ssDNA	#1M ssDNA	#1M ssDNA	#1M ssDNA		
15961011	2842937	252392	25316	2281	2475
200 ng	20 ng	2 ng	200 pg	-	-
#7C ssDNA	#7C ssDNA	#7C ssDNA	#7C ssDNA		
66536	37955	10621	2736	1955	1934

#7 probe SumAboveBG signals

Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Sample 6
66392	17071	118109	153372	166211	666853
200 ng	20 ng	2 ng	200 pg	-	-
#1M ssDNA	#1M ssDNA	#1M ssDNA	#1M ssDNA		
80873	23256	6371	2399	1831	1857
200 ng	20 ng	2 ng	200 pg	-	-
#7C ssDNA	#7C ssDNA	#7C ssDNA	#7C ssDNA		
10619299	2052027	207827	15470	2055	1844


Recovery

Elution conditions	#1 signal (% input ssDNA)	#7 signal (% input ssDNA)	#1/#7 signal ratio
1 M NaCl	21.6 %	10.0 %	2.2
PBS	3.7 %	2.6 %	1.4
1 M urea	26.4 %	17.7 %	1.5
1 % SDS	59.9 %	23.0 %	2.6
10 mM NaOH	70.7 %	24.9 %	2.8

10

Recovery figures are plotted below for the various different elution conditions:

- 149 -

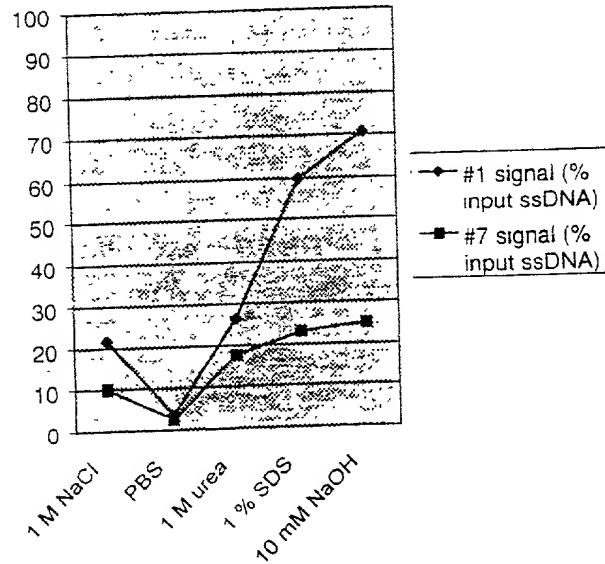

5 Example #3b

A single cycle of inter-population perfectly matched duplex depletion wherein bacteriophage T4 endonuclease VIII protein containing a cleavage-inactivating N62D point mutation is used to capture an A/A mismatch-containing duplex

10 In this example, PCR fragments are again prepared and used to demonstrate each of the individual steps for a single cycle of inter-population perfectly matched duplex depletion using bacteriophage T4 endonuclease VIII protein containing a cleavage-inactivating N62D point mutation.

Clone design and DNA preparation

Clone insert design was exactly as described in example #3a (see clone #1: mutant sequence (#1M), control sequence (#1C) and clone

20 #7: control sequence (#7C)).

- 150 -

Preparation of upper strand 5'-biotinylated and lower strand 5'-biotinylated #1C double-stranded PCR product, and upper strand 5'-biotinylated and lower strand 5'-biotinylated #7C double-stranded PCR product were as described in example #3a.

- 5 Preparation of non-biotinylated #1M single-stranded DNA and non-biotinylated #7C single-stranded DNA were again as described in example #3a.

Denaturation and annealing

- 10 Denaturation and annealing reactions were prepared as described in example #3a.

- Reannealing will therefore again give rise to an A/A mismatch-containing duplex, 5'-biotinylated on the upper strand, for clone insert #1 and a perfectly matched duplex, 5'-biotinylated on the upper strand, for clone insert #7.

15 One 50 µl sample (the pre-enrichment control) was adjusted to 100 µl with TE buffer and was used directly for capture of biotinylated PCR product strands and release of non-biotinylated strands.

20 Mismatch capture

- 50 µl samples of the annealing reaction were mixed with an equal volume of 200 mM sodium phosphate (pH 6.5), 100 mM KCl. 10 µg of GST-tagged T4 endonuclease VII N62D mutant (obtained from Prof. Börries Kemper, Univ. Cologne) were added to the annealed DNA and the samples were incubated for 15 minutes at 16°C.

- Samples were then mixed with 200 µl of a 50 % slurry of Glutathione Sepharose 4B (Amersham Pharmacia Biotech, lot #279991) in 100 mM sodium phosphate (pH 6.5), 50 mM KCl. The tubes were incubated at 16°C for 30 minutes with regular mixing. The mixture was then transferred to a spin column to separate solid from liquid phases.

Samples were finally eluted from the Glutathione Sepharose

- 151 -

4B matrix for 10 minutes at room temperature in 100 μ l of 10 mM reduced glutathione in 50 mM tris-HCl (pH 8.0).

Capture of biotinylated PCR product strands and release of non-biotinylated strands

All separations were carried out using an Amersham magnetic separator (RPN1682, batch #1).

100 μ l of the eluate from the N62D T4 endonuclease VII mismatch-capture reaction, and the pre-enrichment control were each mixed with an equal volume of 4 mg/ml streptavidin-coated colloidal Fe_3O_4 particles in 20 mM tris-HCl (pH 7.4), 2 mM EDTA (pH 8.0), 2 M NaCl. The tubes were incubated at room temperature for 30 minutes with regular mixing.

The streptavidin-coated colloidal Fe_3O_4 particles were then washed twice with 500 μ l of 10 mM tris-HCl (pH 7.4), 1 mM EDTA (pH 8.0), 1 M NaCl at room temperature.

The washed streptavidin-coated colloidal Fe_3O_4 particles were incubated in 10 μ l of 0.1 M NaOH for 10 minutes at room temperature. The supernatant was removed and added to 2.5 μ l of 2 M HEPES (free acid).

5 μ l fractions were finally spotted onto Hybond N+ nylon membranes along with 200 ng, 20 ng, 2 ng, and 200 pg amounts of the following:

C = non-biotinylated #1M single-stranded DNA (the lower of the two #1M strands shown previously)
D = non-biotinylated #7C single-stranded DNA (the lower of the two #7C strands shown previously)

25

Probe labelling

Probe labelling was exactly as described in example #3a.

Hybridization

Each membrane was placed in a 55 mm x 35 mm x 21 mm plastic box and 2.5 ml of pre-hybridization solution (5x SSC; Denhardt's solution; 1 % SDS; 10 % dextran sulphate [MW 500,000]; 0.3 % tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA – pre-warmed to 42°C) was added. Each box was closed and incubated at 42°C for 10 minutes on a rocking platform. The pre-hybridization solution was removed and replaced with hybridization solution (5x SSC; Denhardt's solution; 1 % SDS; 10 % dextran sulphate [MW 500,000]; 0.3 % tetrasodium pyrophosphate; 100 µg/ml denatured, sonicated DNA – containing 2.5 µl of the appropriate ³³P-labelled probe) and the box was incubated at 42°C overnight on a rocking platform.

Washing

The membranes were drained and transferred to 200 ml of 2x SSC, 0.1 % SDS at 42°C for 10 minutes. A further wash was carried out in 0.2x SSC, 0.1 % SDS at 42°C for 10 minutes. The membranes were rinsed in 2x SSC at room temperature and laid out on blotting paper to remove excess liquid. Once dry, the membrane was covered in Saran Wrap and exposed to a Kodak Phosphor Screen for 1 hour. The phosphor screen was subsequently imaged using a Molecular Dynamics Storm 860 Phosphorimager.

Signal intensities from the spots above were again reported using ImageQuant 5.0 software (Molecular Dynamics), with the SumAboveBG figures being used after drawing a grid over the array of spots.

- 153 -

#1 probe SumAboveBG signals

200 ng #1M ssDNA 22418970	20 ng #1M ssDNA 3317389	2 ng #1M ssDNA 270664	200 pg #1M ssDNA 25763	GSH eluted sample 144398
200 ng #7C ssDNA 229951	20 ng #7C ssDNA 105343	2 ng #7C ssDNA 16354	200 pg #7C ssDNA 3632	Pre- enrichment control 997077

#7 probe SumAboveBG signals


200 ng #1M ssDNA 117267	20 ng #1M ssDNA 42226	2 ng #1M ssDNA 9421	200 pg #1M ssDNA 2812	GSH eluted sample 113608
200 ng #7C ssDNA 3665967	20 ng #7C ssDNA 22248096	2 ng #7C ssDNA 203869	200 pg #7C ssDNA 19406	Pre- enrichment control 1135786

Recovery

Elution conditions	#1 signal (% input ssDNA)	#7 signal (% input ssDNA)	#1/#7 signal ratio
10 mM reduced glutathione in 50 mM tris-HCl (pH 8.0)	14.5 %	10.0 %	1.45

10 The enrichment results are presented graphically below:

- 154 -

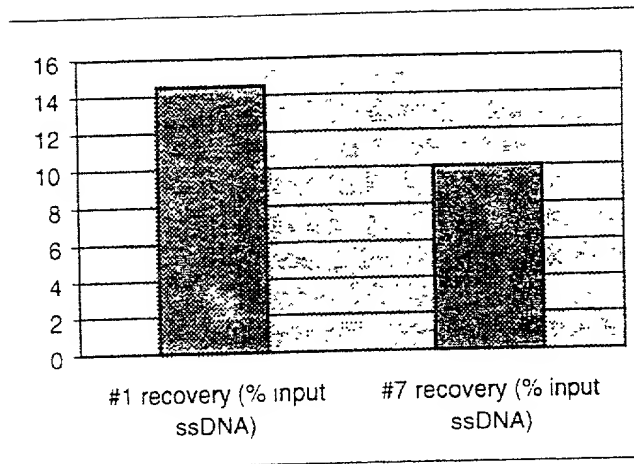

5

10 *Dedicated to the memory of Chris Griffin and Richard Beer*

201604091660

CLAIMS

- 5 1. A method of providing a mixture of DNA fragments enriched in fragments that are characteristic of a phenotype of interest, by providing affected DNA in fragmented form and providing unaffected DNA in fragmented form, which method comprises:
- a) mixing the fragments of the affected DNA and the fragments
 - 10 of the unaffected DNA under hybridising conditions;
 - b) recovering a mixture of hybrids that contain mismatches;
 - c) recovering fragments of the affected DNA from the mixture of hybrids that contain mismatches;
- 15 and optionally repeating steps a), b) and c) one or more times.
2. The method of claim 1 wherein the affected DNA is pooled DNA of individuals who show the phenotype of interest, and the unaffected DNA is pooled DNA of individuals who do not show the phenotype of
- 20 interest.
3. The method of claim 1, wherein the affected DNA is DNA of one individual who shows the phenotype of interest, and the unaffected DNA is pooled DNA of individuals who do not show the phenotype of
- 25 interest.
4. The method of claim 1, wherein the affected DNA is DNA of one individual who shows the phenotype of interest, and the unaffected DNA is pooled DNA of a complete set of ancestors who do not show the
- 30 phenotype of interest.

5

10

15

20

25

10. The method of any one of claims 1 to 9, wherein the mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest, is mixed with an excess of the fragments of the affected DNA under hybridisation conditions, followed by recovery of perfectly matched duplexes.

- 157 -

11. The method of any one of claims 1 to 10, wherein each of the affected DNA and the unaffected DNA is provided in fragmented form by digestion with from 4 to 7 six-cutter restriction endonuclease enzymes together with from 0 to 5 four-cutter restriction endonuclease enzymes.

5

12. A mixture of DNA fragments enriched in fragments that are characteristic of the phenotype of interest, provided by the method of any one of claims 1 to 11.

- 10 13. A method of making a set of arrays of fragments of DNA of interest, which method comprises:
- a) selecting, from a set of n restriction endonuclease enzymes, a subset of r restriction endonuclease enzymes;
 - b) digesting genomic DNA with the subset of r enzymes;
 - 15 c) ligating to the resulting fragments restriction-enzyme-cutting-site-specific adapters with unique polymerase chain reaction amplifiable sequences;
 - d) splitting the resulting fragments into r^2 aliquots;
 - e) amplifying each aliquot with two restriction-enzyme-specific
 - 20 primers of which one is tagged;
 - f) forming an array of the r^2 aliquots of non-tagged amplicon strands; and
 - g) repeating steps a) to f) using one or more different subsets of r restriction endonuclease enzymes.

25

09914604-013402

14. A method of making a set of arrays of fragments of DNA of interest, which method comprises:
- a) selecting, from a set of n restriction endonuclease enzymes, a subset of r restriction endonuclease enzymes;
 - 5 b) digesting genomic DNA with the subset of r enzymes;
 - c) ligating to the resulting fragments restriction-enzyme-cutting-site-specific adapters with unique polymerase chain reaction amplifiable sequences;
 - d) splitting the resulting fragments into r^2 aliquots;
 - 10 e) amplifying each aliquot with two restriction-enzyme-specific primers;
 - f) forming an array of the r^2 aliquots of the amplimer strands; and
 - g) repeating steps a) to f) using one or more different subsets of
 - 15 r restriction endonuclease enzymes.
15. The method of claim 13 or claim 14, wherein steps a to f) are repeated using each different subset of r restriction endonuclease enzymes to give $(n!)/[(n-r)!r!]$ different arrays.
- 20 16. The method of any one of claims 13 to 15, wherein the n restriction endonuclease enzymes are selected from 4-cutters and 5-cutters and 6-cutters.
- 25 17. The method of any one of claims 13 to 16, wherein n is 3 to 10 and r is 2 to 4.
18. The method of claim 17, where $n = 6$ and $r = 3$.
- 30 19. A set of arrays of fragments of DNA of interest, which set results from performance of the method of any one of claims 13 to 18.

20. The set of arrays of claim 19, which set results from performance of the method of claim 13 and claim 14 and claim 15.

5 21. The set of arrays of claim 19 or claim 20, derived from a set of $n = 6$ six-cutter restriction endonuclease enzymes which are *BamHI*; *Bsr GI*; *Hind III*; *NcoI*; *SpeI*; and *AflII*.

22. The set of arrays of claim 19 or claim 20, derived from the set
10 of $n = 6$ six-cutter restriction endonuclease enzymes which are *EcoRI*; *BspHI*; *BglII*; *XbaI*; *Acc65I*; and *ApaLI*.

23. A nucleic acid characterisation method which comprises
presenting to the set of arrays of any one of claims 19 to 22 a nucleic acid
15 fragment of interest under hybridisation conditions, and observing a pattern of hybridisation.

24. The method of claim 23, wherein a plurality of nucleic acid
fragments of interest are separately presented to the set of arrays, and the
20 resulting patterns of hybridisation are compared.

25. The method of claim 24, wherein the plurality of nucleic acid
fragments of interest are drawn from the mixture of DNA fragments,
enriched in fragments that are characteristic of a phenotype of interest, of
25 claim 13.

26. A method of identifying fragments of DNA that are
characteristic of a phenotype of interest, which method comprises
recovering, cloning and amplifying individual DNA fragments from the
30 mixture of DNA fragments of claim 12, presenting the individual DNA
fragments to the set of arrays of any one of claims 19 to 22 under

hybridisation conditions, observing a pattern of hybridisation generated by each individual DNA fragment, and subjecting to further investigation any two or more individual DNA fragments whose hybridisation patterns are similar or identical, or near to each other in a genome of interest.

5

27. A double-stranded DNA molecule having the sequence a-A-b-B...X-y-Y-z where A, B...X and Y are unique restriction sites for n different restriction endonuclease enzymes, and a, b...y, z denotes distances in base pairs, characterised in that each fragment, obtainable by cutting the
- 10 DNA molecule by means of any one or more up to n of the restriction enzymes, has a different length from every other fragment.

28. The double-stranded DNA molecule of claim 27, wherein the following criteria are satisfied:
- 15 a) inter-fragment length differences are greater for larger fragments;
- b) all possible fragments are unambiguously resolvable by electrophoresis from one another;
- c) size gaps between bands comprising different numbers of
- 20 inter-restriction-site units are larger than size gaps between bands comprising the same number of inter-restriction-site units;
- d) the size gaps and size spread from the largest to the smallest fragment are electrophoretically compatible.

2007-04-04 09:44:00

Figure 1


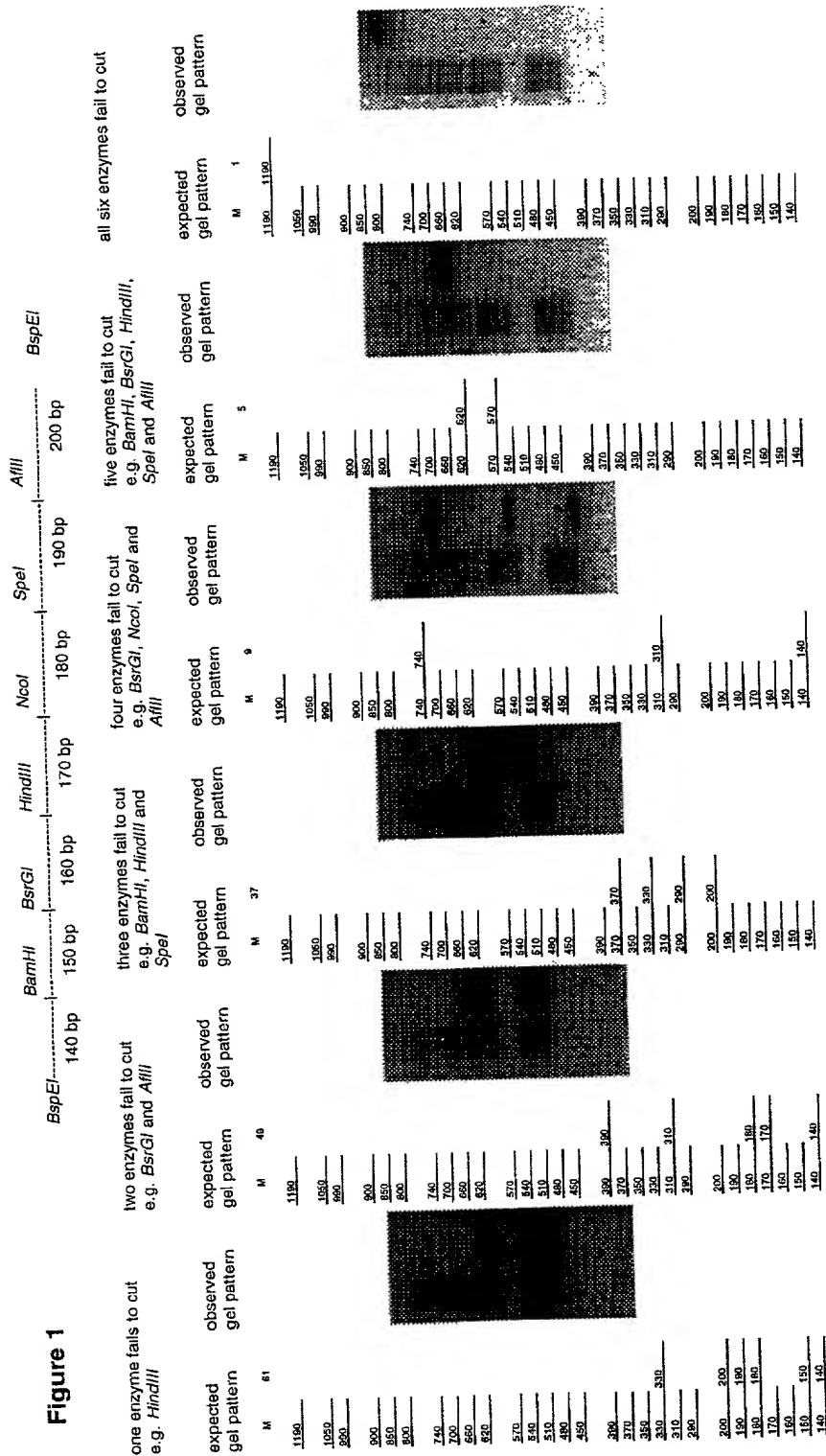


Figure 2

Vistra Green fluorescence

1 µg of Stratagene Kb DNA Ladder (band sizes: 250 bp, 500 bp, 750 bp, 1 kb, 1.5 kb, 2 kb, 3 kb, 4 kb, 5 kb, 6 kb, 7 kb, 8 kb, 9 kb, 10 kb and 12 kb), 3x gaps, 2 µg of canine genomic DNA digested with 0.25, 0.025, 0.0025, 0.00025, 0.000025, 0.0000025 and 0 U/µl *Bam*H1 / *Bsr*GI / *Hind*III / *Nco*I / *Spe*I / *Afl*II, gap, 0.4 µg of canine genomic DNA and 1 µl of *Bsp*EI-released internal control DNA digested with 0.25, 0.025, 0.0025, 0.00025, 0.000025, 0.0000025 and 0 U/µl *Bam*H1 / *Bsr*GI / *Hind*III / *Nco*I / *Spe*I / *Afl*II.

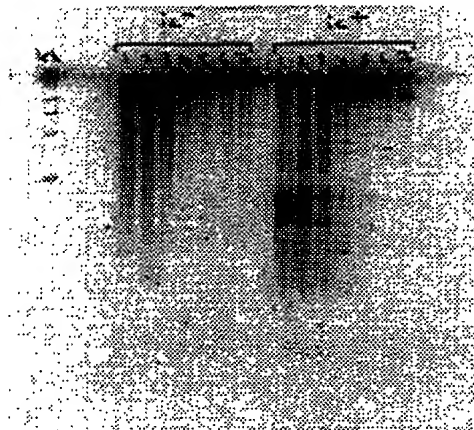

From the results above, the fragment pattern from the internal control DNA is diagnostic of both the degree of digestion and the nature of any partial restriction at less than limit digestion.

Very little digestion of canine genomic DNA and *Bsp*EI-released internal control DNA is observed for 0.00025 U/µl *Bam*H1 / *Bsr*GI / *Hind*III / *Nco*I / *Spe*I / *Afl*II and fewer units. With 0.0025 U/µl, *Bam*HI and *Afl*II are both failing to cut more than the other enzymes (the 140 bp and 200 bp bands in 3ic are diagnostic of this). With 0.025 U/µl, most enzymes have cut – though not quite to completion. 0.25 U/µl thus appears to be necessary in order to achieve limit digestion of 20 µg of canine genomic DNA in a volume of 200 µl.

Figure 3

Figure 3 shows the products of the restriction digests after electrophoresis on an 8 % polyacrylamide gel. The numbers 1–5 indicate the dilution of enzyme mix used. 1=0.5 U/ μ l, 2=0.1 U/ μ l, 3=0.02 U/ μ l, 4=0.004 U/ μ l, and 5=0 U/ μ l. 'ic' indicates the samples that were spiked with 4 bp cutter internal control DNA. PM indicates that 180 ng of PCR molecular weight markers were added as a size standard.



All five of the partial digestion products are visible in lane ic4, between the top 130 bp band and the bottom 40 bp band, indicating that all three enzymes have failed to cut to completion at 0.004 U/ μ l. In lane ic3, the four complete digestion products are all visible but the 65 bp band is also present. *HaeIII* and *MboI* have digested to completion but the *MseI* digestion is still partial at 0.2 U/ μ l. At 0.1 U/ μ l and at 0.5 U/ μ l all three enzymes have produced limit digests.

Using the internal control spike, it is a great deal easier to determine the point at which limit digestion has occurred. The undigested internal control DNA is 130 bp. Partial digestion products generate bands of 105 bp, 90 bp, 75 bp, 65 bp, and 55 bp. The products of complete digestion generate bands of 40 bp, 35 bp, 30 bp, and 25 bp. The fragment pattern from the internal control DNA is therefore diagnostic of both the degree of digestion and the nature of any partial restriction at less than limit digestion.

Figure 4

Expected TRSPA-2 pattern for the pNW33 *Bam*HI, *Hind*III, and *Afl*III matrix probed with the 140 bp *Bsp*EI to *Bam*HI fragment from pNW33

The TRSPA-2 pattern expected for matrix #7 (*Bam*HI, *Hind*III, and *Afl*III) hybridized with a 33 P-labelled probe derived from the 140 bp *Bsp*EI – *Bam*HI fragment of pNW33 is shown below:

#7	<i>Bam</i> HI	<i>Hind</i> III	<i>Afl</i> III
<i>Bam</i> HI	-	HYB	-
<i>Hind</i> III	HYB	-	-
<i>Afl</i> III	-	-	-

z

Observed TRSPA-2 pattern for the pNW33 *Bam*HI, *Hind*III, and *Afl*III matrix probed with the 140 bp *Bsp*EI to *Bam*HI fragment from pNW33

The TRSPA-2 pattern observed for matrix #7 (*Bam*HI, *Hind*III, and *Afl*III) hybridized with a 33 P-labelled probe derived from the 140 bp *Bsp*EI – *Bam*HI fragment of pNW33 is shown below:










#7	<i>Bam</i> HI	<i>Hind</i> III	<i>Afl</i> III
<i>Bam</i> HI			
<i>Hind</i> III			
<i>Afl</i> III			

Figure 5










Expected TRSPA-2 pattern for the pNW33 *HindIII*, *NcoI*, and *SpeI* matrix probed with the 140 bp *BspEI* to *BamHI* fragment from pNW33

The TRSPA-2 pattern expected for matrix #17 (*HindIII*, *NcoI*, and *SpeI*) hybridized with a ^{33}P -labelled probe derived from the 140 bp *BspEI* – *BamHI* fragment of pNW33 is shown below:

#17	<i>HindIII</i>	<i>NcoI</i>	<i>SpeI</i>
<i>HindIII</i>	HYB	HYB	HYB
<i>NcoI</i>	HYB	-	-
<i>SpeI</i>	HYB	-	-

Observed TRSPA-2 pattern for the pNW33 *HindIII*, *NcoI*, and *SpeI* matrix probed with the 140 bp *BspEI* to *BamHI* fragment from pNW33

The TRSPA-2 pattern expected for matrix #17 (*HindIII*, *NcoI*, and *SpeI*) hybridized with a ^{33}P -labelled probe derived from the 140 bp *BspEI* – *BamHI* fragment of pNW33 is shown below:

#17	<i>HindIII</i>	<i>NcoI</i>	<i>SpeI</i>
<i>HindIII</i>			
<i>NcoI</i>			
<i>SpeI</i>			

Please type a plus sign (+) inside this box → ☐

PTO/SB/01 (12-97)
Approved for use through 9/30/00. OMB 0651-0032
Patent and Trademark Office: U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it contains a valid OMB control number.

**DECLARATION FOR UTILITY OR
DESIGN
PATENT APPLICATION
(37 CFR 1.63)**

☐ Declaration Submitted with Initial Filing OR ☒ Declaration Submitted after Initial Filing (surcharge (37 CFR 1.16 (e)) required)

Attorney Docket Number	PA-9912
First Named Inventor	M. Reeve
COMPLETE IF KNOWN	
Application Number	09 / 914,604
Filing Date	28-Aug-2001
Group Art Unit	To be assigned
Examiner Name	To be assigned

As a below named inventor, I hereby declare that:

My residence, post office address, and citizenship are as stated below next to my name.

I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter which is claimed and for which a patent is sought on the invention entitled:

Genetic Analysis

the specification of which (Title of the Invention)

☐ is attached hereto
OR

☒ was filed on (MM/DD/YYYY) 08/28/2001 as United States Application Number or PCT International

Application Number 09/914,604 and was amended on (MM/DD/YYYY) (if applicable).

I hereby state that I have reviewed and understand the contents of the above identified specification, including the claims, as amended by any amendment specifically referred to above.

I acknowledge the duty to disclose information which is material to patentability as defined in 37 CFR 1.56.

I hereby claim foreign priority benefits under 35 U.S.C. 119(a)-(d) or 365(b) of any foreign application(s) for patent or inventor's certificate, or 365(a) of any PCT international application which designated at least one country other than the United States of America, listed below and have also identified below, by checking the box, any foreign application for patent or inventor's certificate, or of any PCT international application having a filing date before that of the application on which priority is claimed.

Prior Foreign Application Number(s)	Country	Foreign Filing Date (MM/DD/YYYY)	Priority Not Claimed	Certified Copy Attached?	
				YES	NO
99301933.0	EPO	03/12/1999	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

☐ Additional foreign application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto:

I hereby claim the benefit under 35 U.S.C. 119(e) of any United States provisional application(s) listed below.

Application Number(s)	Filing Date (MM/DD/YYYY)	<input type="checkbox"/> Additional provisional application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto.

[Page 1 of 2]

Burden Hour Statement: This form is estimated to take 0.4 hours to complete. Time will vary depending upon the needs of the individual case. Any comments on the amount of time you are required to complete this form should be sent to the Chief Information Officer, Patent and Trademark Office, Washington, DC 20231. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Assistant Commissioner for Patents, Washington, DC 20231.

Please type a plus sign (+) inside this box → ☐

PTO/SB/02A (3-97)
Approved for use through 9/30/98. OMB 0651-0032
Patent and Trademark Office: U.S. DEPARTMENT OF COMMERCE
Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it contains a valid OMB control number.

+

DECLARATION

ADDITIONAL INVENTOR(S)
Supplemental Sheet
Page 1 of 1

Name of Additional Joint Inventor, if any:

☐ A petition has been filed for this unsigned inventor

Given Name (first and middle (if any))

Family Name or Surname

Nicholas Ian

Workman

Inventor's
Signature

N. Z. Workman

Date

16 OCT
2001

Residence: City

State

Country

GB

Citizenship

GB

Post Office Address

4 Cranford Court, Aylesbury

Post Office Address

Buckinghamshire, Great Britain HP21 7NZ GBN

City

State

ZIP

Country

Name of Additional Joint Inventor, if any:

☐ A petition has been filed for this unsigned inventor

Given Name (first and middle (if any))

Family Name or Surname

Luis

Martin-Parras

Inventor's
Signature

Luis Martin-Parras

Date

16 OCT
2001

Residence: City

State

Country

GB

Citizenship

ES

Post Office Address

7 Gladstone Rise, High Wycombe

Post Office Address

Buckinghamshire, Great Britain HP13 7NW GBN

City

State

ZIP

Country

Name of Additional Joint Inventor, if any:

☐ A petition has been filed for this unsigned inventor

Given Name (first and middle (if any))

Family Name or Surname

Inventor's
Signature

Date

Residence: City

State

Country

Citizenship

Post Office Address

Post Office Address

City

State

ZIP

Country

Burden Hour Statement: This form is estimated to take 0.4 hours to complete. Time will vary depending upon the needs of the individual case. Any comments on the amount of time you are required to complete this form should be sent to the Chief Information Officer, Patent and Trademark Office, Washington, DC 20231. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Assistant Commissioner for Patents, Washington, DC 20231.

+

Please type a plus sign (+) inside this box → ☐

PTO/SB/01 (12-97)
Approved for use through 9/30/00. OMB 0651-0032
Patent and Trademark Office: U.S. DEPARTMENT OF COMMERCE
Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it contains a valid OMB control number.

DECLARATION — Utility or Design Patent Application

I hereby claim the benefit under 35 U.S.C. 120 of any United States application(s), or 365(c) of any PCT international application designating the United States of America, listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in the prior United States or PCT International application in the manner provided by the first paragraph of 35 U.S.C. 112, I acknowledge the duty to disclose information which is material to patentability as defined in 37 CFR 1.56 which became available between the filing date of the prior application and the national or PCT international filing date of this application.

U.S. Parent Application or PCT Parent Number	Parent Filing Date (MM/DD/YYYY)	Parent Patent Number (if applicable)
PCT/GB00/00916	03/10/2000	

☐ Additional U.S. or PCT international application numbers are listed on a supplemental priority data sheet PTO/SB/01A attached hereto.
As a named inventor, I hereby appoint the following registered practitioner(s) to prosecute this application and to transact all business in the Patent and Trademark Office connected therewith: ☒ Customer Number 22840 OR ☐ Registered practitioner(s) name/registration number listed below

Name	Registration Number	Name	PATENT TRADEMARK OFFICE Registration Number

☐ Additional registered practitioner(s) named on supplemental Registered Practitioner Information sheet PTO/SB/02C attached hereto.
Direct all correspondence to: ☒ Customer Number or Bar Code Label 22840 OR ☐ Correspondence address below

Name					
Address					
Address					
City		State		ZIP	
Country		Telephone		Fax	

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under 18 U.S.C. 1001 and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

1-00

Name of Sole or First Inventor:		<input type="checkbox"/> A petition has been filed for this unsigned inventor			
Given Name (first and middle (if any))		Family Name or Surname			
Michael Alan		Reeve			
Inventor's Signature				Date	16 Oct 2001
Residence: City	State	Country	GB	Citizenship	GB
Post Office Address	118 St. Andrews Road, Henley-on-Thames				
Post Office Address	Oxfordshire, Great Britain RG9 1PL GBN				
City	State	ZIP		Country	

☒ Additional inventors are being named on the 1 supplemental Additional Inventor(s) sheet(s) PTO/SB/02A attached hereto